

COMPREHENSIVE SERVICES

We offer competitive repair and calibration services, as well as easily accessible documentation and free downloadable resources.

SELL YOUR SURPLUS

We buy new, used, decommissioned, and surplus parts from every NI series. We work out the best solution to suit your individual needs.

 Sell For Cash  Get Credit  Receive a Trade-In Deal

OBSOLETE NI HARDWARE IN STOCK & READY TO SHIP

We stock **New**, **New Surplus**, **Refurbished**, and **Reconditioned** NI Hardware.



Bridging the gap between the manufacturer and your legacy test system.

 1-800-915-6216

 www.apexwaves.com

 sales@apexwaves.com

All trademarks, brands, and brand names are the property of their respective owners.

Request a Quote

 **CLICK HERE**


GPIB-1284CT

The LabVIEW logo is displayed vertically within a rectangular frame. The word "LabVIEW" is written in a bold, serif font, with the 'L' and 'V' being significantly larger than the other letters. The frame has a thin black border and a white background.

LabVIEW

LabVIEW[®] Instrument I/O VI Reference Manual

January 1996 Edition
Part Number 320537C-01

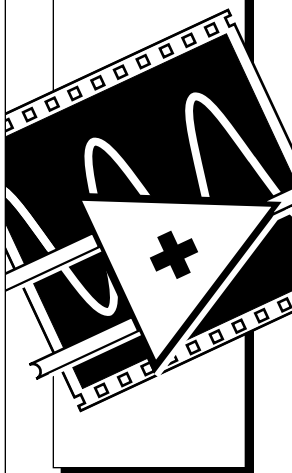

© Copyright 1992, 1995 National Instruments Corporation. All Rights Reserved.



Internet Support

GPIB: gpib.support@natinst.com

DAQ: daq.support@natinst.com

VXI: vxi.support@natinst.com

LabVIEW: lv.support@natinst.com

LabWindows: lw.support@natinst.com

HiQ: hiq.support@natinst.com

E-mail: info@natinst.com

FTP Site: [ftp.natinst.com](ftp://ftp.natinst.com)

Web Address: <http://www.natinst.com>



Bulletin Board Support

BBS United States: (512) 794-5422 or (800) 327-3077

BBS United Kingdom: 01635 551422

BBS France: 1 48 65 15 59



FaxBack Support

(512) 418-1111 or (800) 329-7177



Telephone Support (U.S.)

Tel: (512) 795-8248

Fax: (512) 794-5678 or (800) 328-2203



International Offices

Australia 03 9 879 9422, Austria 0662 45 79 90 0, Belgium 02 757 00 20,
Canada (Ontario) 519 622 9310, Canada (Québec) 514 694 8521, Denmark 45 76 26 00,
Finland 90 527 2321, France 1 48 14 24 24, Germany 089 741 31 30, Hong Kong 2645 3186,
Italy 02 48301892, Japan 03 5472 2970, Korea 02 596 7456, Mexico 95 800 010 0793,
Netherlands 0348 433466, Norway 32 84 84 00, Singapore 2265886, Spain 91 640 0085,
Sweden 08 730 49 70, Switzerland 056 200 51 51, Taiwan 02 377 1200, U.K. 01635 523545

National Instruments Corporate Headquarters

6504 Bridge Point Parkway Austin, TX 78730-5039 Tel: (512) 794-0100

Important Information

Warranty

The media on which you receive National Instruments software are warranted not to fail to execute programming instructions, due to defects in materials and workmanship, for a period of 90 days from date of shipment, as evidenced by receipts or other documentation. National Instruments will, at its option, repair or replace software media that do not execute programming instructions if National Instruments receives notice of such defects during the warranty period. National Instruments does not warrant that the operation of the software shall be uninterrupted or error free.

A Return Material Authorization (RMA) number must be obtained from the factory and clearly marked on the outside of the package before any equipment will be accepted for warranty work. National Instruments will pay the shipping costs of returning to the owner parts which are covered by warranty.

National Instruments believes that the information in this manual is accurate. The document has been carefully reviewed for technical accuracy. In the event that technical or typographical errors exist, National Instruments reserves the right to make changes to subsequent editions of this document without prior notice to holders of this edition. The reader should consult National Instruments if errors are suspected. In no event shall National Instruments be liable for any damages arising out of or related to this document or the information contained in it.

EXCEPT AS SPECIFIED HEREIN, NATIONAL INSTRUMENTS MAKES NO WARRANTIES, EXPRESS OR IMPLIED, AND SPECIFICALLY DISCLAIMS ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE. CUSTOMER'S RIGHT TO RECOVER DAMAGES CAUSED BY FAULT OR NEGLIGENCE ON THE PART OF NATIONAL INSTRUMENTS SHALL BE LIMITED TO THE AMOUNT THEREFORE PAID BY THE CUSTOMER. NATIONAL INSTRUMENTS WILL NOT BE LIABLE FOR DAMAGES RESULTING FROM LOSS OF DATA, PROFITS, USE OF PRODUCTS, OR INCIDENTAL OR CONSEQUENTIAL DAMAGES, EVEN IF ADVISED OF THE POSSIBILITY THEREOF. This limitation of the liability of National Instruments will apply regardless of the form of action, whether in contract or tort, including negligence. Any action against National Instruments must be brought within one year after the cause of action accrues. National Instruments shall not be liable for any delay in performance due to causes beyond its reasonable control. The warranty provided herein does not cover damages, defects, malfunctions, or service failures caused by owner's failure to follow the National Instruments installation, operation, or maintenance instructions; owner's modification of the product; owner's abuse, misuse, or negligent acts; and power failure or surges, fire, flood, accident, actions of third parties, or other events outside reasonable control.

Copyright

Under the copyright laws, this publication may not be reproduced or transmitted in any form, electronic or mechanical, including photocopying, recording, storing in an information retrieval system, or translating, in whole or in part, without the prior written consent of National Instruments Corporation.


Trademarks

LabVIEW[®], NAT4882[®], NI-488[®], and NI-488.2[™] are trademarks of National Instruments Corporation.

Product and company names listed are trademarks or trade names of their respective companies.

WARNING REGARDING MEDICAL AND CLINICAL USE OF NATIONAL INSTRUMENTS PRODUCTS

National Instruments products are not designed with components and testing intended to ensure a level of reliability suitable for use in treatment and diagnosis of humans. Applications of National Instruments products involving medical or clinical treatment can create a potential for accidental injury caused by product failure, or by errors on the part of the user or application designer. Any use or application of National Instruments products for or involving medical or clinical treatment must be performed by properly trained and qualified medical personnel, and all traditional medical safeguards, equipment, and procedures that are appropriate in the particular situation to prevent serious injury or death should always continue to be used when National Instruments products are being used. National Instruments products are NOT intended to be a substitute for any form of established process, procedure, or equipment used to monitor or safeguard human health and safety in medical or clinical treatment.

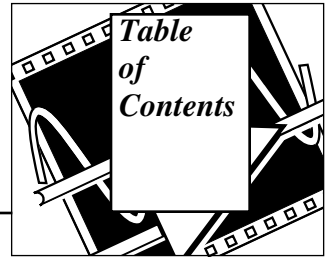

About This Manual

Organization of This Manual	xv
Conventions Used in This Manual	xvi
Related Documentation	xviii
Customer Communication	xix

Chapter 1 Introduction

Instrument Driver Overview	1-1
Introduction to Instrument Drivers	1-1
Instrument Driver Distribution	1-1
Instrument Driver Bulletin Board Access	1-2
Instrument Driver Internet Access	1-2
CD-ROM Instrument Driver Distribution	1-2
Instrument Driver Template VIs	1-3
VISA Library Reference Overview	1-4
VISA Functions	1-4
VISA Attribute Node	1-6
GPIOB Function Overview	1-9
Introduction to the GPIOB	1-9
History of the GPIOB	1-9
The IEEE 488.2 Standard	1-9
Compatible GPIOB Hardware	1-10
LabVIEW for Windows 3.x	1-10
LabVIEW for Windows NT	1-11
LabVIEW for HP-UX	1-11
LabVIEW for Sun	1-11
LabVIEW for Macintosh	1-11
Traditional GPIOB Functions	1-12

LabVIEW GPIB 488.2 Functions	1-13
Single-Device Functions	1-14
Multiple-Device Functions	1-15
Bus Management Functions	1-15
Low-Level Functions	1-16
General Functions	1-17
Serial Port VI Overview	1-17

Chapter 2

Overview of LabVIEW Instrument Drivers

Introduction	2-1
LabVIEW Instrument Driver External Interface Model	2-2
Functional Body	2-2
Interactive Developer Interface	2-3
Programmatic Developer Interface	2-3
I/O Interface	2-3
Subroutine Interface	2-4
LabVIEW Instrument Driver Internal Design Model	2-4
Instrument Driver Application VIs	2-6
Instrument Driver Component VIs	2-6
Error Reporting	2-9

Chapter 3

Developing a LabVIEW Instrument Driver

Development Procedure	3-1
Designing the Instrument Driver Structure	3-1
Instrument Driver Structure and VI Hierarchy	3-2
Instrument Driver VI Organization	3-3
Guidelines and Recommendations	3-4
Design Example	3-6
Modifying the Instrument Driver Templates	3-8
Adding Instrument Driver Component VI VIs	3-9
Tips for Developing a LabVIEW Instrument Driver	3-10
Loop Termination Conditions	3-10
Assembling Command Strings	3-11
Data Dependency	3-13
Programmatically Checking Parameter Ranges	3-13
Guidelines	3-16
Front Panel	3-16
Required Front Panel Controls	3-16
Control Guidelines	3-17

Block Diagram	3-19
Icon	3-20
Connector Pane	3-20
Error Reporting	3-21
Online Help Information	3-21
VI Descriptions	3-22
Control and Indicator Descriptions	3-22
Application VIs	3-23
LabVIEW Instrument Driver Standards Checklist	3-24
General Issues	3-24
Files Created	3-25
LabVIEW Instrument Driver Archive - .LLB	3-25
VI Front Panels	3-25
Icon/Connector Pane/VI Setup	3-26
Block Diagram	3-26

Chapter 4

Instrument Driver Template VIs

Instrument Driver Template VI Overview	4-1
Initialize	4-2
Close	4-3
Reset	4-3
Self-Test	4-3
Error Query and Error Message	4-3
Revision Query	4-3
Message-Based Template and Register-Based Template	4-4
Example	4-4
Instrument Driver Template VI Error Reporting	4-5
Instrument Driver Template VI Descriptions	4-6
PREFIX Close	4-6
PREFIX Error Message	4-7
PREFIX Error Query	4-9
PREFIX Initialize	4-10
PREFIX Message-Based Template	4-13
PREFIX Register-Based Template	4-14
PREFIX Reset	4-16
PREFIX Revision Query	4-17
PREFIX Self-Test	4-19

Chapter 5

VISA Library Reference

Operations	5-2
VISA Library Reference Parameters	5-2
VISA Operation Descriptions	5-3
VISA Assert Trigger	5-3
VISA Clear	5-4
VISA Close	5-4
VISA Find Resource	5-5
VISA Lock	5-6
VISA Open	5-7
VISA Read	5-8
VISA Read STB	5-9
VISA Status Description	5-10
VISA Unlock	5-11
VISA Write	5-11
Event Handling Functions	5-12
VISA Disable Event	5-12
VISA Discard Events	5-13
VISA Enable Event	5-14
VISA Wait On Event	5-15
High Level Register Access Functions	5-16
VISA In8 / In16 / In32	5-16
VISA Memory Allocation	5-18
VISA Memory Free	5-18
VISA Move In8 / Move In16 / Move In32	5-19
VISA Move Out8 / Move Out16 / Move Out32	5-20
VISA Out8 / Out16 / Out32	5-21
Low Level Register Access Functions	5-22
VISA Map Address	5-23
VISA Memory Allocation	5-23
VISA Memory Free	5-24
VISA Peek8 / Peek16 / Peek32	5-24
VISA Poke8 / Poke16 / Poke32	5-25
VISA Unmap Address	5-25
VISA Attribute Node	5-26
VISA Attribute Descriptions	5-27
Fast Data Channel Mode	5-27
Fast Data Channel Number	5-27
Fast Data Channel Pairs	5-27
Fast Data Channel Signal Enable	5-28
GPIO Primary Address	5-28

GPIB Secondary Address	5-28
IO Protocol	5-28
Immediate Servant	5-29
Increment Destination Count	5-29
Increment Source Count	5-29
Interface Number	5-30
InterfaceType	5-30
Mainframe Logical Address	5-30
Manufacturer ID	5-30
Maximum Queue Length	5-31
Model Code	5-31
Resource Lock State	5-31
Resource Manufacturer Identification	5-31
Resource Manufacturer Name	5-32
Resource Name	5-32
Send End Enable	5-33
Slot	5-33
Suppress End Enable	5-34
Termination Character	5-34
Termination Character Enable	5-34
Timeout Value	5-34
Trigger Identifier	5-35
User Data	5-35
Version of Implementation	5-36
Version of Specification	5-36
VXI Commander Logical Address	5-36
VXI Logical Address	5-37
VXI Memory Address Space	5-37
VXI Memory Base Address	5-37
VXI Memory Size	5-37
Window Access	5-38
Window Base Address	5-38
Window Size	5-38

Chapter 6

Traditional GPIB Functions

Traditional GPIB Function Parameters	6-2
Traditional GPIB Function Behavior	6-3
Traditional GPIB Function Descriptions	6-4
GPIB Clear	6-4
GPIB Initialization	6-4
GPIB Misc	6-5

GPIB Read	6-7
GPIB Serial Poll	6-9
GPIB Status	6-9
GPIB Trigger	6-12
GPIB Wait	6-13
Wait for GPIB RQS	6-15
GPIB Write	6-16
GPIB Device and Controller Functions	6-17
Device Functions	6-17
loc – Go to local	6-17
off– Take device offline	6-18
pct – Pass control	6-18
ppc – Parallel poll configure	6-18
Controller Functions	6-19
cac – Become active Controller	6-19
cmd – Send IEEE 488 commands	6-20
dma – Set DMA mode or programmed I/O mode	6-20
gts – Go from active Controller to standby	6-20
ist – Set individual status bit	6-21
llo – Local lockout	6-22
loc – Place Controller in local state	6-22
off– Take controller offline	6-22
ppc – Parallel poll configure (enable and disable)	6-22
ppu – Parallel poll unconfigure	6-23
rpp – Conduct parallel poll	6-23
rsc – Release or request system control	6-23
rsv – Request service and/or set the serial poll status byte	6-24
sic – Send interface clear	6-24
sre – Unassert or assert remote enable	6-25

Chapter 7

GPIB 488.2 Functions

GPIB 488.2 Common Function Parameters	7-2
GPIB 488.2 Function Descriptions	7-3
Single-Device Functions	7-3
DevClear	7-3
PPollConfig	7-4
PassControl	7-5
ReadStatus	7-5
Receive	7-6
Send	7-7
Trigger	7-8

Multiple-Device Functions	7-9
AllSPoll	7-9
DevClearList	7-10
EnableLocal	7-11
EnableRemote	7-11
FindRQS	7-12
PPoll	7-13
PPollUnconfig	7-14
SendList	7-14
TriggerList	7-15
Bus Management Functions	7-16
FindLstn	7-16
ResetSys	7-17
SendIFC	7-18
SendLLO	7-18
SetRWLS	7-19
TestSRQ	7-20
TestSys	7-21
WaitSRQ	7-22
Low-Level I/O Functions	7-23
RcvRespMsg	7-23
ReceiveSetup	7-24
SendCmds	7-25
SendDataBytes	7-25
SendSetup	7-26
General Functions	7-27
MakeAddr	7-27
SetTimeOut	7-28

Chapter 8

Serial Port VIs

Common Serial Port VI Parameters	8-2
Port Number	8-2
Handshaking Modes	8-2
Software Handshaking–XON/XOFF	8-3
Error Codes	8-3
Serial Port VI Descriptions	8-3
Bytes at Serial Port	8-3
Serial Port Break	8-4
Serial Port Init	8-5
Serial Port Read	8-7
Serial Port Write	8-7

Appendix A

Multiline Interface Messages

Appendix B

Operation of the GPIB

Types of Messages	B-1
Talkers, Listeners, and Controllers	B-1
The Controller-In-Charge and System Controller	B-2
GPIB Signals and Lines	B-2
Physical and Electrical Characteristics	B-4
Configuration Requirements	B-7

Appendix C

Instrument I/O VI Common Questions

GPIB	C-1
Serial I/O	C-4

Appendix D

Error Codes

VISA Error Codes	D-1
Instrument Driver Error Codes	D-4
GPIB Error Codes	D-5
Serial Port Error Codes	D-7

Appendix E

Customer Communication

Glossary

Index


Figures

Figure 2-1. VIs in Tek VX4790 Example Diagram	2-8
Figure 2-2. Tek VX4790 Config Std Wave Diagram	2-9
Figure 3-1. Range Test VI (Front Panel and Block Diagram)	3-14
Figure 3-2. Simple Trigger VI with Programmatic Range Testing	3-15
Figure 3-3. Simple Trigger VI without Programmatic Range Testing	3-15
Figure 3-4. Simple Trigger VI Front Panel (See Figure 3-3 for Diagram)	3-19

Figure B-1. GPIB Connector Showing Signal Assignment	B-5
Figure B-2. Linear Configuration	B-6
Figure B-3. Star Configuration	B-7

Tables

Table 6-1. Command String Functions	6-6
Table 6-2. Status Bits	6-10
Table 6-3. GPIB Error Codes	6-11
Table 6-4. Wait State Vector Bits	6-14
Table D-1. VISA Error Codes	D-1
Table D-2. Instrument Driver Error Codes	D-4
Table D-3. GPIB Error Codes	D-5
Table D-4. Serial Port Error Codes	D-7

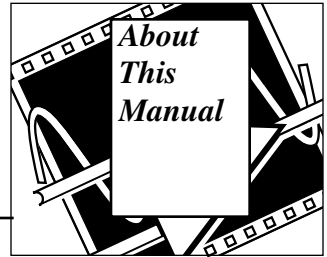

The *LabVIEW Instrument I/O VI Reference Manual* contains descriptions for the use and development of LabVIEW instrument drivers and LabVIEW virtual instruments (VIs) for GPIB, VXIbus, and serial port I/O operation.

This manual is a supplement to your LabVIEW User Manual and assumes that you are familiar with that material. You should also know how to operate LabVIEW, and your computer and operating system.

Organization of This Manual

The *LabVIEW Instrument I/O VI Reference Manual* is organized as follows:

- Chapter 1, *Introduction*, contains an overview of LabVIEW instrument drivers and the GPIB, serial port, instrument driver template and VISA transition library VIs. It also contains a history of the GPIB, and an explanation of GPIB improvements and standards.
- Chapter 2, *Overview of LabVIEW Instrument Drivers*, contains an overview of the LabVIEW instrument driver model, VIs, and VI templates.
- Chapter 3, *Developing a LabVIEW Instrument Driver*, describes the procedure for developing a LabVIEW instrument driver.
- Chapter 4, *Instrument Driver Template VIs*, describes the Instrument Driver Template VIs.
- Chapter 5, *VISA Library Reference*, describes the VISA (Virtual Instrument Software Architecture) Transition Library.
- Chapter 6, *Traditional GPIB Functions*, describes the traditional GPIB Functions. You use these functions to handle situations outside the scope of the IEEE 488.2 standard.
- Chapter 7, *GPIB 488.2 Functions*, describes the IEEE 488.2 (GPIB) Functions.

- Chapter 8, *Serial Port VIs*, describes the VIs for serial port operations.
- Appendix A, *Multiline Interface Messages*, lists multiline interface messages, which are commands that IEEE 488 defines.
- Appendix B, *Operation of the GPIB*, describes basic concepts you need to understand to operate the GPIB.
- Appendix C, *Instrument I/O VI Common Questions*, lists answers to questions frequently asked by LabVIEW users.
- Appendix D, *Error Codes*, lists error codes returned by LabVIEW functions, including the name and description of each code.
- Appendix E, *Customer Communication*, contains forms you can use to request help from National Instruments or to comment on our products and manuals.
- The *Glossary* contains an alphabetical list and description of terms used in this manual, including abbreviations, acronyms, metric prefixes, mnemonics, and symbols.
- The *Index* contains an alphabetical list of key terms and topics in this manual, including the page where you can find each one.

Conventions Used in This Manual

	The following conventions are used in this manual:
bold	Bold text denotes a parameter, menu item, return value, function panel item, or dialog box button or option.
<i>italic</i>	Italic text denotes emphasis, a cross reference, or an introduction to a key concept.
<i>bold italic</i>	Bold italic text denotes a note, caution, or warning.
monospace	Text in this font denotes text or characters that you should literally enter from the keyboard. Sections of code, programming examples, and syntax examples also appear in this font. This font also is used for the proper names of disk drives, paths, directories, programs, subprograms, subroutines, device names, variables, filenames, and extensions, and for statements and comments taken from program code.
<i>italic monospace</i>	Italic text in this font denotes that you must supply the appropriate words or values in the place of these items.
<>	Angle brackets enclose the name of a key. A hyphen between two or more key names enclosed in angle brackets denotes that you should

simultaneously press the named keys—for example, <Ctrl-Alt-Delete>.

<enter> Key names are lowercase

» The » symbol leads you through nested menu items and dialog box options to a final action. The sequence

File»Page Setup»Options»Substitute Fonts

directs you to pull down the **File** menu, select the **Page Setup** item, select **Options**, and finally select the **Substitute Fonts** option from the last dialog box.

paths Paths in this manual are denoted using backslashes (\) to separate drive names, directories, and files, as in
drivename\dir1name\dir2name\myfile

IEEE 488.1 IEEE 488.1 and IEEE 488.2 refer to the ANSI/IEEE Standard
and 488.1-1987 and the ANSI/IEEE Standard 488.2-1987, respectively,
IEEE 488.2 which define the GPIB.



Warning: *This icon to the left of bold italicized text denotes a warning, which alerts you to the possibility of damage to you or your equipment.*



Caution: *This icon to the left of bold italicized text denotes a caution, which alerts you to the possibility of data loss or a system crash.*



Note: *This icon to the left of bold italicized text denotes a note, which alerts you to important information.*

Data Types Each VI description gives a data type picture for each control and indicator, as illustrated in the following table.

Control	Indicator	Data Type
		Boolean
		String
		Unsigned 32-Bit Integer
		Unsigned 16-Bit Integer
		Unsigned 8-Bit Integer

Control	Indicator	Data Type
		Signed 32-Bit Integer
		16-Bit Integer (Word)
		Cluster of Clusters or Cluster of Arrays
		Array of 16-Bit Integers (Word)
		Array of Unsigned 8-bit integers
		Boolean Array
		String array
		Refnum
		Double-Precision Floating Point Number
		Array of Double-Precision Floating Point Numbers

Abbreviations, acronyms, metric prefixes, mnemonics, symbols, and terms are listed in the *Glossary*.

Related Documentation


The following documents contain information that you may find helpful as you read this manual:

- *LabVIEW Tutorial*
- *LabVIEW User Manual*
- *LabVIEW Function Reference Manual*
- *LabVIEW VXI VI Reference Manual*
- *IEEE Standard Digital Interface for Programmable Instrumentation*, ANSI/IEEE Standard 488.1-1987
- *IEEE Standard Codes, Formats, Protocols, and Common Commands*, ANSI/IEEE Standard 488.2-1987

Customer Communication

National Instruments wants to receive your comments on our products and manuals. We are interested in the applications you develop with our products, and we want to help if you have problems with them. To make it easy for you to contact us, this manual contains comment and configuration forms for you to complete. These forms are in Appendix E, *Customer Communication*, at the end of this manual.

Introduction

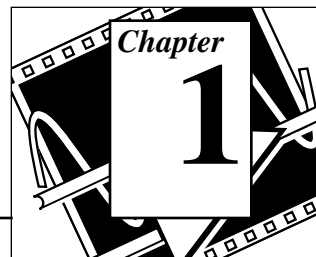

This chapter contains an overview of LabVIEW instrument drivers and the GPIB, serial port, instrument driver template and VISA VIs and functions. It also contains a history of the GPIB, and an explanation of GPIB improvements and standards.

Instrument Driver Overview

Introduction to Instrument Drivers

A LabVIEW instrument driver is a set of VIs that control a programmable instrument. Each VI corresponds to a programmatic operation such as configuring, reading from, writing to, or triggering the instrument. LabVIEW instrument drivers simplify instrument control and reduce test program development time by eliminating the need to learn the low-level programming protocol for each instrument.

The LabVIEW instrument driver library contains instrument drivers for a variety of programmable instrumentation, including GPIB, VXI, and serial. If a driver for your instrument is in the library, you can use it as is to control your instrument. Instrument drivers are distributed with their block diagram source code, so you can customize them for your specific application. If a driver for your particular instrument does not exist, you can:

- Try using a driver for a similar instrument. Often similar instruments from the same manufacturer have similar if not identical instrument drivers.
- Modify the Instrument Driver Template VIs to create a new driver for your instrument.
- Use either the GPIB, VXI, Serial, or VISA I/O libraries provided with LabVIEW to send commands directly to your instrument.

Instrument Driver Distribution

LabVIEW instrument drivers are distributed in a variety of media including electronic via bulletin board and internet and CD-ROM.

Instrument Driver Bulletin Board Access

You can download the latest versions of the LabVIEW instrument drivers from one of the National Instruments bulletin boards. The bulletin boards are located in National Instruments offices in the United States, United Kingdom, and France. All of the bulletin boards can auto baud to 9,600 baud and use 8 data bits, 1 stop bit, and no parity. The United States bulletin board handles up to 14,400 baud. Telephone numbers for the bulletin boards are:

United States	(512) 794-5422
United Kingdom	(0635) 551422
Outside the United Kingdom	(44) 635 551422
France	(1) 48 65 15 59
Outside France	(33) 1 48 65 15 59

Instrument Driver Internet Access

If you have internet access, you can download the latest instrument driver files from the National Instrument File Transfer Protocol (FTP) site. Use FTP to log on to host `ftp.natinst.com`. Use anonymous as the user name and your e-mail address as the password. The instrument drivers are in the following locations:

LabVIEW for Windows

`/support/labview/windows/instruments/`

LabVIEW for Sun

`/support/labview/sun/instruments/`

LabVIEW for Macintosh

`/support/labview/mac/instruments/`

CD-ROM Instrument Driver Distribution

The entire library of LabVIEW instrument drivers is available on CD-ROM. The instrument driver CD-ROM is available from National Instruments at no charge.

You can retrieve the latest instrument driver list on a touch-tone phone by calling the National Instruments automated fax system, Fax Back, at (512) 418-1111 or by calling National Instruments.

Instrument Driver Template VIs

The LabVIEW instrument driver templates are the foundation for all LabVIEW instrument driver development. The templates have a simple, flexible structure and a common set of instrument driver VIs that you can use for driver development. The VIs establish a standard format for all LabVIEW drivers and each has instructions for modifying it for a particular instrument.

The LabVIEW instrument driver templates are predefined instrument driver VIs that perform common operations such as initialization, self- test, reset, error query, and so on. Instead of developing your own VIs to accomplish these tasks, you should use the LabVIEW instrument driver template VIs, which already conform to the LabVIEW standards for instrument drivers.

The following table lists the Instrument Driver Template VIs:

Instrument Driver Template VI	Description
PREFIX Init	Establishes communication with the instrument and performs user selectable identification query and reset operations.
PREFIX Close	Terminates communication with the instrument.
PREFIX Revision	Queries the instrument driver and instrument firmware revisions.
PREFIX Self-Test	Performs an instrument self-test and returns the result.
PREFIX Reset	Places the instrument in a known state.
PREFIX Error Query	Queries the instrument for errors.
PREFIX Error Message	Converts error numbers from other instrument driver functions to user-readable strings.
PREFIX Message-Based Template	A template for creating instrument driver VIs for a message-based instrument.
PREFIX Register-Based Template	A template for creating instrument driver VIs for a register-based instrument.

VISA Library Reference Overview

VISA Functions

VISA (Virtual Instrument Software Architecture) is a single interface library for controlling VXI, GPIB, RS-232, and other types of instruments. The VISA Library provides a standard set of I/O routines used by all LabVIEW instrument drivers. Using the VISA functions, you can construct a single instrument driver VI which controls a particular instrument model across different I/O interfaces.

An instrument descriptor string is passed to the VISA Open function in order to select which kind of I/O will be used to communicate with the instrument. Once the session with the instrument is open, functions such as VISA Read and VISA Write perform the instrument I/O activities in a generic manner such that the program is not tied to any specific GPIB or VXI functions. Such an instrument driver is considered to be interface independent and can be used as is in different systems.

Instrument drivers which use the VISA functions perform activities specific to the instrument, not to the communication interface. This creates more opportunities for using the instrument driver in many diverse situations.

The following table lists the VISA functions:

VISA Function	Description
VISA Assert Trigger	Asserts a software or hardware trigger.
VISA Clear	Performs an IEEE 488.1-style clear of the device.
VISA Close	Closes a specified device session.
VISA Disable Event	Disables servicing of an event. Event occurrences already queued are not lost.
VISA Discard Events	Discards all pending occurrences of the specified event types and mechanisms from the specified session.
VISA Enable Event	Enables notification of a specified event.

VISA Function	Description
VISA Find Resource	Queries the system to locate the devices associated with a specified interface.
VISA In 8/16/32	Reads in 8-bits, 16-bits, or 32-bits of data, respectively, from the specified memory space.
VISA Lock	Establishes exclusive access to the specified resource.
VISA Map Address	Maps in a specified memory space.
VISA Memory Allocation	Returns an offset into a device's region that has been allocated for use by the session.
VISA Memory Free	Frees the memory previously allocated by the VISA Memory Allocation function.
VISA Move In 8/16/32	Moves a block of data from device memory to local memory in accesses of 8-bits, 16-bits, or 32-bits, respectively.
VISA Move Out 8/16/32	Moves a block of data from local memory to device memory in accesses of 8-bits, 16-bits, or 32-bits, respectively.
VISA Open	Opens a session to the specified device and returns a session identifier.
VISA Out 8/16/32	Writes 8-bits, 16-bits, or 32-bits of data, respectively, to the specified memory space.
VISA Peek 8/16/32	Reads an 8-bit, 16-bit, or 32-bit value, respectively, from the mapped memory region.
VISA Poke 8/16/32	Writes an 8-bit, 16-bit, or 32-bit value, respectively, to the mapped memory region.
VISA Read	Reads data from a device.
VISA Read STB	Reads a service request status from a message-based device.
VISA Status Description	Retrieves a user-readable string that describes the status code.
VISA Unlock	Relinquishes the lock previously obtained using the VISA Lock function.

VISA Function	Description
VISA Unmap Address	Unmaps memory space previously mapped by VISA Map Address.
VISA Wait On Event	Suspends execution of a thread of application and waits for an event of the specified type.
VISA Write	Writes data to the device.

VISA Attribute Node

The VISA Attribute Node gets and/or sets the indicated attributes. The following table lists the VISA attributes:

VISA Attribute	Description
FDC Mode	Specifies which fast data channel (FDC) transfer mode to use (VXI only).
FDC Number	Determines which FDC will be used to transfer the buffer (VXI only).
FDC Pairs	Specifies use of a channel pair for transferring data.
FDC Signal Enable	Lets the servant send a signal when control of the FDC channel is passed back to the commander.
GPIB Primary Address	Specifies the primary address of the GPIB device used by the given session.
GPIB Secondary Address	Specifies the secondary address of the GPIB device used by the given session.
IO Protocol	Specifies which protocol to use.
Immediate Servant	Determines if the VXI device is an immediate servant of the local controller.
Increment Destination Count	Specifies the number of elements by which to increment the destination address on block move operations.

VISA Attribute	Description
Increment Source Count	Specifies the number of elements by which to increment the source address on block move operations.
Interface Number	Specifies the board number for the given interface.
Interface Type	Specifies the interface type of the given session.
Mainframe Logical Address	Specifies the lowest logical address in the mainframe.
Manufacturer ID	The manufacturer identification number of the VXIbus device.
Maximum Queue Length	Specifies the maximum number of events that can be queued at any time on the given session.
Model Code	Specifies the model code for the VXIbus device.
Resource Lock State	Reflects the current locking state of the resource that is associated with the given session.
Resource Manufacturer ID	A value corresponding to the VXI manufacturer ID of the manufacturer that created the implementation.
Resource Manufacturer Name	A string that corresponds to the VXI manufacturer name of the manufacturer that created the implementation.
Resource Name	Unique identifier for a resource.
Send End Enable	Specifies whether to assert END during the transfer of the last byte of the buffer.
Slot	Specifies the physical slot location of the VXIbus device.
Suppress End Enable	Specifies whether to suppress the END bit termination.
Termination Character	The termination character.
Termination Character Enable	Determines whether the read operation should terminate when a termination character is received.
Timeout Value	Specifies the timeout value to use when accessing the device associated with the given session.

VISA Attribute	Description
Trigger Identifier	Identifier for the current triggering mechanism.
User Data	Used privately by the application for a particular session.
Version of Implementation	Uniquely identifies the current revision or implementation of a resource.
Version of Specification	Uniquely identifies the version of the VISA specification to which the implementation is compliant.
VXI Commander LA	The logical address of the commander of the VXI device.
VXI Logical Address	Specifies the logical address of the VXI device used by the given session.
VXI Memory Address Space	Specifies the VXIbus address space used by the device.
VXI Memory Base Address	Specifies the base address of the device in VXIbus memory address space.
VXI Memory Size	Specifies the size of memory requested by the device in VXIbus address space.
Window Access	Specifies the modes in which the current window may be accessed.
Window Base Address	Specifies the base address of the interface bus to which this window is mapped.
Window Size	Specifies the size of the region mapped to this window.

GPIB Function Overview

Introduction to the GPIB

The General Purpose Interface Bus (GPIB) is a link, or interface system, through which interconnected electronic devices communicate.

History of the GPIB

Hewlett-Packard designed the GPIB (originally called the HP-IB) to interconnect and control its line of programmable instruments. The GPIB was soon applied to other applications such as intercomputer communication and peripheral control because of its 1 Mbytes/s maximum data transfer rates. It was later accepted as IEEE Standard 488-1975 and has since evolved into ANSI/IEEE Standard 488.1-1987. The versatility of the system prompted the name General Purpose Interface Bus. For a basic description of the GPIB, see Appendix B, *Operation of the GPIB*.

National Instruments brought the GPIB to users of non-Hewlett-Packard computers and devices, specializing in both high-performance, high-speed hardware interfaces and comprehensive, full-function software. The GPIB functions for LabVIEW follow the IEEE 488.1 specification.

The IEEE 488.2 Standard

The ANSI/IEEE Standard 488.2-1987 expanded on the earlier IEEE 488.1 standard to describe exactly how the Controller should manage the GPIB, including the standard messages that compliant devices should understand, the mechanisms for reporting device errors and other status information, and the various protocols that discover and configure compliant devices connected to the bus.

The original standard, renamed IEEE 488.1, addressed only the hardware specifications of the GPIB cable and basic protocols. Its main shortcoming was that it left the interpretation of the standard as it applied to GPIB devices up to the instrument manufacturers. Thus, each GPIB instrument had a unique command set. To integrate each instrument into a particular GPIB system, programmers had to learn programming particulars for each device, a time-consuming and frustrating process. IEEE 488.2 specifically states how compliant

devices must communicate. This standard, along with Standard Commands for Programmable Instruments (SCPI), which defines specific function-dependent command sets, makes instrument programming more uniform.

The IEEE 488.2 standard also addresses Controller issues, such as the capabilities a compatible Controller must have. For example, the ability to monitor any of the bus lines at any time is crucial for detecting active devices (Talkers and Listeners) on the GPIB. IEEE 488.2 also defines the bus commands and protocols a Controller must use. The new standard also lists minimum functionality requirements, which directly influence the style of the NI-488.2 software in general and the GPIB 488.2 functions for LabVIEW in particular. See Appendix B, *Operation of the GPIB*, for more information on Talkers, Listeners, and Controllers.

Compatible GPIB Hardware

The following National Instruments GPIB hardware products are compatible with LabVIEW:

LabVIEW for Windows 3.x

AT-GPIB/TNT, AT-GPIB	IEEE 488.2 interface board for the IBM PC AT and compatible computers that have 16-bit plug-in slots
PCMCIA-GPIB	IEEE 488.2 interface card for computers with Type II PCMCIA slots
MC-GPIB	IEEE 488.2 interface board for the IBM PS/2 and compatible computers that have MicroChannel plug-in slots
GPIB-PCII/IIA	IEEE 488.2 interface board for the IBM PC/XT/AT
GD-GPIB	IEEE 488.2 interface board for the GRiDCASE 1500 series computer
GPIB-232CT-A	External RS-232/IEEE 488.2 Controller
GPIB-1284CT	External parallel/IEEE 488.2 Controller

LabVIEW for Windows NT

AT-GPIB/TNT, AT-GPIB	IEEE 488.2 interface board for the IBM PC AT and compatible computers that have 16-bit plug-in slots
----------------------	--

LabVIEW for HP-UX

AT-GPIB/TNT	Plug-in ISA IEEE 488.2 interface board for the HP series 700
EISA-GPIB	Plug-in EISA IEEE 488.2 interface board for the HP series 700
GPIB-ENET	External Ethernet/IEEE 488.2 Controller for the HP series 700

LabVIEW for Sun

SB-GPIB/TNT, SB-GPIB	Plug-in IEEE 488.2 Controller for Sun SPARCstation series computers
GPIB-SCSI-A	External IEEE 488.2/SCSI Controller for Sun SPARCstation series computers
GPIB-ENET	External Ethernet/IEEE 488.2 Controller for Sun SPARCstation series computers

LabVIEW for Macintosh

NB-GPIB/TNT, NB-GPIB	IEEE 488.2 interface board for Macintosh computers with full-size NuBus plug-in slots
NB-GPIB-P/TNT, NB-GPIB-P	IEEE 488.2 interface board for Macintosh computers with half-size NuBus plug-in slots
LC-GPIB	IEEE 488.2 interface board for the Macintosh LC series computers

NB-DMA2800	IEEE 488 DMA interface board for Macintosh NuBus computers
NB-DMA-8-G	IEEE 488 DMA interface board for Macintosh NuBus computers
GPIB-SCSI-A	External SCSI-to-IEEE 488.2 Controller
GPIB-SCSI	External SCSI-to-IEEE 488 Controller
GPIB-422CT	External RS-422/IEEE 488 Controller

Traditional GPIB Functions

The traditional GPIB functions are compatible with the GPIB boards listed in the *Compatible GPIB Hardware* section of this chapter.

These traditional GPIB functions are compatible with existing LabVIEW programs. You can also use traditional GPIB functions with GPIB 488.2 functions to handle situations outside the scope of the IEEE 488.2 standard, such as communicating with devices that do not comply with the IEEE 488.2 standard and managing the bus in unusual ways.

The following table lists the traditional GPIB functions:

Traditional GPIB Function	Description
GPIB Clear	Clears a single device.
GPIB Initialization	Configures characteristics of the GPIB interface.
GPIB Misc	Performs specialized GPIB operations.
GPIB Read	Reads data bytes from a GPIB device.
GPIB Serial Poll	Serial polls a single device.
GPIB Status	Indicates the status of the GPIB Controller.
GPIB Trigger	Triggers a single device.
GPIB Wait	Waits for a GPIB event to occur.

Traditional GPIB Function	Description
Wait for GPIB RQS	Waits for the assertion of the SRQ line Service Request.
GPIB Write	Sends data bytes to a single GPIB device.

LabVIEW GPIB 488.2 Functions

Using LabVIEW GPIB 488.2 functions together with IEEE 488.2-compatible devices improves the predictability of instrument and software behavior and lessens programming differences between instruments of different manufacturers.

The latest revisions of many National Instruments GPIB boards are fully compatible with the IEEE 488.2 specification for Controllers. The LabVIEW package also contains functions that make use of IEEE 488.2. By using these functions, your programming interface will strictly adhere to the IEEE 488.2 standard for command and data sequences.

The GPIB 488.2 functions contain the same basic functionality as the traditional GPIB functions, and include the following enhancements and additions:

- You specify the GPIB device address with an integer instead of a string. Further, you specify the bus number with an additional numeric control, which makes dealing with multiple GPIB interfaces easier.
- You can determine the GPIB status, error, and/or byte count immediately from the connector pane of each GPIB 488.2 function. You no longer need to use the GPIB Status Function to obtain error and other information.
- The FindLstn Function implements the IEEE 488.2 Find All Listeners protocol. You can use this function at the beginning of an application to determine which devices are present on the bus without knowing their addresses.
- The GPIB Misc Function is still available, but it is no longer necessary in most cases. IEEE 488.2 specifies routines for most GPIB application needs, which are implemented as functions. However, you can mix the GPIB Misc Function, as well as other GPIB functions, with the GPIB 488.2 functions if you need to.

- There are GPIB 488.2 functions with low-level as well as high-level functionality, to suit any GPIB application. You can use the low-level functions in Non-Controller situations or when you need additional flexibility.
- Although you must use an IEEE 488.2-compatible Controller to use these functions, they can control both IEEE 488.1 and IEEE 488.2 devices. The GPIB 488.2 functions are divided into five functional categories: single-device, multiple-device, bus management, low-level, and general.

Single-Device Functions

The single-device functions perform GPIB I/O and control operations with a single GPIB device. In general, each function accepts a single-device address as one of its inputs. The following table lists the single-device functions:

Single-Device Function	Description
DevClear	Clears a single device.
PPollConfig	Configures a device for parallel polls.
PassControl	Passes control to another device with Controller capability.
ReadStatusByte	Serial polls a single device to get its status byte.
Receive	Reads data bytes from a GPIB device.
Send	Sends data bytes to a single GPIB device.
Trigger	Triggers a single device.

Multiple-Device Functions

The multiple-device functions perform GPIB I/O and control operations with several GPIB devices at once. In general, each function accepts an array of addresses as one of its inputs.

The following table lists the multiple-device functions:

Multiple-Device Function	Description
AllSpoll	Serial polls all devices.
DevClearList	Clears multiple devices simultaneously.
EnableLocal	Enables operations from the devices' physical front panels.
EnableRemote	Enables remote GPIB programming of devices.
FindRQS	Determines which device is requesting service.
PPoll	Performs a parallel poll.
PPollUnconfig	Unconfigures a device for parallel polls.
SendList	Sends data bytes to multiple GPIB devices.
TriggerList	Triggers multiple devices simultaneously.

Bus Management Functions

The bus management functions perform system-wide functions or report system-wide status. The following table lists the bus management functions:

Bus Management Function	Description
FindLstn	Finds all Listeners on the GPIB.
ResetSys	Performs an IEEE 488.2 three-step system initialization.
SendIFC	Clears the GPIB functions with Interface Clear.
SendLLO	Sends the Local Lockout message to all devices.

Bus Management Function	Description
SetRWLS	Places devices in the Remote Mode With Lockout State.
TestSRQ	Determines the state of the SRQ line.
TestSys	Instructs IEEE 488.2 devices to conduct self-tests.
WaitSRQ	Suspends operation until SRQ is asserted.

Low-Level Functions

The low-level functions let you create a more specific, detailed program than higher-level functions. You use low-level functions for unusual situations or for situations requiring additional flexibility.

The following table lists the low-level functions:

Low-Level Function	Description
SendDataBytes	Sends data bytes to previously addressed devices.
RcvRespMsg	Reads data bytes from a previously addressed device.
ReceiveSetup	Prepares a device to send data bytes and prepares the GPIB interface to read them.
SendCmds	Sends GPIB command bytes.
SendSetup	Prepares particular devices to receive data bytes.

General Functions

The general functions are useful for special situations. The following table lists the general functions:


General Function	Description
MakeAddr	Formats a primary and secondary address for use with the GPIB 488.2 Functions.
SetTimeOut	Changes the global timeout period for all GPIB 488.2 Functions and sets the default timeout period for all GPIB functions.

Serial Port VI Overview

The serial port VIs configure the serial port of your computer and conduct I/O using that port. The following table lists the serial port VIs:

Serial Port VI	Description
Bytes at Serial Port	Reports the number of bytes available in the serial port input buffer.
Serial Port Break	Sends a break on the specified output port.
Serial Port Init	Initializes the selected serial port by specifying the protocol, handshaking, and buffer size.
Serial Port Read	Reads characters from the specified serial port.
Serial Port Write	Writes a data string to the specified serial port.

Overview of LabVIEW Instrument Drivers

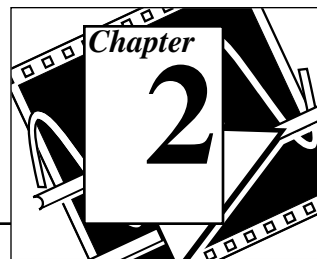

This chapter contains an overview of the LabVIEW instrument driver model, VIs, and VI templates.


Introduction

A LabVIEW instrument driver is a set of VIs that control a programmable instrument which each correspond to a programmatic operation. LabVIEW instrument drivers simplify instrument control and reduce test program development time by eliminating the need to learn the low-level programming protocol for each instrument.

The following two conceptual models help define a standard for LabVIEW instrument driver software design, development and use. The first model, the instrument driver external interface model, shows how the instrument driver interfaces with other system components. The second model, the instrument driver internal design model, defines the internal organization of an instrument driver software module.

LabVIEW Instrument Driver External Interface Model

The following figure shows a general model of how a LabVIEW instrument driver interfaces to the rest of the system:

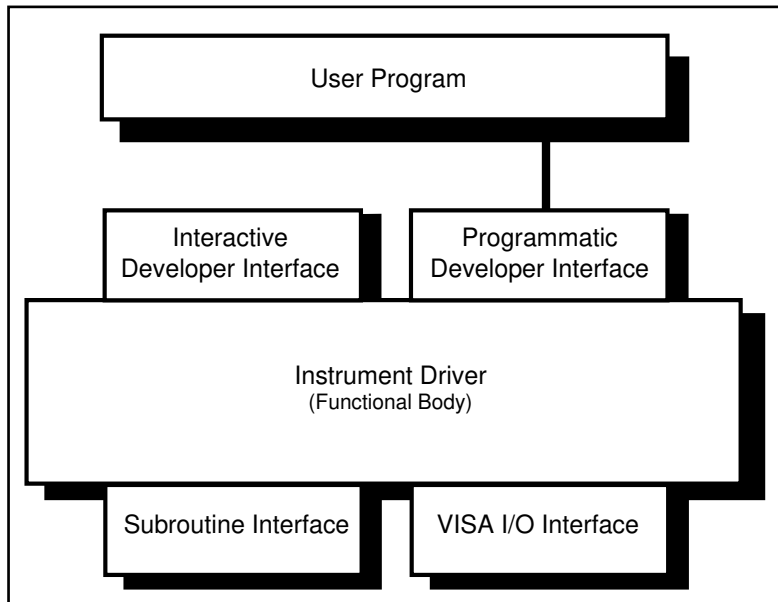

Functional Body

The *functional body* is the actual code for the instrument driver. Refer to the *LabVIEW Instrument Driver Internal Design Model* section of this chapter, for more information.

The most successful instrument driver products have, historically, been developed by using a standard programming language for the functional body. This is the approach LabVIEW instrument drivers use. The advantages include greater developer control over the driver, more robust drivers, and increased functionality. LabVIEW instrument drivers are written using the standard LabVIEW graphical programming environment.

The functional body of a LabVIEW instrument driver is a set of VIs that control a specific instrument. The source code for these VIs are block diagrams consisting of executable icons connected by data flow wires. Because the functional body is developed with the standard tools

provided in LabVIEW, users can easily view instrument driver source code and optimize it for their application.

Interactive Developer Interface

The interactive developer interface of a LabVIEW instrument driver is the front panel. It is analogous to a physical instrument panel and is the interactive user interface of the VI. On the panel, controls and indicators graphically represent the inputs and outputs of the VI. With the LabVIEW front panel, users can interactively operate individual instrument driver VIs and verify communication.

Programmatic Developer Interface

The icon/connector is the programmatic interface of the LabVIEW instrument driver VI. It consists of a graphical representation of the VI (icon) and a definition of the input and output terminals for the VI (connector). When you call or execute a VI from another VI, you place a copy of the subVI icon/connector in the block diagram of the calling VI. Information passes between the two VIs through the connector terminals. There are several benefits to this approach. You can easily assemble test systems utilizing LabVIEW instrument drivers by combining a few instrument driver VIs, each using multiple parameters. The instrument driver interface in the user program is modular and easy to identify, and you can recall the VI front panels during debugging to understand how the program uses the instrument driver.

I/O Interface

An important consideration for instrument drivers is how they perform their I/O to and from instruments. The I/O interfaces for LabVIEW instrument drivers are the VISA and GPIB function libraries, and the VXI and Serial VI libraries. These libraries contain sets of functions and VIs that cover the capabilities of GPIB, VXIbus, and Serial bus capabilities, including both message-based and register-based programming, interrupt and event handling, and direct access to the VXI backplane.

VISA, an acronym for Virtual Interface Software Architecture, is a single interface library for controlling VXI, GPIB, RS-232, Ethernet, and other types of instruments. Refer to Chapter 5, *VISA Library Reference*, for further information.


Subroutine Interface

Because you write LabVIEW instrument drivers in standard LabVIEW graphical code, an instrument driver is a software program with the same capabilities as any other LabVIEW VI. While some VIs (such as instrument drivers) perform only simple I/O to and from an instrument, other VIs may control multiple instruments or use support libraries to integrate data analysis or other measurement-specific operations. With LabVIEW, you can build virtual instruments that combine hardware and software capabilities. You can develop and package complete, high-level tests as single VIs, which other test developers can reuse.

By ensuring compatibility with the virtual instrument concept, the LabVIEW instrument driver standard has unlimited potential for delivering baseline as well as sophisticated application-specific instrument drivers. The LabVIEW instrument driver standard defined in this document applies both to instrument drivers that control only a single instrument, and to virtual instrument drivers that combine features of multiple instruments and add software processing.

LabVIEW Instrument Driver Internal Design Model

The LabVIEW instrument driver internal design model, shown in the following figure, defines the organization of the LabVIEW instrument driver *functional body*. Because development guidelines and all LabVIEW instrument drivers are based on this model, it is important to both developers and end users of instrument drivers. When you understand the model and how to use one instrument driver, you can use that knowledge across numerous instrument drivers.

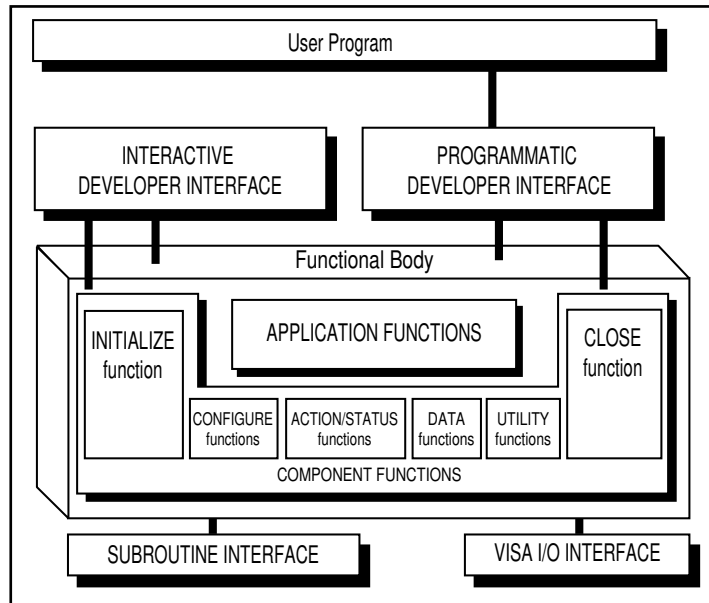


The functional body of a LabVIEW instrument driver consists of two main categories of VIs. The first category is a collection of *component VIs*, which are individual software modules that each control a specific type of instrument function. The second category is a collection of higher-level *application VIs* that show how to combine the component VIs to perform basic test and measurement operations with the instrument.

The internal design model of LabVIEW instrument drivers is built on a proven methodology. With this model, you have the necessary granularity to control instruments properly in your software applications. You can, for example, initialize all instruments once at the start, configure multiple instruments, and then trigger several instruments simultaneously. As another example, you can initialize and configure an instrument once, and then trigger and read from the instrument several times.

Instrument Driver Application VIs

The *application VIs* are at the highest level of the instrument driver hierarchy. They are written in LabVIEW block diagram source code and control the most commonly used instrument configurations and measurements. These VIs serve as a code example for how to configure the instrument for a common operation, trigger the instrument, and take measurements. Because the application VIs are standard VIs, with icons and connector panes, you can call them from any higher-level application when you want a single, measurement-oriented, interface to the driver. For many developers, the application VIs are the only instrument driver VIs needed for instrument control. The Tek VX4790 Example VI, shown in the following figure, demonstrates an application VI front panel.

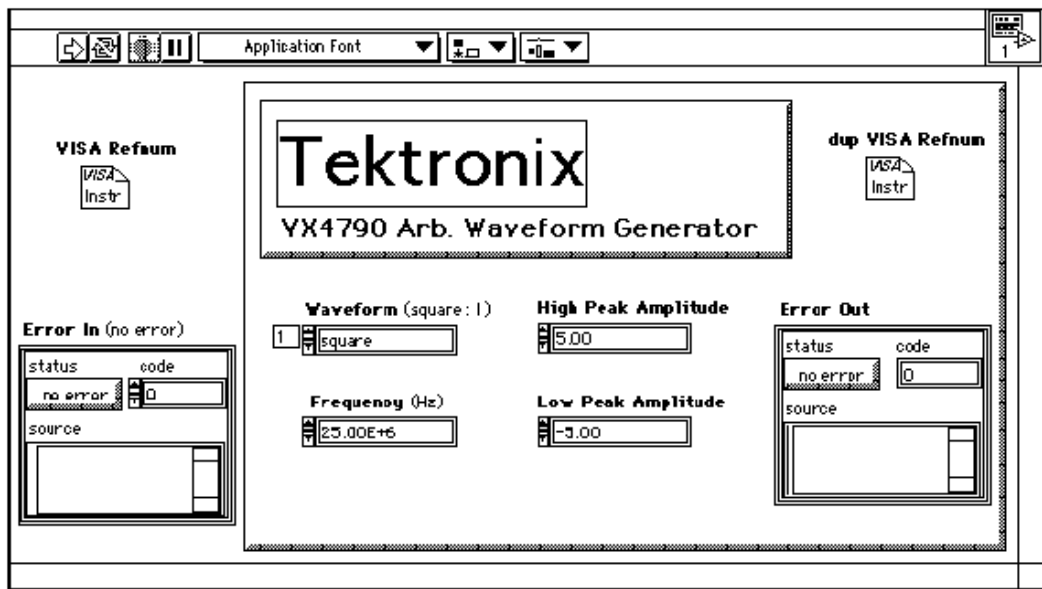

The application VIs are built from a lower-level set of instrument driver *component VIs*.

Instrument Driver Component VIs

LabVIEW instrument drivers have *component VIs*, which are a modular set of VIs that contain all of the instrument configuration and measurement capabilities. The component VIs fit into six categories: initialize, configuration, action/status, data, utility, and close.

All LabVIEW instrument drivers should have an *initialize VI*. It is the first instrument driver VI called, and establishes communication with the instrument. Additionally, it can perform any necessary actions to place the instrument either in its default power on state or in some other specific state.

The *configuration VIs* are a collection of software routines that configure the instrument to perform the desired operation. There may be numerous configuration VIs, depending on the particular instrument. After these VIs are called, the instrument is ready to take measurements or stimulate a system.

The *action/status* category contains two types of VIs. *Action VIs* cause the instrument to initiate or terminate test and measurement operations. These operations can include arming the trigger system or generating a stimulus. These VIs are different from the configuration VIs because they do not change the instrument settings, but only order the instrument to carry out an action based on its current configuration.

The *status VIs* obtain the current status of the instrument or the status of pending operations. The specific routines in this category and the actual operations they perform are left up to you.

The *data VIs* transfer data to or from the instrument. Examples include VIs for reading a measured value or waveform from a measurement instrument, VIs for downloading waveforms or digital patterns to a source instrument, and so on. The specific routines in this category and the actual operations performed by those routines are left up to you.

The *utility VIs* can perform a variety of operations that are auxiliary to the most often used instrument driver VIs. These VIs include the majority of the *instrument driver template VIs* such as reset, self-test, revision, error query, and error message and may include other custom instrument driver VIs, such as calibration or storing and recalling setups.

All LabVIEW instrument drivers should include a *close VI*. The close VI terminates the software connection to the instrument and deallocates system resources.

Each of these categories, with the exception of the initialize and close VIs, consists of several modular VIs. Most of the critical work in developing an instrument driver lies in the initial design and organization of the instrument driver component VIs. The specific

routines in each category are further categorized as either *template VIs* or *developer-specified VIs*.

The *template VIs* are instrument driver VIs that you can use as templates or examples. These VIs perform common operations such as initialize, close, reset, self-test, and revision query. The template VIs contain modification instructions for their use in a specific instrument driver for a particular instrument. For more information, refer to Chapter 4, *Instrument Driver Template VIs*.

The remainder of instrument driver VIs are known as *developer-specified VIs*, and the actual operations performed by those routines are left up to you. Although all instruments will have configuration VIs, some instruments can have a different number of configuration VIs depending on the unique capabilities of the instrument.

Figure 2-1 shows how the Tek VX4790 Example application VI diagram uses the instrument driver component VIs:

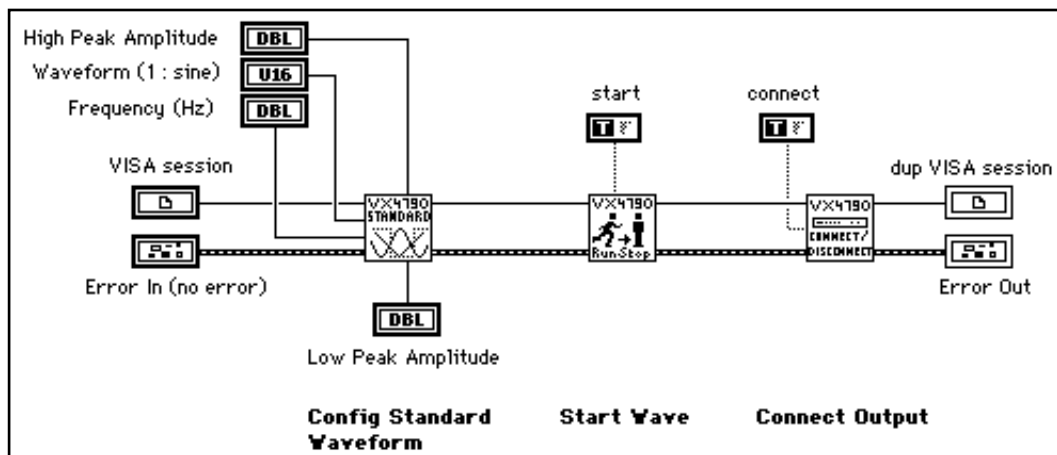


Figure 2-1. VIs in Tek VX4790 Example Diagram

The block diagram of the instrument driver component VIs uses standard LabVIEW VIs, as well as VISA VIs to build command strings and send them to the instrument. In the following figure, the Tek VX4790 Config Std Wave component VI block diagram assembles the command string and wires it into the VISA Write function. This function performs the necessary I/O, checks for errors, and updates the appropriate error indicators.

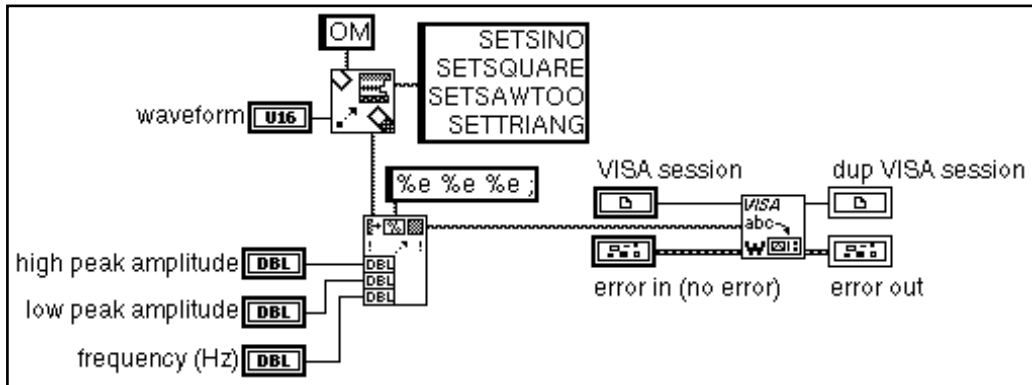



Figure 2-2. Tek VX4790 Config Std Wave Diagram

Error Reporting

LabVIEW instrument drivers use error clusters to report all errors, as shown in the following figure. Inside the cluster, a Boolean error indicator, a numeric error code, and an error source string indicator report if there is an error, the specific error condition, and the source (name) of the VI in which the error occurred. Additional comments may also be included. Each instrument driver VI has an error in and an error out terminal defined on its connector pane in the lower left and lower right terminals respectively. By wiring the error out cluster of one VI to the error in cluster of another VI, you can pass error information all the way through your instrument driver and out to your full application.

Another benefit of error input/output is that data dependency is added to VIs that are not otherwise data dependent.

Developing a LabVIEW Instrument Driver

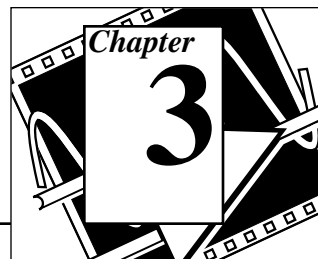

This chapter describes the procedure for developing a LabVIEW instrument driver. The ideal LabVIEW instrument driver has full function control of the instrument. Rather than mandate the required functionality of all instrument types, such as DMMs, counter/timers, and so on, this chapter focuses on the architectural guidelines of all drivers. With this information, driver developers can implement functionality unique to a particular instrument, and still organize, package and use all drivers in the same way.

Development Procedure

The best way to develop a LabVIEW Instrument Driver is to follow a three-step process. In step one, you design the instrument driver structure. In step two, you modify the instrument driver templates VIs. In step three, you add developer defined VIs.

Designing the Instrument Driver Structure

The ideal instrument driver does what the user needs—no more and no less. No particular type of driver design is perfect for everyone, but by carefully studying the instrument and grouping controls into modular VIs, you can satisfy most users.

When the number of programmable controls in an instrument increases, so does the need for modular instrument driver design since a single VI cannot access all features. However, when an instrument driver contains hundreds of VIs, each controlling a single instrument feature, more instrument rules regarding command order and interaction apply. Modular design simplifies the tasks of controlling the instrument and modifying VIs to meet special requirements.

Ideally, you should devise the overall structure of your instrument driver before you build the individual VIs. A useful instrument driver is more than a series of VIs; it is a tool to help users develop application programs. You should design an instrument driver with the application and end user in mind.

You must create some instrument driver VIs that control unique instrument features. However, you can use template VIs for common operations. For more information about template VIs see Chapter 4, *Instrument Driver Template VIs*.

Instrument Driver Structure and VI Hierarchy

When you develop a LabVIEW instrument driver, it is important to clearly define the structure and VI hierarchy of the driver. First, define the primary VIs and develop a modular VI hierarchy. This hierarchy is the design document for a LabVIEW instrument driver.

Useful instrument drivers come from an in-depth knowledge of the instrument operation and use in test applications. The following steps outline one approach to developing the structure for the LabVIEW instrument drivers:

1. Familiarize yourself with the instrument operation. Read the operating manual thoroughly. Typically the foundation of the driver hierarchy is in the instrument programming manual. Learn how to use the instrument interactively before you attempt any programming.
2. Use the instrument in an actual test set-up to get practical experience. (The operating manual may explain how to set up a simple test.)
3. Study the programming section of the manual. Skim the instruction set to see which controls and functions are available and how the features are organized. Decide which features are best suited for programmatic use.
4. Examine instrument drivers for similar instruments. Often instruments from the same family have the same programming command set and you can easily modify their corresponding instrument drivers.
5. Determine which LabVIEW template VIs are suitable for use with your instrument.
6. Develop a structure for the driver by looking for controls that are used together to perform a single task or function. The sections of a well organized manual often correspond to the functional groupings of an instrument driver.

Instrument Driver VI Organization

After you have developed your Instrument Driver structure, you can develop a VI hierarchy to organize the VIs that will be necessary to create the driver.

The VI organization of an instrument driver defines the hierarchy and overall relationship of the instrument driver component VIs.

You define the majority of instrument driver VIs and design them to access the unique capabilities of a particular instrument. However, many operations common to all types of instrumentation are performed by the template instrument driver VIs: initialize, close, reset, self-test, revision, error query, and error message.

The template VIs for LabVIEW instrument drivers include prewritten VIs to perform these common instrument operations. The command strings are based on IEEE 488.2 Common Commands. To include these VIs in your instrument driver, modify the command strings as required for your instrument. If the instrument is IEEE 488.2 compliant, little or no modifications are needed. If you are developing a driver for a non-IEEE 488.2 compliant or a register-based device, you will develop equivalent VIs for your instrument.

A class is a group of VIs that perform similar operations. Common classes of VIs are configuration, action, data, and utility.

The following table shows an example instrument driver organization for an oscilloscope. At the highest level of the hierarchy, you see the template VIs, initialize and close and the typical classes of VIs.

VI Hierarchy	Type
Initialize VI	(Template)
Application VIs <ul style="list-style-type: none"> • Autosetup and Read Waveform • Rise-Time/Fall-Time Measurement 	(Developer Defined) (Developer Defined)

VI Hierarchy	Type
Configuration VIs <ul style="list-style-type: none"> • Configure Vertical • Configure Horizontal • Configure Trigger • Configure Acquisition Mode • Autosetup 	(Developer Defined) (Developer Defined) (Developer Defined) (Developer Defined) (Developer Defined)
Action VIs <ul style="list-style-type: none"> • Acquire Data 	(Developer Defined)
Data VIs <ul style="list-style-type: none"> • Read Waveform • Voltmeter Measurement • Counter/Timer Measurement 	(Developer Defined) (Developer Defined) (Developer Defined)
Utilities VIs <ul style="list-style-type: none"> • Reset • Self-Test • Revision • Error Query • Error Message 	(Template) (Template) (Template) (Template) (Template)
Close VI	(Template)

Guidelines and Recommendations

- Design an instrument driver VI front panel that contains all the controls required to perform the VI task.

For example, a configure measurement VI would contain only the necessary controls to configure the instrument to take the measurement. It would not take the measurement or configure any other features. Additional VIs included in the instrument driver perform these tasks.

- Design a modular instrument driver that contains a set of VIs, each performing a logical task or function such as configuring the instrument or taking a measurement.

A modular instrument driver is flexible and easy to use. For example, consider a digital multimeter driver design that uses a single VI to both configure the instrument and read a measurement. The user cannot read multiple measurements without reconfiguring the meter each time the VI executes. A

better approach is to build two VIs: one to configure the instrument, and one to read a measurement. Then the user can configure the meter once and take multiple measurements.

- Concentrate on the correct level of granularity of driver VIs and how these VIs will be used in a system.

An instrument driver with a few very high-level VIs may not give the user enough control of the instrument operation. Conversely, an instrument driver with many low-level VIs is difficult for users unfamiliar with instrument rules regarding command order and interaction. For example, when using a measurement device such as an oscilloscope, the user typically configures the instrument once and takes many measurements. In this case, you should write high-level configuration VIs for the device. On the other hand, when using a stimulus device such as a pulse generator, the user may want to vary individual parameters of the pulse to test the boundary conditions of his system, or perform frequency response tests. In this case, you should write lower-level VIs, so that users can access individual instrument capabilities instead of reconfiguring each time they want to change one component of the output.

- Consider the relationship of the driver with other instrument drivers in the system.

Typically, test designers want to initialize all of the instruments in a system at once, then configure them, take measurements, and finally close them at the end of the test. Good driver design includes logical division of operations.

- Create an instrument driver design (both in appearance and functional structure) that is similar to other instruments of the same type.

Instrument drivers across a family of similar instruments should be consistent in appearance, structure, and style. For example, all oscilloscope drivers should resemble each other, as should all multimeters, scanners, and sources. If possible, modify a copy of an existing driver of a similar instrument.

- Design an instrument driver that optimizes the programming capability of the instrument.

You can sometimes exclude documented functions that are not well-suited for programmatic use.

- Design each VI to be independent of other VIs.

If two or more VIs must always be used together, consolidate them into one VI.

- Minimize redundant parameters.

For example, the parameters for each channel of a multi-channel oscilloscope are similar or identical. Rather than duplicate the programming controls for each channel, you can include a VI control for selecting which channel to configure. The user can use this VI to change the settings for an individual channel, rather than configuring every channel each time the VI is called.

Design Example

Deciding which parameters to include in an instrument driver VI is one of the greatest challenges facing the instrument driver developer. Fortunately, organizational information is often available in the instrument's manuals. In particular, the programming section of the manual may group the commands into sections such as configuring a measurement, triggering, reading measurements, and so on. These groupings can serve as a model for a driver hierarchy. Begin to develop a structure for the driver by looking for controls that are used together to perform a single task or function. A modular driver will contain individual VIs for each of the control groups.

The following table shows how the command summary from the *Tektronix VX4790 Arbitrary Waveform Generator Operating Manual* relates to developer specified instrument driver VIs.


Instrument Manual Section	Instrument Driver VI
Setup Commands <ul style="list-style-type: none"> • External clock input enable • External trigger source • Sync pulse control • Isolation relay control 	TKVX4790 Setup
Pre-Programmed Waveform Commands <ul style="list-style-type: none"> • Sine wave • Square wave • Triangle wave • Sawtooth wave 	TKVX4790 Config Std. Waveform

Instrument Manual Section	Instrument Driver VI
Frequency Commands <ul style="list-style-type: none"> • Frequency • Period • Divide • Low-Pass filters 	TKVX4790 Config Sample Frequency
Voltage/Attenuator Commands <ul style="list-style-type: none"> • Voltage control • Attenuator enable • Attenuation level 	TKVX4790 Config Volt/Atten.
Arbitrary Waveform Commands <ul style="list-style-type: none"> • Sample voltage • Breakpoint/Last commands 	TKVX4790 Download Arb. Waveform
Trigger Commands <ul style="list-style-type: none"> • Start location • Breakpoint/last commands 	TKVX4790 Run/Stop

While the instrument manual can provide a great deal of information about how to structure the instrument driver, you should not rely on it exclusively. Your knowledge of the instrument and how it is used should be the ultimate guide. The preceding table shows manual sections that map nicely to VIs found in the instrument driver. There are instances when it is more appropriate to place commands from several different command groups in your VI.

Conversely, it is often necessary to take one group of commands and divide it into two or more VIs. Consider how an instrument manual groups the trigger configuration commands with the commands that actually perform the trigger arming and execution. In this case, you should separate the commands into two VIs; one to configure the trigger, and one that arms or triggers the instrument.

The following figure shows the LabVIEW instrument driver VIs for the Tektronix VX4790 Arbitrary Function Generator:

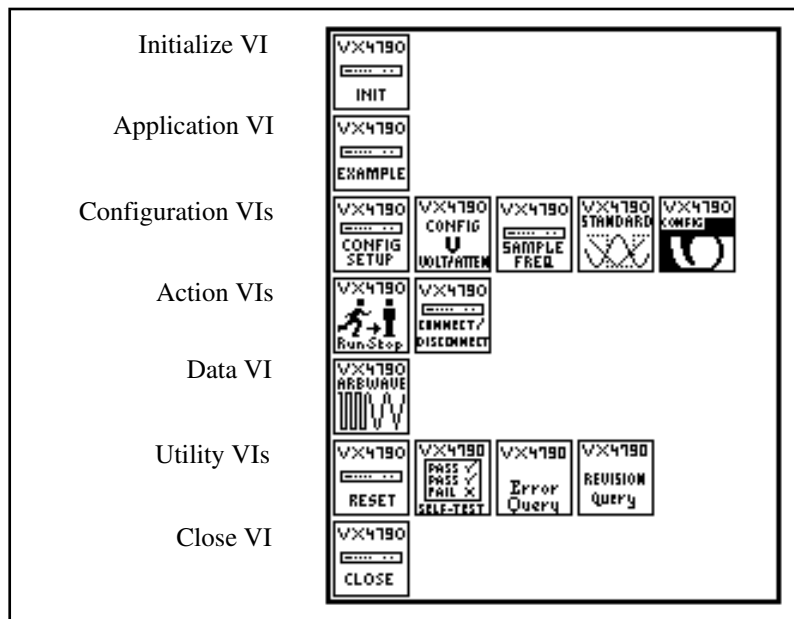

Modifying the Instrument Driver Templates

After you design the LabVIEW instrument driver structure, the next step is to modify the template VIs to represent your instrument. Most of the modifications involve the instrument prefix. The prefix is a unique identifier for the instrument driver, and is used as the filename for all files associated with the driver and as the prefix to all instrument VI names. Typically, the prefix is the combination of an abbreviation for the instrument vendor name and the model number. For example, the instrument prefix for the Tektronix VX4790 instrument driver is `tkvx4790`. As a default, the template instrument drivers use `PREFIX` as the instrument prefix.

Use the following procedure for modifying the LabVIEW instrument driver template:

1. Open the `PREFIX` Initialize VI in the file `insttmpl.llb`.
2. Save the VI into a new VI library file by using the prefix for your instrument as the filename of the `.llb` file. Save the VI replacing `PREFIX` in the VI name with the prefix for your instrument.

3. Follow the instructions in the Instructions to Modify This VI string control on the initialize panel to modify the VI for your particular instrument.
4. Edit all **Show VI Info...** and control and indicator descriptions.
5. Edit the icon. Create an icon for each of the color modes of the icon: Black and White, 16-Color, and 256-Color.
6. Delete the Instructions to Modify This VI string control after you have completed the modifications.
7. Resize the front panel and save the VI.
8. Repeat steps 1 through 7 for PREFIX Close VI and the remaining template VIs that your instrument uses. All LabVIEW instrument drivers should have initialize, close, reset, revision, and error message VIs.

After completing this procedure, you will have a base level driver that implements all template instrument driver VIs and is a good framework from which you can create your driver.

Adding Instrument Driver Component VI VIs

The final step in developing a LabVIEW instrument driver is to add the developer defined component VIs that define the functionality of the instrument driver and access the unique capabilities of your instrument. The VIs that you create will be added to the source code along with the template VIs in the file `prefix.llb`.

You can use the following procedure to add your new VIs:

1. Open either the PREFIX Message-Based or PREFIX Register-Based templates VI in the file `insttmpl.llb`. Use the PREFIX Message-Based template VI for message-based operations. Use the PREFIX Register-Based template VI for register-based operations.
2. Edit the VI front panel. Create the controls and indicators for the VI.
3. Edit all control and indicator Help information. Edit the **Show VI Info...** description.
4. Edit the icon. Create an icon for each of the color modes of the icon: Black and White, 16-Color, and 256-Color.
5. Edit the connector pane. Select an appropriate connector pattern and wire all controls and indicators to the terminals.

6. Edit the block diagram. Program all operations necessary to carry out the functionality of the instrument driver VI.
7. Save the VI.
8. Test the instrument driver VI.
9. Repeat these steps for every instrument driver component VI and application VI that you define for your instrument.
10. Edit the instrument driver .11b by selecting **File»Edit VI Library...** from the menu. Edit the **Functions** and **Controls** names. Edit the arrangement of icons in the Functions and Controls palettes.

Editing the block diagram source code is the most difficult step in adding a component VI to the instrument driver. Defining a block diagram structure makes it easier to edit the block diagram source code. You can divide this process into the following steps:

1. Place the appropriate I/O routines in the block diagram.
2. Wire the **error in** cluster terminal to the first I/O VI error input connector. Then wire the **error out** connector of that VI to the **error in** connector of the next VI. Continue this process for all of the I/O VIs. Then wire the **error out** connector of the last VI to the **error out** terminal of the icon.
3. Wire the **VISA session** to every I/O VI.
4. Use the LabVIEW string VIs to assemble a command string based on the VI inputs.
5. Wire the command string to the viWrite I/O VI.
6. Use the viRead I/O VI to read the response if an instrument response is generated.
7. Use the string VIs to parse the response and wire it to the appropriate indicator terminals.



Tips for Developing a LabVIEW Instrument Driver

Loop Termination Conditions

When you use looping structures in instrument driver block diagrams, you must include a way to escape from While Loops if an error occurs. This escape method is important if you are using a While Loop

containing I/O routines and the loop termination depends on the result of the I/O.

If there is an error, the I/O routines automatically shut down and LabVIEW may be stuck in an endless loop. Therefore, always test the error cluster status in conjunction with your normal loop termination condition to determine when to terminate the loop. The following two figures show incorrect and correct mechanisms for terminating a While Loop, respectively:

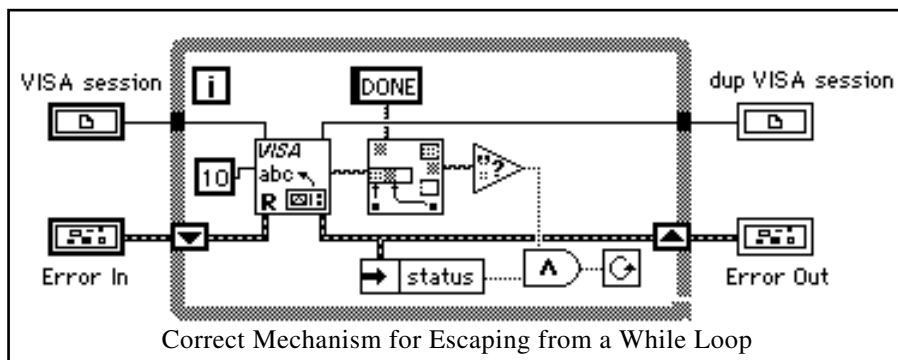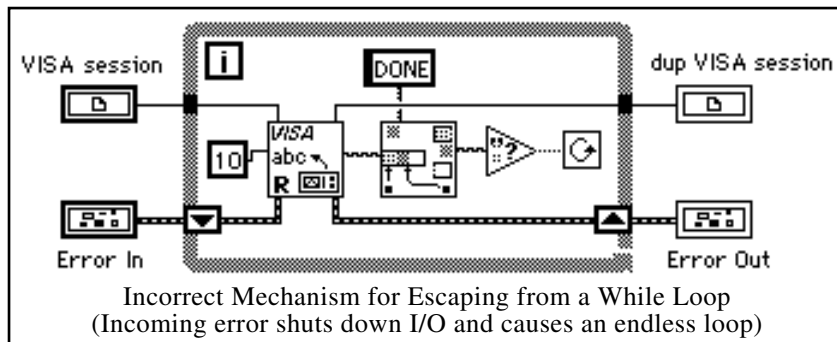


Assembling Command Strings

After you develop your front panel, the next step is to create the block diagram which performs the function required by the VI. Each type of front panel control has a corresponding block diagram string VI that simplifies the task of building command strings.

You can use Pick Line & Append to choose from a selection of strings and concatenate it to another string in a single step. This procedure is easier than using a Case structure and Concatenate Strings.

Use this block diagram...

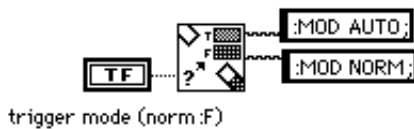


rather than this one...




You can use Format & Append to format and concatenate simple numeric values. This procedure is easier than using one of the To Decimal or To Exponential type conversion VIs in conjunction with Concatenate Strings.

Use this block diagram...

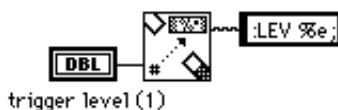


rather than this one...

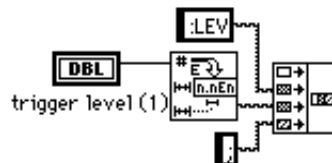


By using Select & Append you can select a string constant and concatenate it to another string in a single step. This procedure is easier than using Select and Concatenate Strings.

Use this block diagram...

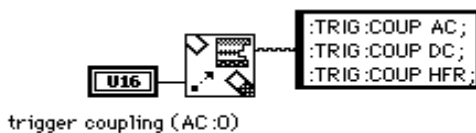

rather than this one...

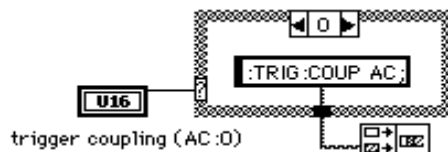

Data Dependency

Carefully consider the control flow when you build your diagrams. LabVIEW does not necessarily execute in a left-to-right, top-to-bottom fashion. Data dependency automatically determines execution order. Add artificial data dependency wherever appropriate (see the *LabVIEW User Manual* for more information). By using the clusters to chain I/O VIs together, you can define the execution order without using Case or Sequence structures, as illustrated in Figure 2-1, *VIs in Tek VX4790 Example Diagram*. Sequence structures, which hide parts of the diagram, are also effective at controlling execution order. Whichever method you use, make sure that you clearly define control flow so that the correct branch of the diagram executes first.

Programmatically Checking Parameter Ranges

Be prepared to handle unexpected values entered into controls. The front panel **Data Range...** item is a deterrent to many range errors. When circumstances prevent its use, you can either build a block diagram routine based on the Range Test VI shown in Figure 3-1 to report range errors to the error cluster, or you can build block diagram routines to coerce the data into range. Also, you must often account for control interaction by using case statements to alter the block diagram routines based on the setting of a particular control. Ideally, you can run the VI with any combination of settings and values without crashing the instrument. If not, add more error detection/correction to your block diagrams. Even though programmatic range checking can simplify the VI use, it can also complicate the development of the VI.

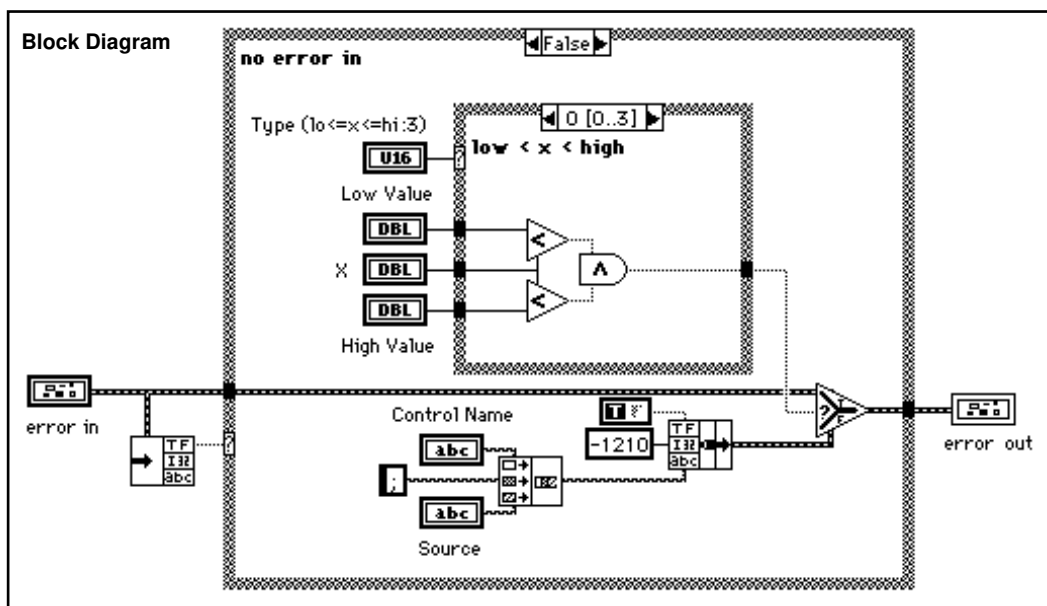
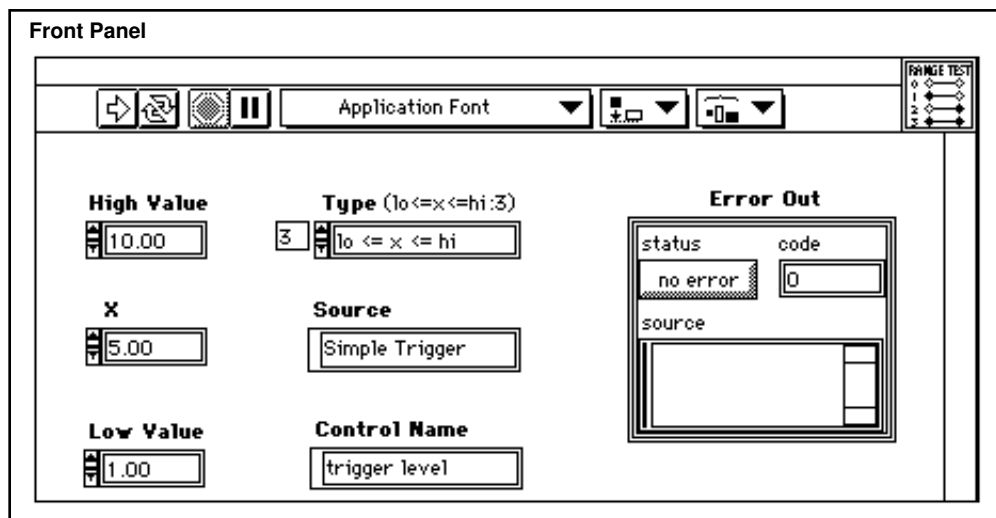


Figure 3-1. Range Test VI (Front Panel and Block Diagram)

Programmatic range checking can easily double the size of your VI and add some execution speed penalties. Figures 3-2 and 3-3 show the changes made to the Simple Trigger VI to programmatically check the ranges of the numeric inputs.

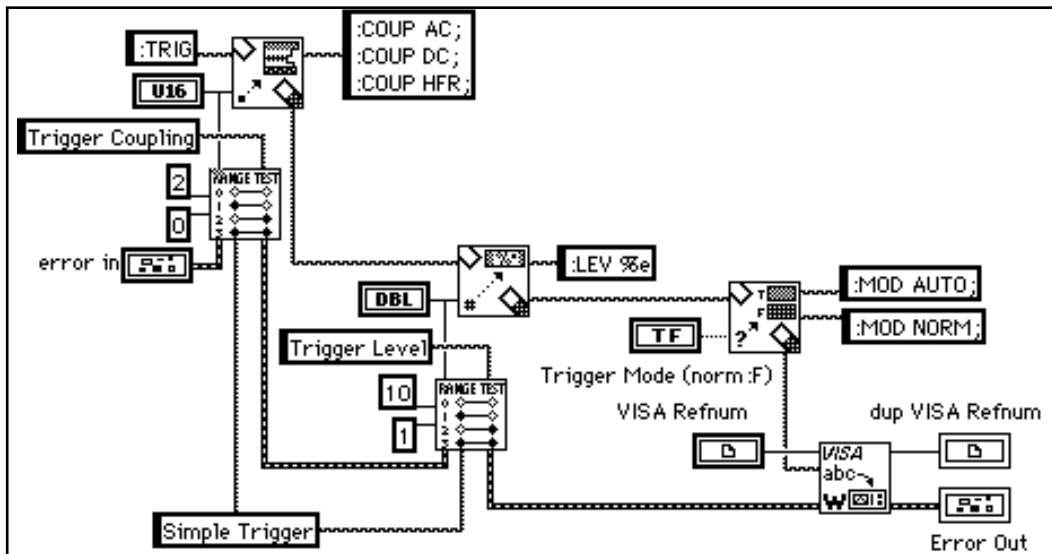


Figure 3-2. Simple Trigger VI with Programmatic Range Testing

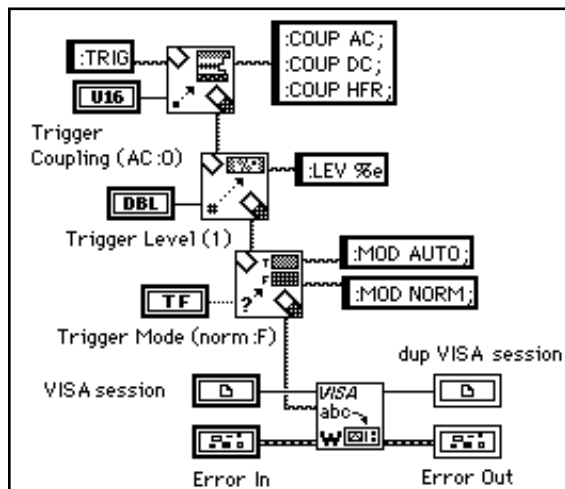


Figure 3-3. Simple Trigger VI without Programmatic Range Testing

Guidelines

Like the LabVIEW VI, the standard components of an instrument driver VI are the front panel, block diagram, and icon/connector pane. Special guidelines concerning these components, as well as error reporting and on-line help information, are described in the following sections.

Front Panel


Each VI in your instrument driver should contain a front panel that groups all the necessary controls together to perform the function of the VI. When you develop an instrument driver VI, decide which control styles best represent the instrument commands and options. Typically, you can categorize instrument commands into three types of control styles: Boolean, digital numeric, and text or ring numeric.

For example, you can represent any instrument command that has two options (such as TRIG:MODE:AUTO | NORMAL) on the front panel with a Boolean switch. In this case, label the switch **Trigger Mode** and add a free label showing the options: auto or normal. For commands that have a discrete number of options (such as TRIG:COUP:AC | DC | HFREJ), use a text ring rather than a digital numeric because the text ring labels each numeric value with the command it represents. Any command requiring a numeric parameter whose value varies over a wide range and cannot be better represented by a ring can be represented with a digital numeric.

The three control types of Boolean, numeric, and text ring represent most instrument commands on the front panel of your VIs. In addition, block diagram string icons specifically designed for use with these controls exist. These features can simplify string formatting and appending instrument commands into command messages, as discussed in the *Assembling Command Strings* and *Block Diagram* sections of this chapter.

Required Front Panel Controls

In addition to the controls required to operate the instrument, your front panel must also have the following controls:

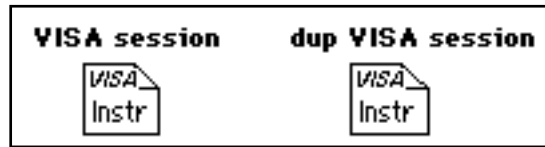

VISA session (except for the initialize VI) input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the Instrument Driver Template VI Error Reporting in chapter four for more information on errors in clusters.



error out is a cluster containing error information. If **error in** indicates an error, the **status**, **code**, and **source** elements of **error out** have the same values as the corresponding elements of **error in**. If **error in** does not indicate an error, **error out** describes the error encountered by the VI. Refer to Appendix D, *Error Codes*, for a description of the possible error codes. See the Instrument Driver Template VI Error Reporting in chapter four for more information on errors in clusters.

To gain consistency with other LabVIEW instrument drivers, place the **VISA session** control and **dup VISA session** indicator in the upper left and upper right corners of the front panel, and the **error in** and **error out** clusters in the lower left and lower right corners, respectively.



Note:

You can place the error in cluster outside the panel window because it has no interactive meaning and is only needed for programmatic use.

Control Guidelines

When placing controls on your front panels, use the following style guidelines to ensure uniformity with other LabVIEW VI front panels:

- Use the default font (application) for all LabVIEW instrument driver front panel control labels. The application font is available on all platforms that are using LabVIEW.
- Use **bold** text for control name labels that denote important or primary controls, and reserve plain text for secondary controls.



Note:

In most cases, all instrument driver controls are primary and require bold text. If you are finding yourself placing many secondary or auxiliary controls on panels, this may indicate the need to subdivide your VI into two or more VIs.

- Use lowercase letters in all words, except abbreviations or acronyms, which require caps (such as ID or GPIB)
- Enclose control default information in parentheses in the control name.

By including default information in the control name, you can later access that information through the help window. This feature is helpful when you are using the VI in higher-level applications.

For example, label a function selector ring control whose default is DC volts as item zero **Function** (DCV:0), and label a Boolean mode switch that defaults to true indicating automatic **Mode** (auto:T). (Notice that the default information is in plain text).

The following figure shows the simple trigger VI after modification to meet the style guidelines:

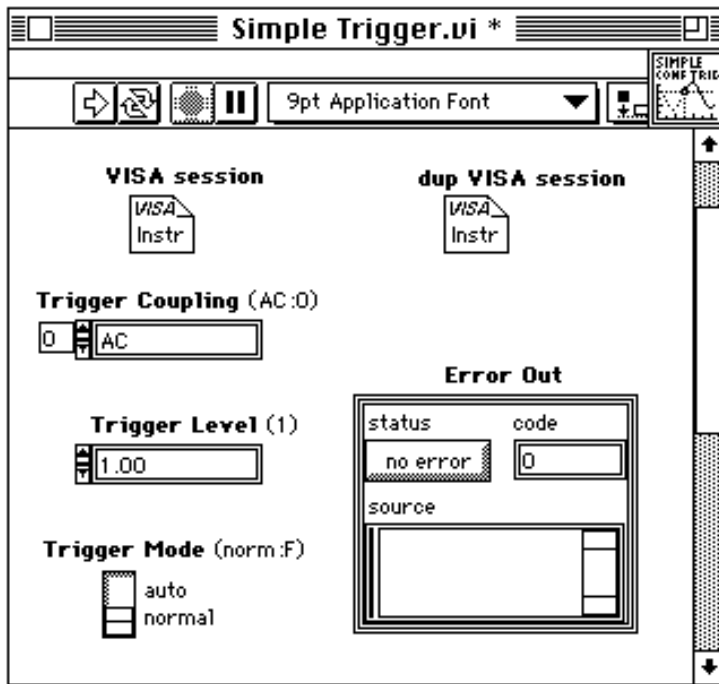


Figure 3-4. Simple Trigger VI Front Panel (See Figure 3-3 for Diagram)

Block Diagram


Proper wiring style improves the diagram appearance and eases understanding. The following are recommendations for developing your instrument driver block diagrams:

- Add text labels to each frame of Case and Sequence structures.
- Label control and indicator nodes with normal text.
- Use bold text to make your free label comments stand out.
- Leave room for labels and wires. Do not crowd the diagram. Do not cover wires with loops, cases, labels, or other diagram objects.
- Reduce the number of bends in the wires by aligning data terminals whenever possible. You can use the cursor keys to move objects one pixel width at a time. Use the **Align** and **Distribute** options in the **Edit** menu to add symmetry and straight lines to your diagram.
- Label long wires and complex operations to increase understandability.

Icon

When you use an instrument driver VI programmatically, the icon graphically represents the function (much like the function name of a C library call). Use meaningful icons for every VI. Include text in the icon that identifies the instrument model controlled by the VI. If you are unable to create an icon to express the function of the VI, you can use text only.

You can borrow icons from similar VIs in other instrument drivers. The following figure shows some examples of icons. These sample icons are available in the file `insticon.llb`.

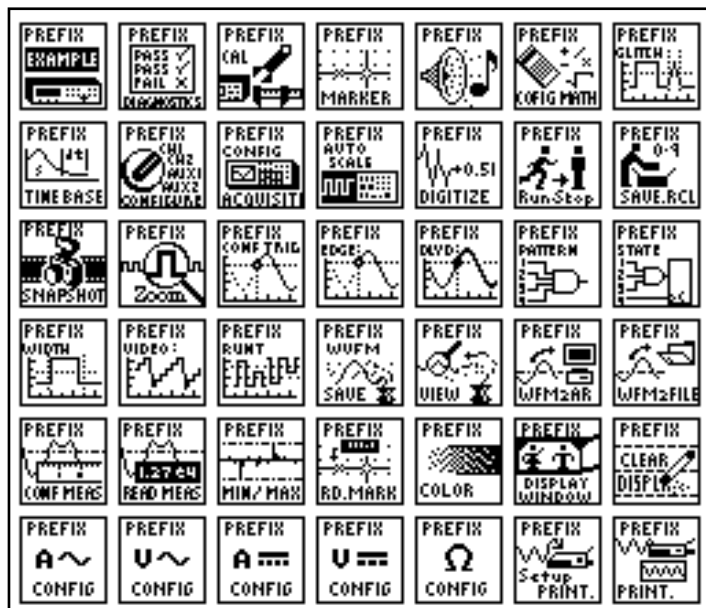

Connector Pane

When you use an instrument driver programmatically, the connector pane defines how to pass parameters to and from the VI. Use the following rules when creating your instrument driver connector panes:

- Place the **VISA session** input and **dup VISA session** output in the upper left and upper right terminals of the LabVIEW instrument driver connector pane.

- Place the **error in** and **error out** clusters in the lower left and lower right terminals of the LabVIEW instrument driver connector pane respectively.
- Place inputs on the left and outputs on the right of the connector pane whenever possible. This promotes a left-to-right data flow when the VI is used in a block diagram.

**Note:**

It is acceptable to choose a connector pane pattern that has extra terminals in case you make unforeseen control or indicator additions to your instrument driver VIs in the future. This procedure prevents you from having to change the pattern and replace all instances of calls to a modified VI.

Error Reporting


Refer to Appendix D, *Error Codes*, for a list of error codes reserved for LabVIEW instrument drivers.

Online Help Information

LabVIEW has two types of help mechanisms available to users: *VI Descriptions* and *Control Descriptions*. You should implement both *VI Descriptions* and *Control Descriptions* for all LabVIEW instrument driver VIs and controls that you develop.

VI Descriptions

Users can access VI Description help from the description box of the information window by selecting **Windows»Show VI Info...**, as shown in the following figure.

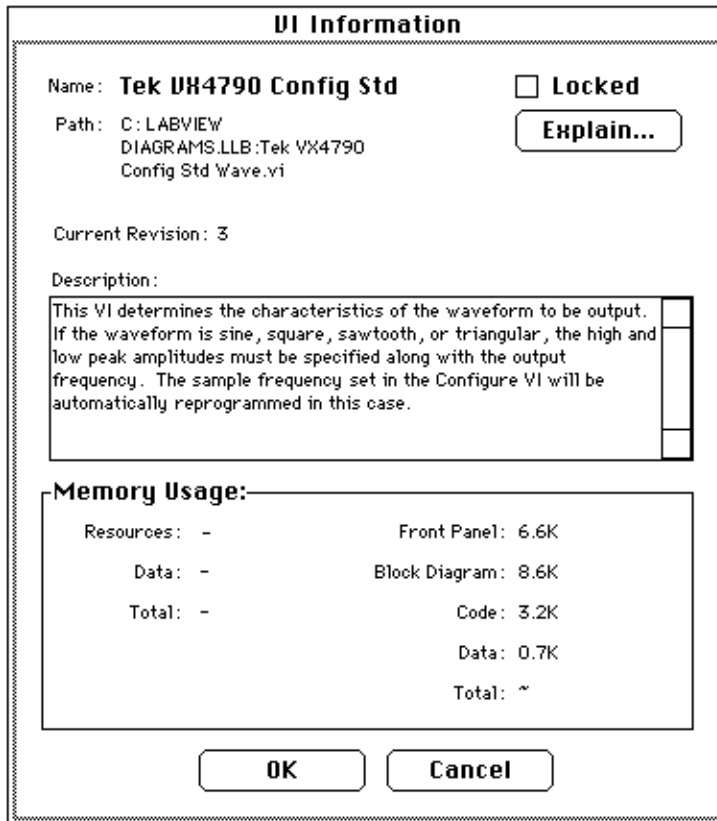


This dialog box should contain the following information:

- A general description of the instrument driver VI
- Control usage rules
- VI interaction with other instrument driver VIs
- Important information concerning the use of the VI

Control and Indicator Descriptions

Control and indicator help is the information most frequently viewed by the user. You can obtain control or indicator help by selecting **Data**

Operations»Description... from the control or indicator pop-up menu, as shown in the following figure:

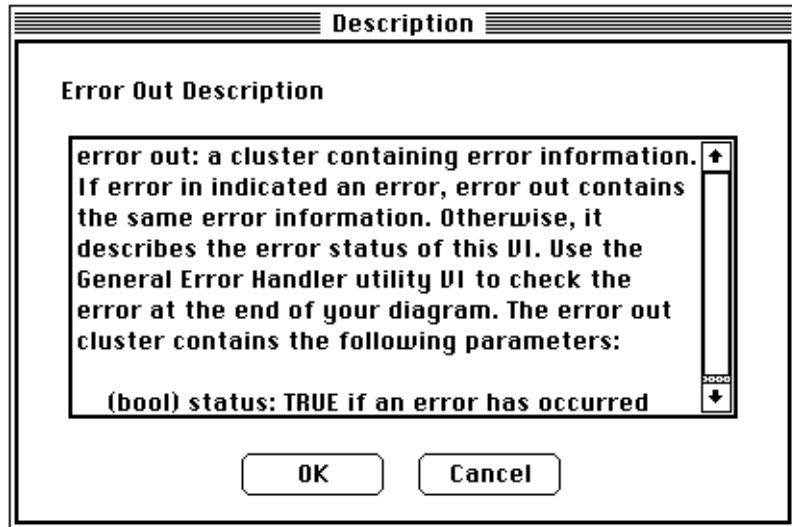

The control and indicator help information should contain the following:

- Name of the parameter
- Brief description of the parameter
- Valid range
- Default value
- Interaction with other controls

Be sure to include information showing index numbers and corresponding settings for all ring and slide controls, and settings corresponding to True/False positions on Boolean controls.

Application VIs

The application VIs demonstrate a common use of the instrument and show how the component VIs are used programmatically to perform a task. For example, an oscilloscope application VI would configure the vertical and horizontal amplifiers, trigger the instrument, acquire a waveform, and report errors. Consider the following points when developing application VIs for your instrument driver:

- Concentrate on building simple, quality examples that can serve as general models for users. It is not necessary to make your application VIs perform every function found in your instrument driver.
- Build the instrument driver top-level examples from the instrument driver component VIs, and perform common test- and measurement-oriented operations for this particular instrument.
- Do not use the instrument driver application VIs to call the initialize or close instrument driver VIs, because doing so will make the application VIs less useful to higher level applications.

LabVIEW Instrument Driver Standards Checklist

All LabVIEW instrument drivers should conform to recommendations for programming style, error handling, front panels, block diagrams, and online help described in this section. Use the following checklist to verify that your instrument driver complies with library standards:

General Issues

- _____ All VIs designed for programmatic use; no pop-up VIs or dialog boxes; no interactive input; all controls wired to connector pane. ('Soft front panel' VIs are allowed only as additional examples built upon the driver VIs; the driver VIs themselves must be for programmatic use.)
- _____ All VIs are multi-instance: no uninitialized shift registers; no global storage VIs unless specifically designed to work with multiple instruments simultaneously.
- _____ All VIs fully documented including **Show VI Info...** and control descriptions.
- _____ Driver uses VISA for all instrument I/O.
- _____ Driver follows instrument internal driver model: Initialize, Configure, Action/Status, Data, Utility, and Close VIs; also Getting Started, Application, and VI Tree VIs.
- _____ All VIs use error I/O clusters.

Files Created

- _____ All VIs in one or more LabVIEW instrument driver archives (* .llb).
 - _____ All .llb files must have instrument prefix in file name (abbreviate as necessary, for example: FL45conf.llb, FL45trig.llb, or FL45meas.llb).
 - _____ Text file allowed with name in the form 'prefix.txt' (example: FL45.txt).
- _____ All .llb and .txt files placed in directory named with instrument prefix.
- _____ Instrument directory and files compressed into 'prefix.zip' using PKzip (for example, FL45.zip).

LabVIEW Instrument Driver Archive - .LLB

- _____ All VIs saved with meaningful names including instrument prefix and description; no special characters; use Initial Cap form (for example, Fluke 45 Read Measurement).
- _____ Function and control menu palettes well organized, and names entered.
 - _____ Getting Started, Application, and VI Tree VIs given top-level status.
 - _____ Support VIs hidden from the function menu palette.

VI Front Panels

- _____ Contain VISA session, dup VISA session, error in, and error out controls.
- _____ Front panel **Show VI Info...** description is complete.
- _____ VI History is updated with comments as needed.
- _____ Control and Indicators:
 - _____ All descriptions are complete. Include defaults, range, items in rings, etc.
 - _____ Place labels at upper left of controls; color labels transparent; use the **Size to Text** feature.

- _____ Capitalize initial letters of control names; use the Application font; bold for primary controls (most) and plain for secondary controls (rarely needed).
- _____ Proper defaults are set for each control and default included in name (example: ring control with default DC Volts selection at value zero: Function Select (DCV:0)).
- _____ The proper display format is used, such as hexadecimal for status registers, etc.
- _____ Use color sparingly or standard gray.

Icon/Connector Pane/VI Setup


- _____ Create meaningful icons for all VIs.
 - _____ Place instrument prefix in upper 1/3 and description text in lower 1/3 of icon.
 - _____ Try to keep a common theme for all VIs of a particular driver (if nothing else, use an image of the instrument).
 - _____ Create icons for Black & White (required), 16-Color and 256-Color (optional).
- _____ Select appropriate connector pane.
 - _____ Try to keep inputs on left & top; outputs on right & bottom.
 - _____ Instrument handle in/out on upper left/right terminal; Error in/out on lower left/right terminal.
 - _____ Extra terminals are acceptable for future modifications.
- _____ Use VI Setups and window options carefully; do not create pop-up panels.

Block Diagram

- _____ Use bold text labels (Application font) to describe each frame of Case and Sequence structures; left-justify; color borders transparent.
- _____ Use plain text labels (Application font) for controls/indicators; place beneath control; left-justify; color borders transparent.
- _____ Don't crowd diagram or cover wires with labels, objects, structures, etc.

- _____ Try for general left-right, top-down data flow.
- _____ Use proper Error I/O wiring techniques; use correct error codes for error reporting.
- _____ Save diagrams with first (or most important) frames or cases showing.

Instrument Driver Template VIs

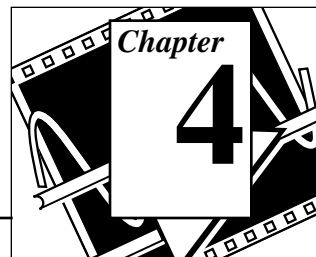

This chapter describes the Instrument Driver Template VIs. These VIs are located in `examples\instr\insttmpl.llb`.

Instrument Driver Template VI Overview

The LabVIEW instrument driver templates are the foundation for all LabVIEW instrument driver development. The templates have a simple, flexible structure and a common set of instrument driver VIs that you can use for driver development. The templates establish a standard format for all LabVIEW drivers and each has instructions for modifying it for a particular instrument. The LabVIEW instrument driver templates contain the following nine predefined template component VIs:

- PREFIX Initialize
- PREFIX Close
- PREFIX Reset
- PREFIX Self Test
- PREFIX Error Query
- PREFIX Error Message
- PREFIX Revision Query
- PREFIX Message-Based Template
- PREFIX Register-Based Template

Rather than developing your own VIs to accomplish these tasks, you should use the LabVIEW instrument driver template VIs which already conform to the LabVIEW standards for instrument drivers. The template VIs are IEEE 488.2 compatible and work with IEEE 488.2 instruments with minimal modifications. For non-IEEE 488.2 instruments, use the template VIs as a shell or pattern, which you can modify by substituting your corresponding instrument-specific commands where applicable. After modifying the VIs, you will have the base level driver that implements all of the template instrument driver VIs for your particular instrument.

Additionally, LabVIEW instrument drivers developed from the template VIs will be similar to other instrument drivers in the library. Therefore, you will have a higher level of familiarity and understanding when you work with multiple instrument drivers.

Initialize

The *Initialize VI* is the first VI called when you are accessing an instrument driver. It configures the communications interface, manages handles, and sends a default command to the instrument. Typically, the default setup configures the instrument operation for the rest of the driver (including turning headers on or off, or using long or short form for queries). After successful operation, the Initialize VI returns a VISA session that addresses the instrument in all subsequent instrument driver VIs.

The VI has an instrument descriptor string as an input. Based on the syntax of this input, the VI configures the I/O interface and generates an instrument handle for all other instrument driver VIs. The following table shows the grammar for the instrument descriptor. Optional parameters are shown in square brackets ([]).

Interface	Syntax
GPIB	GPIB[<i>board</i>][: <i>primary address</i>][: <i>secondary address</i>][:INSTR]
VXI	VXI::VXI <i>logical address</i> [:INSTR]
GPIB-VXI	GPIB-VXI[<i>board</i>][:GPIB-VXI <i>primary address</i>][:VXI <i>logical address</i>][:INSTR]

The GPIB keyword is used with GPIB instruments. The VXI keyword is used for either embedded or MXIbus controllers. The GPIB-VXI keyword is used for a National Instruments GPIB-VXI controller.

Additionally, the Initialize VI can perform selectable ID query and reset operations. In other words, you can disable the ID query when you are attempting to use the driver with a similar but different instrument without modifying the driver source code. Also, you can enable or disable the reset operation. This feature is useful for debugging when resetting would take the instrument out of the state you were trying to test.

Close

All LabVIEW instrument drivers should include a *Close VI*. The Close VI is the last VI called when controlling an instrument. It terminates the software connection to the instrument and deallocates system resources. Additionally, you can choose to place the instrument in an idle state. For example, if you are developing a switch driver, you can disconnect all switches when closing the instrument driver.

Reset

All LabVIEW instrument drivers have a *Reset VI* that places the instrument in a default state. The default state that the Reset VI places the instrument in should be documented in the help information for the Reset VI. In an IEEE 488.2 instrument, this VI sends the command string *RST to the instrument. When you reset the instrument from the Initialize VI, this VI is called. Also, you can call the Reset VI separately.

Self-Test

If an instrument has self-test capability, the LabVIEW instrument driver should contain a *Self-test VI* to instruct the instrument to perform a self-test and return the result of that self-test.

Error Query and Error Message

If an instrument has error query capability, the LabVIEW instrument driver has *Error Query* and *Error Message* VIs. The Error Query VI queries the instrument and returns the instrument-specific error information. The Error Message VI translates the error status information returned from a LabVIEW instrument driver VI into a user-readable string.

Revision Query

LabVIEW instrument drivers have a *Revision Query VI*. This VI outputs the following:

- The revision of the instrument driver.
- The firmware revision of the instrument being used (If the instrument firmware revision cannot be queried, the Revision Query VI should return the literal string `Not Available`.)

Message-Based Template and Register-Based Template

The *Message-Based* and *Register-Based* template VIs are the starting point for developing your own instrument driver VIs. The template VIs have all required instrument driver controls, and instructions for modification for a particular instrument.

Example


The following table shows the LabVIEW VI organization of the instrument driver for the Tektronix VX4790 Arbitrary Function Generator. The hierarchy shows the developer specified VIs derived from the instrument manual as well as the template VIs.

Level	Name	Type
Application VIs	Output Standard Waveform	(Developer Specified)
	Output Arbitrary Waveform	(Developer Specified)
Component VIs	Initialize VI	(Template)
	Configuration VIs <ul style="list-style-type: none"> • Configure Setup • Configure Sample Frequency • Configure Voltage / Attenuation • Configure Standard Waveform 	(Developer Specified) (Developer Specified) (Developer Specified) (Developer Specified)
	Actions VIs <ul style="list-style-type: none"> • Run / Stop 	(Developer Specified)
	Data Functions <ul style="list-style-type: none"> • Download Arbitrary Waveform 	(Developer Specified)
Component VIs	Utilities VIs <ul style="list-style-type: none"> • Reset • Self-Test • Revision Query • Error Query • Error Message • Send Message • Receive Message 	(Template) (Template) (Template) (Template) (Template) (Template) (Template)

Level	Name	Type
	Close VI	(Template)

Instrument Driver Template VI Error Reporting

The instrument driver template VIs report errors in clusters, as the following illustration shows.

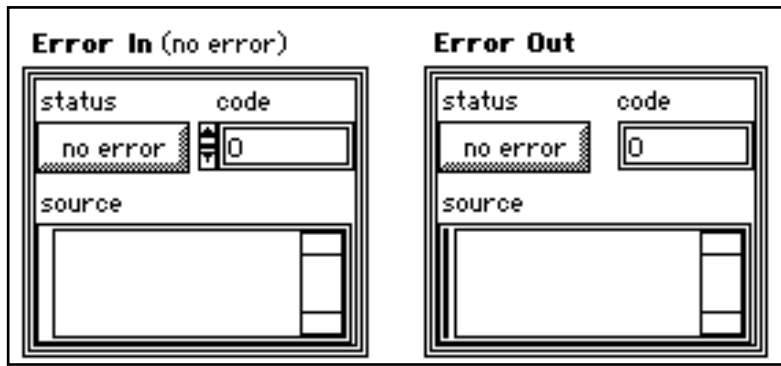

error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. The **error in** cluster contains the following parameters.



status is TRUE if an error occurred. If status is TRUE, this VI does not perform an operation. Instead it places the value of the **error in** cluster in the **error out** cluster.



code is the error code associated with an error. A value of 0 means there is no error. Refer to Appendix D, *Error Codes*, for a description of the possible error codes.



source is the name of the VI that produced the error.



error out is a cluster containing error information. If **error in** indicates an error, the **status**, **code**, and **source** elements of **error out** have the same values as the corresponding elements of **error in**. If **error in** does not indicate an error, **error out** describes the error encountered by the VI. Refer to Appendix D, *Error Codes*, for a description of the possible error codes.


When you use an instrument driver VI in a While Loop, you should stop the loop if the **status** Boolean in the **error out** cluster is TRUE. You can also wire the error cluster to the General Error Handler VI, which can then decipher the error information and inform the user.

The General Error Handler VI is in **Functions»Utilities**. For more information on this VI, refer to the *LabVIEW for Windows User Manual*.

Instrument Driver Template VI Descriptions

PREFIX Close

The PREFIX Close VI is a template for creating a Close VI for your particular instrument. It terminates communication with the instrument and deallocates system resources.

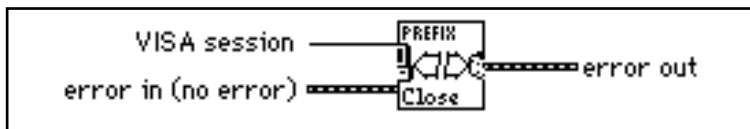

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the close VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Close VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word PREFIX to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram. Replace the word PREFIX in the string constant labeled VI name with the prefix for your instrument.
4. Edit **Show VI Info...** and all control and indicator description information.

5. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O. The Prefix Close VI terminates communication with a device and therefore does not have a **dup VISA session** output.




error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Error Message

The PREFIX Error Message VI is a template for creating an Error Message VI for your particular instrument. It translates the error status information returned from a LabVIEW instrument driver VI to a user-readable string.

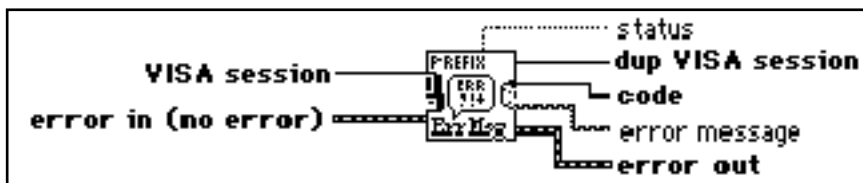

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the error message VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Error Message VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.

2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word PREFIX to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram. Replace the word PREFIX in the string constant labeled **VI name** with the prefix for your instrument.
4. Edit the **user-defined codes** and **user-defined descriptions** array controls that are located on the front panel. Enter instrument specific error codes in the user-defined codes array and enter their corresponding descriptions in the user-defined descriptions array.
5. Select **Operate»Make Current Values Default**.
6. Edit **Show VI Info...** and all control and indicator description information.
7. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.



dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



code contains the error number from the **error in** input cluster.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.




status contains the status from the **error in** input cluster.



error message returns a user-readable string based on the **error in** input cluster.

PREFIX Error Query

The PREFIX Error Query VI is a template for creating an error query VI for your particular instrument. It queries the instrument and returns instrument-specific error information.

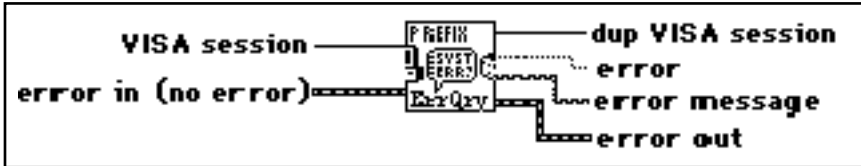

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the error query VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Error Query VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word PREFIX to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram.
 - a. Replace the word PREFIX in the string constant labeled VI name with the prefix for your instrument.
 - b. Edit the string constant labeled error query command by entering the reset command for your instrument.
4. Edit **Show VI Info...** and all control and indicator description information.
5. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.



dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error is the error which has occurred.




error message contains any error messages from the instrument.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Initialize

The PREFIX Initialize VI is a template for creating an Initialize VI for message-based instruments. The function is used to establish communication with a remote instrument. Also the VI can perform user selectable ID query and/or reset operations.

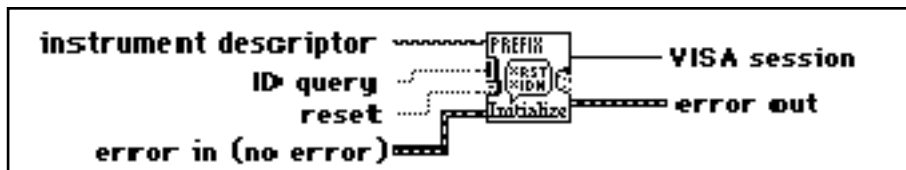

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the Initialize VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Initialize VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.

2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word **PREFIX** to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram.
 - a. Replace the word **PREFIX** in the string constant labeled **VI name** with the prefix for your instrument.
 - b. Set the value of **ID Query command** to be the ID Query for your instrument. Change the constant labeled *ID Query Response* to be the ID Query response for your instrument. If your instrument driver is meant to work with multiple instruments, place all possible ID Query responses in the string constant separated by carriage returns.
4. Replace the **PREFIX Reset VI** with the name of the reset VI for your instrument. If you have not created a reset VI for your instrument, load the **PREFIX Reset VI** in the file `instrtmp.llb` by selecting **File»Open...** and following the instructions on the front panel.
5. Change the string constant labeled **default setup string** by entering the command string you want to send to the instrument.
6. Edit **Show VI Info...** and all control and indicator description information.
7. Delete the **Modification Instructions** text box from the front panel, resize the front panel, and save this VI.



instrument descriptor. Based on the syntax of the **instrument descriptor** input, the **INSTR Open VI** configures the I/O interface and generates an instrument handle that is used by all other instrument driver functions. The grammar for the **instrument descriptor** is shown below. Optional parameters are shown in square brackets ([]).

Interface	Syntax
GPIB	GPIB[<i>board</i>][: <i>primary address</i>][: <i>secondary address</i>][:INSTR]
VXI	VXI::VXI <i>logical address</i> [:INSTR]
GPIB-VXI	GPIB-VXI[<i>board</i>][:GPIB-VXI <i>primary address</i>][:VXI <i>logical address</i>][:INSTR]

As shown in the preceding table, the GPIB keyword is used for GPIB instruments. The VXI keyword is used for VXI instruments through either embedded or MXIbus

controllers. The GPIB-VXI keyword is used for a National Instruments GPIB-VXI controller.

The following table shows the default values for optional parameters:

Optional Parameter	Default Value
board	0
secondary address	none
GPIB-VXI primary address	1



ID query determines if the Initialize VI will attempt to perform identification query.

TRUE: Perform identification query.

FALSE: Do not perform identification query.



reset determines if the Initialize VI will reset the device.

TRUE: Reset the device.

FALSE: Do not reset the device.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.




VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Message-Based Template

The PREFIX Message-Based Template VI is a template for creating a message-based VI for your particular instrument.

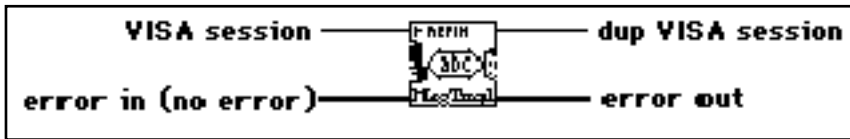

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Rename the VI substituting the prefix for your instrument and a unique name describing the VI's purpose.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the front panel. Add controls and indicators as necessary to perform the VI function.
3. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word **PREFIX** to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
4. Edit the icon connector pane on the front panel.
 - a. Pop-up on the icon pane and select **Show Connector**.
 - b. Pop-up on the connector pane and select **Patterns**. Choose a pattern suitable for the controls and indicators on the front panel. Remember the instrument handle control goes in the upper left-hand corner, and the error in and error out clusters go in the lower left-hand and right-hand corners, respectively.
5. Edit the block diagram. Replace the word **PREFIX** in the string constant labeled **VI name** with the prefix for your instrument.
6. Use the String functions found in the **Functions** palette to build a command string based on the input controls. Wire the resulting command string to the write buffer control of the INSTR Send Message VI.

7. If the instrument generates a response, use the INSTR Receive Message VI to read from the instrument. Use the String functions found in the **Functions** palette to parse the response and wire it to the indicator terminals.
8. Edit **Show VI Info...** and all control and indicator description information.
9. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters




dupVISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Register-Based Template

The PREFIX Register-Based Template VI is a template for creating a register-based VI for your particular instrument.

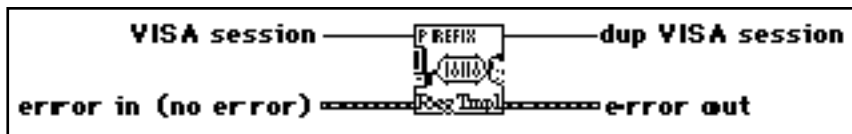

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**

- b. Rename the VI substituting the prefix for your instrument and a unique name describing the VI's purpose.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the front panel. Add controls and indicators as necessary to perform the VI function.
3. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word **PREFIX** to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
4. Edit the icon connector pane on the front panel.
 - a. Pop-up on the icon pane and select **Show Connector**.
 - b. Pop-up on the connector pane and select **Patterns**. Choose a pattern suitable for the controls and indicators on the front panel. Remember the instrument handle control goes in the upper left-hand corner, and the error in and error out clusters go in the lower left-hand and right-hand corners, respectively.
5. Edit the block diagram.
 - a. Replace the word **PREFIX** in the string constant labeled **VI name** with the prefix for your instrument.
 - b. Change the enumerated type labeled **memory space** to be the memory space you plan to access, and change the constant labeled **base address** to be the base address of your instrument.
 - c. Change the constants labeled **register offset** to be the registers that you plan to peek or poke. Add or remove the VISA Peek 16 function or the VISA Poke 16 function, as necessary. If you want to access A24 or A32 addresses, replace these functions with the appropriate VIs from the LabVIEW Instrument Driver Support Library.
6. Edit **Show VI Info...** and all control and indicator description information.
7. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.




dupVISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Reset

The PREFIX Reset VI is a template for creating a reset VI for your particular instrument. The Reset VI places the instrument in a default state.

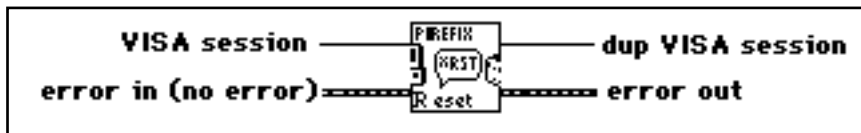

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the reset VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Reset VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word PREFIX to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.

3. Edit the block diagram.
 - a. Replace the word **PREFIX** in the string constant labeled **VI name** with the prefix for your instrument.
 - b. Edit the string constant labeled **Reset Command** by entering the reset command for your instrument.
4. Edit **Show VI Info...** and all control and indicator description information.
5. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters




dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

PREFIX Revision Query

The PREFIX Revision Query VI is a template for creating a Revision Query VI for your particular instrument. It returns the revision of the instrument driver and the firmware revision of the instrument.

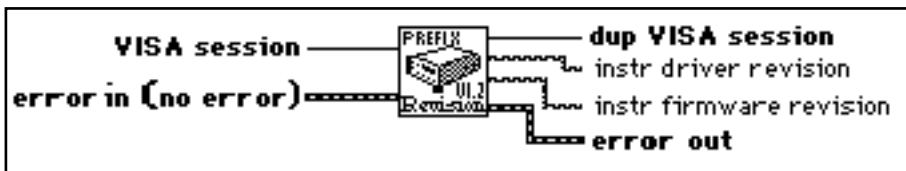

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word **PREFIX** in the name of the VI with the prefix for your instrument. For example, you could name the revision VI for the Tektronix VX4790 Arbitrary Waveform Generator the **TKVX4790 Revision VI**.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word **PREFIX** to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram.
 - a. Replace the word **PREFIX** in the string constant labeled **VI name** with the prefix for your instrument.
 - b. Edit the string constant labeled **revision query command** by entering the revision command for your instrument.
4. Edit **Show VI Info...** and all control and indicator description information.
5. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.



dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.




instr driver revision contains the revision level of the instrument driver source code.



instr firmware revision contains the firmware revision level of the instrument.

PREFIX Self-Test

The PREFIX Self-Test VI is a template for creating a self-test VI for your particular instrument. It tells the instrument to perform a self-test and returns the result.

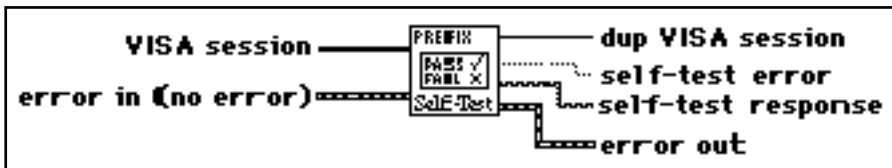

Modify this VI according to the following instructions:

1. Save a copy of this VI.
 - a. Select **File»Save As...**
 - b. Replace the word PREFIX in the name of the VI with the prefix for your instrument. For example, you could name the self-test VI for the Tektronix VX4790 Arbitrary Waveform Generator the TKVX4790 Self-Test VI.
 - c. Save the VI in the file `prefix.llb` where `prefix` is the prefix for your instrument. If you have not already made this `.llb`, you can create it while saving the VI.
2. Edit the icon pane on the front panel.
 - a. Open the icon editor by double-clicking on the icon pane.
 - b. Change the word PREFIX to the prefix for your instrument.
 - c. Edit the icon for all three color modes: Black and White, 16-color, and 256-color.
3. Edit the block diagram.
 - a. Replace the word PREFIX in the string constant labeled **VI name** with the prefix for your instrument.
 - b. Edit the string constant labeled **self-test command** by entering the self-test command for your instrument.
4. Edit **Show VI Info...** and all control and indicator description information.

5. Delete the Modification Instructions text box from the front panel, resize the front panel, and save this VI.



VISA session input is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



error in describes error conditions that occur before this VI executes. The default input of this cluster is `no error`. See the *Instrument Driver Template VI Error Reporting* section of this chapter for a description of the parameters in the **error in** and **error out** clusters.



dup VISA session output is a unique identifier reference to a device I/O session. It identifies the device with which the VI communicates and all configuration information necessary to perform the I/O.



self-test error returns the error code from the instrument.




self-test response contains the result of the instrument driver self-test.




error out contains error information. If **error in** indicates an error, then **error out** contains the same error information. Otherwise it describes the error status that this VI produces.

VISA Library Reference

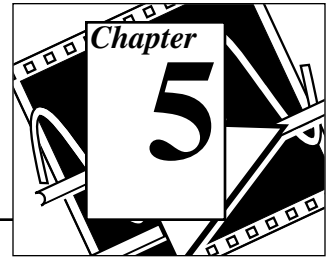

This chapter contains descriptions of the VISA Library Reference operations and attributes.

The following figure shows the **VISA** palette, which you access by selecting **Functions»Instrument I/O»VISA**:

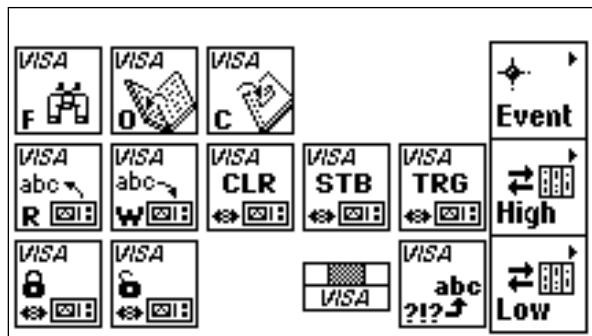

You can find helpful information about individual VIs online by using the LabVIEW Help window (**Help»Show Help**). When you place the cursor on a VI icon, the wiring diagram and parameter names for that VI appear in the Help window. You can also find information for front panel controls or indicators by placing the cursor over the control or indicator with the Help window open. For more information on the LabVIEW Help window, refer to the *Getting Help* section in Chapter 2, *Creating VIs*, of the *LabVIEW User Manual*.



In addition to the Help window, LabVIEW has more extensive online information available. To access this information, select **Help»Online Reference**. For most block diagram objects, you can select **Online Reference** from the object's pop-up menu to access the online description. You can also access this information by pressing the button shown to the left, which is located at the bottom of LabVIEW's Help window. For information on creating your own online reference files, see the *Creating Your Own Help Files* section in Chapter 25, *Managing Your Applications* of the *LabVIEW User Manual*.

Operations

This section describes the VISA Library Reference operations.

VISA Library Reference Parameters

Most of the VISA Library Operations use the following parameters:

- **VISA session** is a unique logical identifier to a session. It is produced by the VISA Open function and used by the VISA primitives. **dup VISA session** is the **VISA session** passed to a function. The dup simplifies dataflow programming and is similar to the dup file refnums provided by file I/O functions.

The **VISA session** drops by default with class *Instr*. You can change the class by popping up on it at edit time. The following classes are currently supported:

- Instr
- GPIB Instr
- VXI/GPIB-VXI RBD Instr
- VXI/GPIB-VXI MBD Instr
- Serial Instr
- Generic Event
- Service Request Event
- Trigger Event
- VXI Signal Event



Note:

The Generic Event, Service Request Event, Trigger Event, and VXI Signal Event classes work only with the VISA Close function and the VISA Attribute Node.

VISA functions vary in the class of **VISA session** which can be wired to them. The valid classes for each function are indicated in the documentation. For example, the functions on the High Level and Low Level Register Access palettes do not accept VISA sessions of class GPIB Instr or Serial Instr. If you wire a **VISA session** to a function that does not accept the class of the session, or if you wire two VISA sessions of differing classes together, your diagram will be broken and the error will be reported as a *Class Conflict*.


- error in** and **error out** terminals comprise the error clusters in each VISA function. The error cluster contains three fields. The **status field** is a Boolean which is TRUE when an error occurs, FALSE when no error occurs. The **code field** will be a VISA error code value if an error occurs during a VISA function. Appendix D lists the VISA Reference Library error codes. The **source field** is a string which describes where the error has occurred. By wiring the **error out** of each function to the **error in** of the next function, the first error condition is recorded and propagated to the end of the diagram where it is reported in only one place.

VISA Operation Descriptions

These functions appear on the main VISA palette. The valid classes for these functions are: Instr (default), GPIB Instr, Serial Instr, VXI/GPIB-VXI RBD Instr, and VXI/GPIB-VXI MBD Instr.

VISA Assert Trigger

Asserts a software or hardware trigger, depending on the interface type.

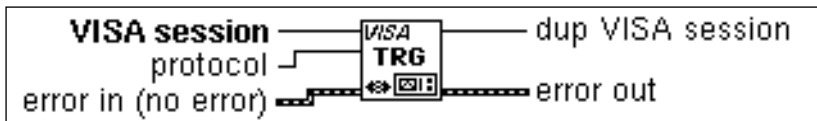

Note: *The Serial Instr class is not valid for VISA Assert Trigger.*



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



protocol is the trigger protocol to use during assertion. Valid protocol values are:

- 0: default
- 1: on
- 2: off
- 5: sync



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Clear

Performs an IEEE 488.1-style clear of the device. For VXI, this is the Word Serial Clear command; for GPIB systems, this is the Selected Device Clear command.

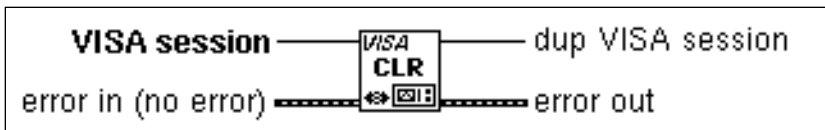

Note: *The Serial Instr class is not valid for VISA Clear.*



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Close

Closes a specified device session or event object. VISA Close accepts all available classes. For a listing of available classes, see the *VISA Operation Parameters* section earlier in this chapter.

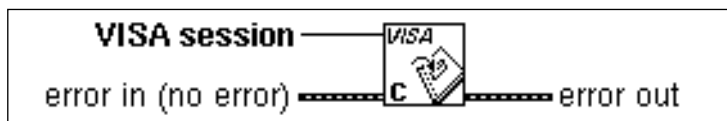



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Find Resource

Queries the system to locate the devices associated with a specified interface.

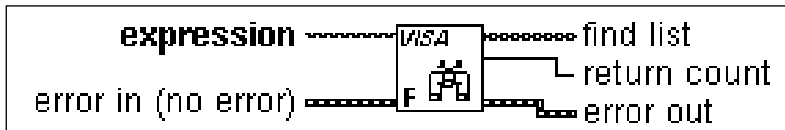

expression matches the value specified with the devices available for a particular interface. The description string specified sets the criteria to search an interface (GPIB, GPIB-VXI, VXI, All VXI, Serial or All) for existing devices. The description string format is shown in the following table:

Interface	Expression
GPIB	GPIB[0-9]*::?*INSTR
GPIB-VXI	GPIB-VXI?*INSTR
GPIB or GPIB-VXI	GPIB?*INSTR
VXI	VXI?*INSTR
All VXI	?*VXI[0-9]*::?*INSTR
Serial	ASRL[0-9]*::?*INSTR
All	?*INSTR



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



find list. An array of strings, each string specifying one interface found by the function.




return count. Specifies the number of matches found; (the number of strings in the **find list** array).



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Lock

Establishes exclusive access to the specified source.

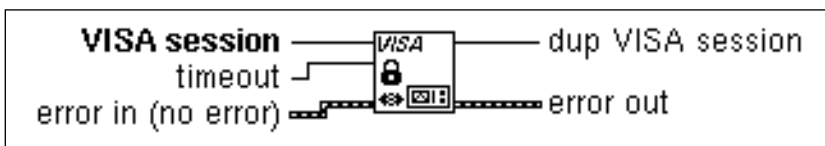

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



timeout. Period of time to wait for access to lock.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Open

Opens a session to the specified device and returns a session identifier that can be used to call any other operations of that device.

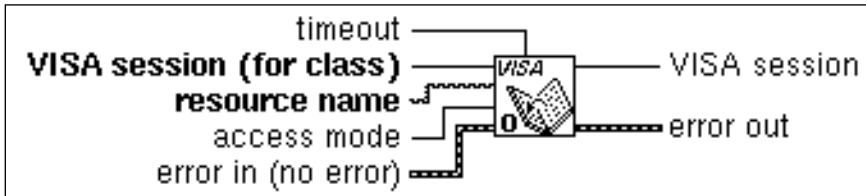

timeout should be set to zero (0).



VISA session is wired to indicate the class of the session being opened. It need not be a valid session, only the class of the session is important.



resource name. Unique symbolic name of a resource.



access mode should be set to zero (0).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



VISA session is the session opened by the function. It must match the class of the input VISA session.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

The following table shows the grammar for the address string. Optional parameters are shown in square brackets ([]).

Interface	Grammar
GPIOB	GPIOB[board]::primary address[::secondary address][::INSTR]
GPIOB-VXI	GPIOB-VXI[board]::VXI logical address[::INSTR]
VXI	VXI[board]::VXI logical address[::INSTR]
Serial	ASRL[board][::INSTR]

The GPIB keyword can be used to establish communication with a GPIB device. The GPIB-VXI keyword is used for a GPIB-VXI controller. The VXI keyword is used for VXI instruments via either embedded or MXIbus controllers. The Serial keyword is used to establish communication with an asynchronous serial (such as RS-232) device.

The INSTR keyword specifies a VISA resource of the type INSTR.

The following table shows the default value for optional parameters:

Optional Parameter	Default Value
board	0
secondary address	none


The following table shows examples of address strings:

Address String	Description
GPIB::1::0::INSTR	A GPIB device at primary address 1 and secondary address 0 in GPIB interface 0.
GPIB-VXI::9::INSTR	A VXI device at logical address 9 in a GPIB-VXI controlled VXI system.
VXI0::1::INSTR	A VXI device at logical address 1 in VXI interface VXI0.
ASRL0::INSTR	A serial device located on port 0.

See the VISA Close description earlier in this chapter.

VISA Read

Reads data from a device. On UNIX platforms data is read synchronously; on all other platforms data is read asynchronously.

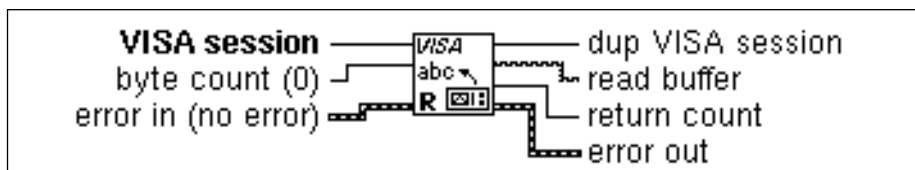



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



byte count (0). Number of bytes to be read.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



read buffer contains the data read from the device.




return count contains the number of bytes actually read.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Read STB

Reads a service request status from a message-based device. For example, on the IEEE 488.2 interface, the message is read by polling devices. For other types of interfaces, a message is sent in response to a service request to retrieve status information. If the status information is only one byte long, the most significant byte is returned with the zero value.

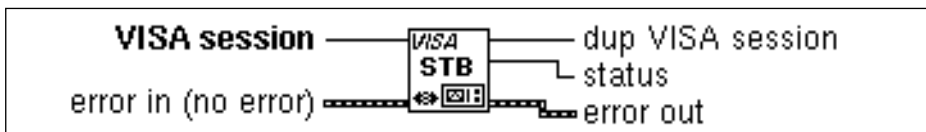

Note: *The Serial Instr class is not valid for VISA Read STB.*



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




status contains the status byte read from the instrument.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Status Description

Retrieves a user-readable string that describes the status code presented in **error in**.

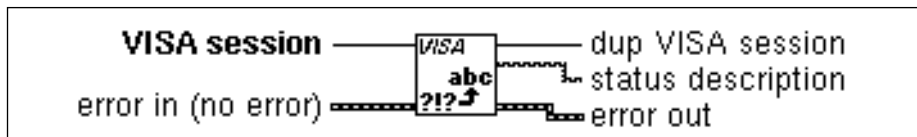

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




status description contains the user-readable string interpretation of the status code passed to the operation.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Unlock

Relinquishes the lock previously obtained using the VISA Lock function.

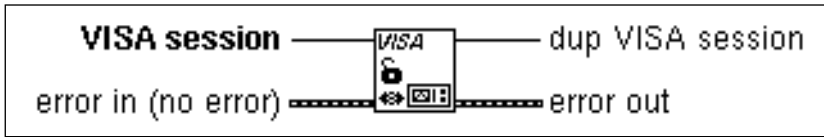

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Write

Writes data to the device. On UNIX platforms data is written synchronously; on all other platforms data is written asynchronously.

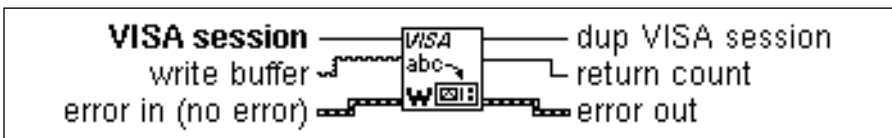

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



write buffer contains the data to be written to the device.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



return count contains the actual number of bytes written.




error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

Event Handling Functions


The following section describes the VISA Event Handling functions. Valid classes for these functions are: Instr (default), GPIB Instr, Serial Instr, VXI/GPIB-VXI RBD Instr, and VXI/GPIB-VXI MBD Instr.

To access the VISA Event Handling functions, pop up on the Event Handling icon on the VISA palette:

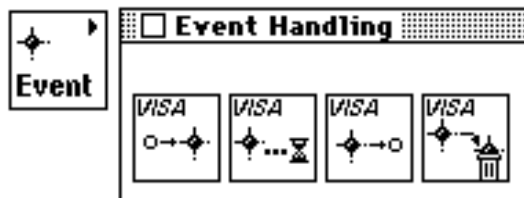

VISA Disable Event

Disables servicing of an event. This operation prevents *new* event occurrences from being queued. However, event occurrences already queued are not lost; use VISA Discard Events if you want to discard queued events.

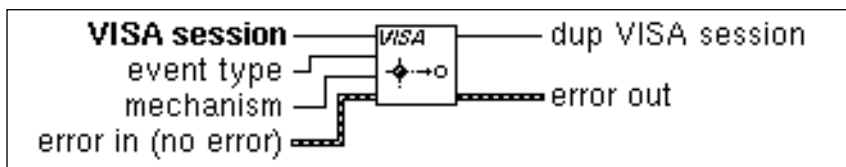

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



event type is the logical event identifier. Special values are:

Trigger:	0xBFFF200A
VXI Signal:	0x3FFF2020
Service Request:	0x3FFF200B
All Events:	0x3FFF7FFF



mechanism. Specifies the event handling mechanism to be disabled. Currently only accepts VI_QUEUE (1).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Discard Events

Discards all pending occurrences of the specified event types and mechanisms from the specified session.

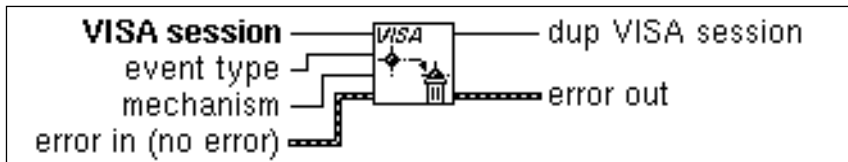

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



event type is the logical event identifier. Special values are:

Trigger:	0xBFFF200A
VXI Signal:	0x3FFF2020
Service Request:	0x3FFF200B
All Events:	0x3FFF7FFF



mechanism. Specifies the event handling mechanism to be disabled. Currently only accepts VI_QUEUE (1).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Enable Event

Enables notification of a specified event.

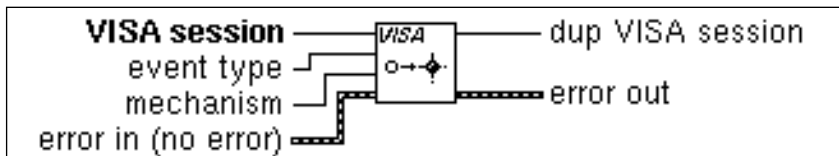

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



event type is the logical event identifier. Special values are:

Trigger:	0xBFFF200A
VXI Signal:	0x3FFF2020
Service Request:	0x3FFF200B



mechanism. Specifies the event handling mechanism to be disabled. Currently only accepts VI_QUEUE (1).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Wait On Event

Suspends execution of a thread of application and waits for an event Event Type for a time period not to exceed that specified by timeout. Refer to individual event descriptions for context definitions. If the specified event type is All Events, the operation waits for any event that is enabled for the given session.

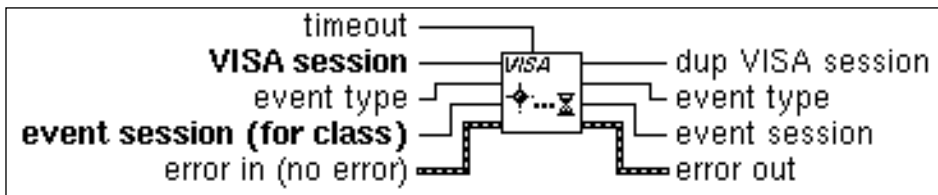

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



event session (for class) specifies the class of event being waited for.



timeout. Period of time to wait for the event.



event type is the logical event identifier. Special values are:

Trigger:	0xBFFF200A
VXI Signal:	0x3FFF2020
Service Request:	0x3FFF200B
All Events:	0x3FFF7FFF



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



event type identifies the event type received if the wait was successful.




event session is valid if the wait was successful. Can be wired to an attribute node to get further information about the event and should be wired to the VISA Close function when examination of the event is complete.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.


High Level Register Access Functions

The following section describes the VISA High Level Register Access functions. Valid classes for these functions are: Instr (default), VXI/GPIB-VXI RBD Instr, and VXI/GPIB-VXI MBD Instr. To access the VISA High Level Register Access functions, pop up on the High Level icon on the VISA palette:

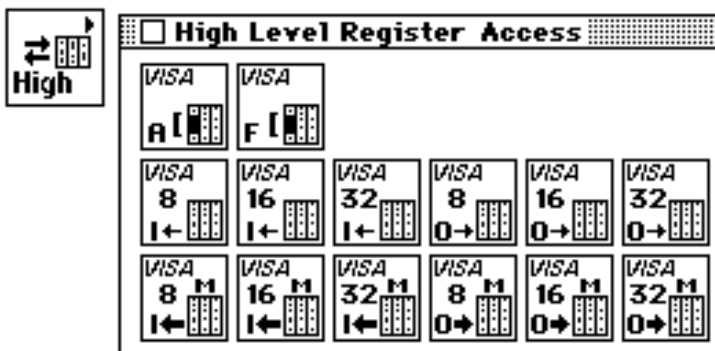

VISA In8 / In16 / In32

Reads in 8-bits, 16-bits, or 32-bits of data, respectively, from the specified memory space (assigned memory base + offset).

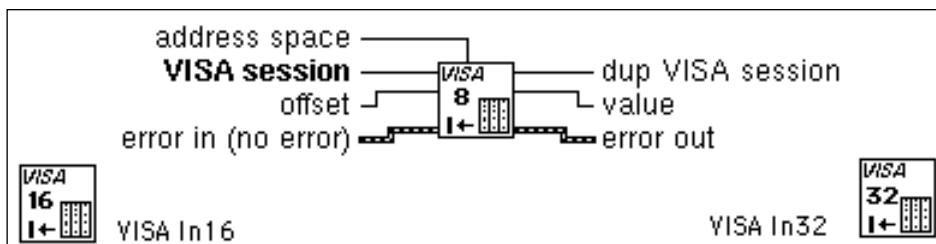



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



address space. Specifies the address space without requiring VISA Map Address to be called. The following table lists the valid entries for specifying address space:

Value	Description
A16 SPACE (1)	Addresses the A16 address space of VXI/MXI bus.
A24 SPACE (2)	Addresses the A24 address space of VXI/MXI bus.
A32 SPACE (3)	Addresses the A32 address space of VXI/MXI bus.



offset. Offset (in bytes) of the device to read from. It is the offset address relative to the device's allocated address base for the corresponding specified address space. For example, if **address space** specifies A16 SPACE, then **offset** specifies the offset from the logical address base address of the VXI device specified. If **address space** specifies A24 or A32 SPACE, **offset** specifies the offset from the base address of the VXI device's memory space allocated by the VXI Resource Manager within VXI A24 or A32 SPACE. See the VISA In description earlier in this section.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




value contains the data read from bus (8-bits for In8, 16-bits for In16, 32-bits for In32; 8-bits is shown).



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Memory Allocation

Returns an offset into a device's region that has been allocated for use by the session. The memory can be allocated on either the device itself or on the computer's system memory.

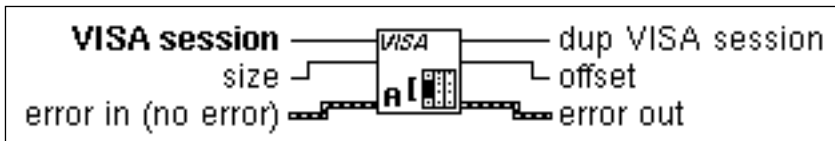

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



size specifies the size of the allocation.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




offset returns the offset which specifies where the memory was allocated.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Memory Free

Frees the memory previously allocated by the VISA Memory Allocation function.

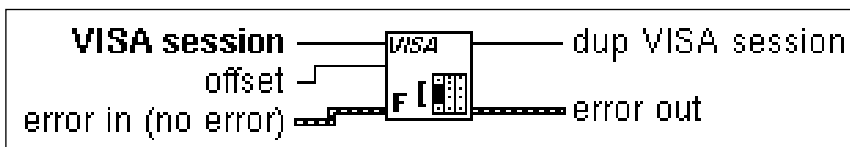

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



offset specifies the memory previously allocated by VISA Memory Allocation.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Move In8 / Move In16 / Move In32

Moves a block of data from device memory to local memory in accesses of 8-bits, 16-bits, or 32-bits, respectively.

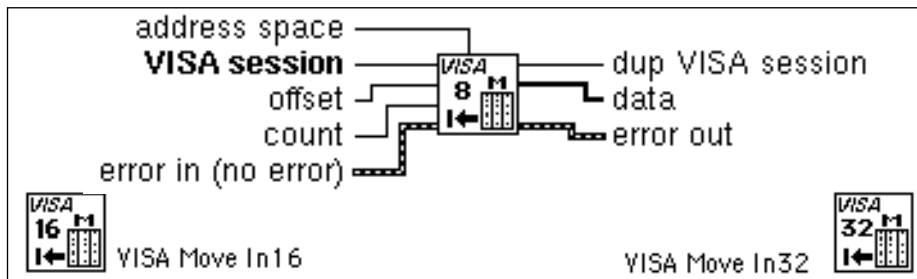

address space. Specifies the address space to map. For more information, see the table listed under **address space**, VISA In.



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



offset. Offset (in bytes) of the device to read from. It is the offset address relative to the device's allocated address base for the corresponding specified address space. For example, if **address space** specifies A16 SPACE, then **offset** specifies the offset from the logical address base address of the VXI device specified. If **address space** specifies A24 or A32 SPACE, **offset** specifies the offset from the base address of the VXI device's memory space allocated by the VXI Resource Manager within VXI A24 or A32 SPACE.



count. Counts the number of data items to move.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




data is the data which is being moved (8-bits for Move In8, 16-bits for Move In16, 32-bits for Move In32; 8-bits is shown).



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Move Out8 / Move Out16 / Move Out32

Moves a block of data from local memory to device memory in accesses of 8-bits, 16-bits, or 32-bits, respectively.

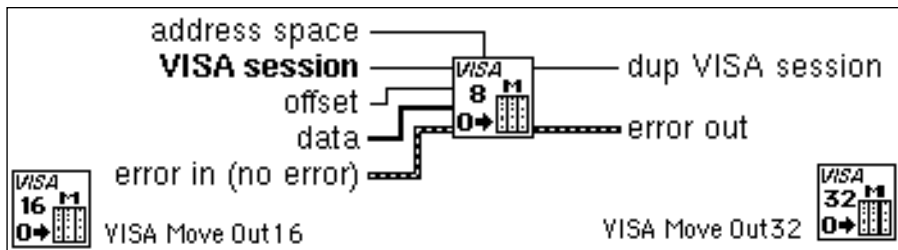

address space. Specifies the address space to map. For more information, see the table listed under **address space**, VISA In.



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



offset. Offset (in bytes) of the device to write to. It is the offset address relative to the device's allocated address base for the corresponding specified address space. For example, if **address space** specifies A16 SPACE, then **offset** specifies the offset from the logical address base address of the VXI device specified. If **address space** specifies A24 or A32 SPACE, **offset** specifies the offset from the base address of the VXI device's memory space allocated by the VXI Resource Manager within VXI A24 or A32 SPACE.



data is the data which is being moved (8-bits for Move Out8, 16-bits for Move Out16, 32-bits for Move Out32; 8-bits is shown).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




count. Counts the number of data items to move.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Out8 / Out16 / Out32

Writes 8-bits, 16-bits, or 32-bits of data, respectively, to the specified memory space (assigned memory base + offset).

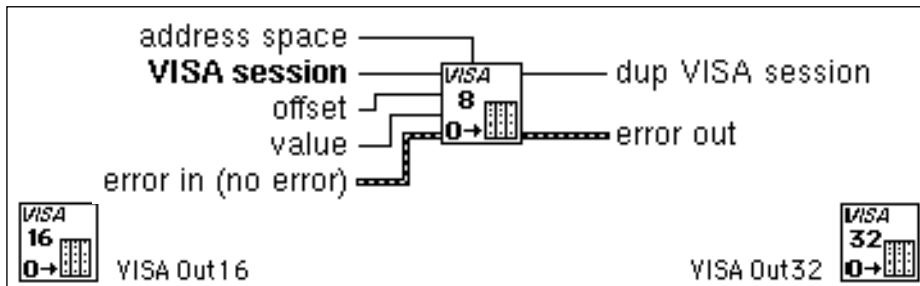

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



address space. Specifies the address space to map. For more information, see the table listed under **address space**, VISA In.



offset. Offset (in bytes) of the device to read from. It is the offset address relative to the device's allocated address base for the corresponding specified address space. For example, if **address space** specifies A16 SPACE, then **offset** specifies the offset from the logical address base address of the VXI device specified. If **address space** specifies A24 or A32 SPACE, **offset** specifies the offset from the base address of the VXI device's memory space allocated by the VXI Resource Manager within VXI A24 or A32 SPACE.



value contains the data to write to bus (8-bits for Out8, 16-bits for Out16, 32-bits for Out32; 8-bits is shown).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.


Low Level Register Access Functions

The following section describes the VISA Low Level Register Access functions. Valid classes for these functions are: Instr (default), VXI/GPIB-VXI RBD Instr, and VXI/GPIB-VXI MBD Instr. To access the VISA Low Level Register Access functions, pop up on the Low Level icon on the VISA palette:

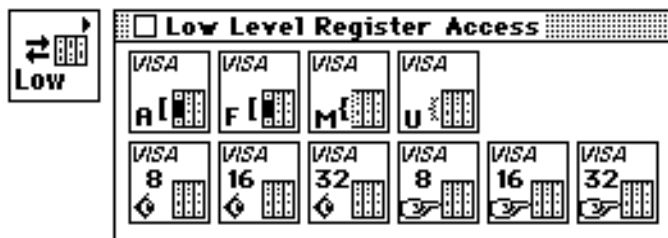

VISA Map Address

Maps in a specified memory space.

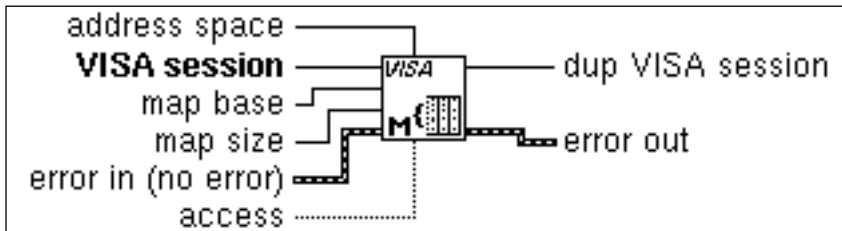

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



address space. Specifies the address space to map. For more information, see the table listed under **address space**, VISA In.



map base contains the offset (in bytes) of the memory to be mapped.



map size contains the amount of memory to map.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



access should be set to `false` (F).



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Memory Allocation


For information on the VISA Memory Allocation function, see the *High Level Register Access Functions* section of this chapter.

VISA Memory Free

For information on the VISA Memory Free function, see the *High Level Register Access Functions* section of this chapter.

VISA Peek8 / Peek16 / Peek32

Reads an 8-bit, 16-bit, or 32-bit value, respectively, from the specified address.

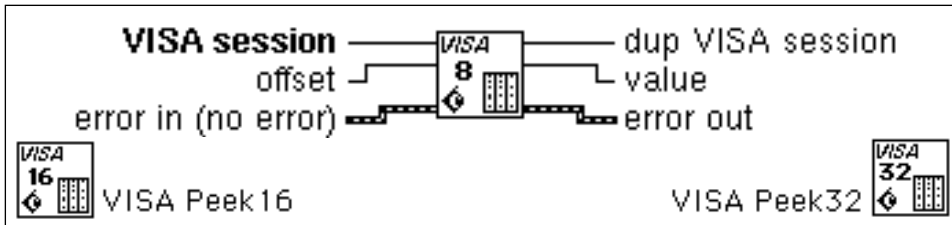

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



offset contains the offset (in bytes) from the register location originally mapped in the mapped memory space.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.



dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.




value contains the read data (8-bits for VISA Peek8, 16-bits for Peek16, 32-bits for Peek32; 8-bits is shown).



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Poke8 / Poke16 / Poke32

Writes an 8-bit, 16-bit, or 32-bit value, respectively, to the specified address.

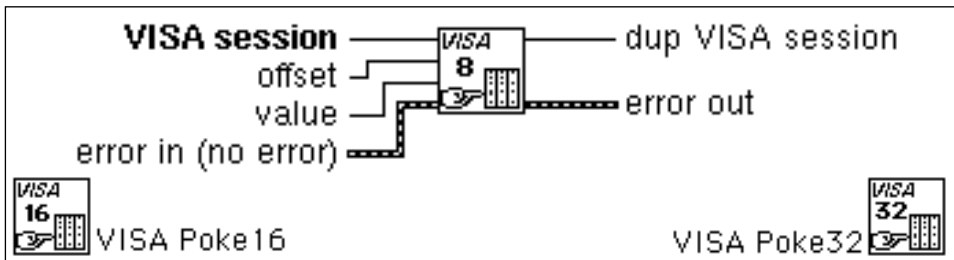

VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



offset contains the offset (in bytes) from the register location originally mapped in the mapped memory space.



value contains the value to write (8-bits for Poke8, 16-bits for Poke16, 32-bits for Poke32; 8-bits is shown).



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Unmap Address

Unmaps memory space previously mapped by VISA Map Address.

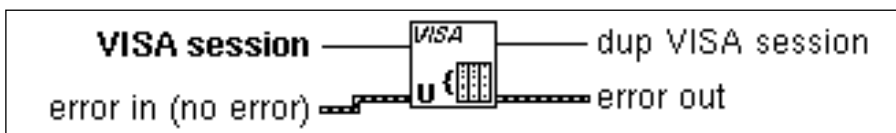



VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error in (no error). See the *VISA Library Reference Parameters* section of this chapter for more information.




dup VISA session. See the *VISA Library Reference Parameters* section of this chapter for more information.



error out. See the *VISA Library Reference Parameters* section of this chapter for more information.

VISA Attribute Node

This section describes the VISA Library attributes. The VISA Attribute Node gets and/or sets the indicated attributes. The node is growable; evaluation starts from the top and proceeds downward until an error, or until the final evaluation, occurs.

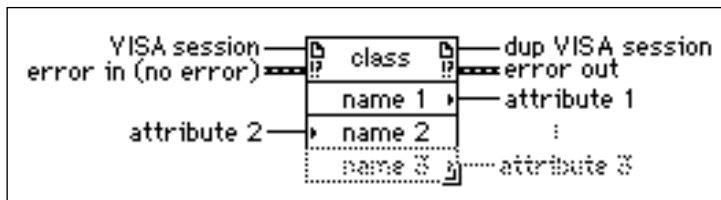

To access the attribute node, select **Functions»Instrument I/O»VISA**. Then select the attribute node icon located on the bottom row of the **VISA** palette.

The VISA Attribute Node only displays attributes for the class of the session that is wired to it. You can change the class of a VISA Attribute Node as long as you have not wired it to a **VISA session**. Once a **VISA session** is wired to a VISA Attribute Node, it adapts to the class of the session and any displayed attributes which are not valid for that class become invalid (this is indicated by turning the attribute item black).

VISA Attribute Descriptions

Each attribute description includes the attribute's range of values, default value, and access privilege. *Local* applies the current session only. *Global* refers to all sessions to the same VISA resource.

Fast Data Channel Mode



Specifies which FDC mode to use (either normal or stream mode).

Range: FDC Normal (1)/FDC Stream (2)
 Default: FDC Normal (1)
 Access Privilege: Read/Write Local

Fast Data Channel Number



Determines which fast data channel (FDC) will be used to transfer the buffer.

Range: 0 to 7
 Default: N/A
 Access Privilege: Read/Write Local

Fast Data Channel Pairs



Specifies use of a channel pair for transferring data; (otherwise, only one channel will be used).

Range: True/False
 Default: False
 Access Privilege: Read/Write Local

Fast Data Channel Signal Enable



Lets the servant send a signal when control of the FDC channel is passed back to the commander. This action frees the commander from having to poll the FDC header while engaging in an FDC transfer.

Range: True/False
 Default: False
 Access Privilege: Read/Write Local

GPIO Primary Address



Specifies the primary address of the GPIO device used by the given session.

Range: 0 to 30
 Default: N/A
 Access Privilege: Read Only Global

GPIO Secondary Address



Specifies the secondary address of the GPIO device used by the given session.

Range: 0 to 31; FFFF for no secondary address
 Default: N/A
 Access Privilege: Read Only Global

IO Protocol



Specifies which protocol to use. In VXI systems you can choose between normal word serial or fast data channel (FDC). In GPIO, you can choose between normal and high-speed (HS488) data transfers.

Range: GPIB: Normal (1)/HS488 (3)
 VXI, GPIB-VXI: Normal (1)/FDC(2)
 Default: Normal (1)
 Access Privilege: Read/Write Local

Immediate Servant



Determines if the VXI device is an immediate servant of the local controller.

Range: True/False
 Default: N/A
 Access Privilege: Read Only Global

Increment Destination Count



Specifies the number of elements by which to increment the destination address on block move operations.

Range: 0, 1
 Default: 1
 Access Privilege: Read/Write Local

Increment Source Count



Specifies the number of elements by which to increment the source address on block move operations.

Range: 0,1
 Default: 1
 Access Privilege: Read/Write Local

Interface Number



Specifies the board number for the given interface.

Range: 0h to FFFFh
 Default: 0
 Access Privilege: Read Only Global

InterfaceType



Specifies the interface type of the given session.

Range: GPIB (1), GPIB-VXI (3), VXI (2), Serial (4)
 Default: N/A
 Access Privilege: Read Only Global

Mainframe Logical Address



Specifies the lowest logical address in the mainframe. If the logical address is not known, UNKNOWN LA is returned.

Range: 0 to 255; UNKNOWN LA (-1)
 Default: N/A
 Access Privilege: Read Only Global

Manufacturer ID



The manufacturer identification number of the VXIbus device.

Range: 0h to FFFh
 Default: N/A
 Access Privilege: Read Only Global

Maximum Queue Length



Specifies the maximum number of events that can be queued at any time on the given session. This attribute is Read/Write until the first time Enable Event is called on a session. Thereafter, this attribute is Read Only.

Range: 1h to FFFFFFFFh
 Default: 50
 Access Privilege: Read/Write Local

Model Code



Specifies the model code for the VXIbus device.

Range: 0h to FFFFh
 Default: N/A
 Access Privilege: Read Only Global

Resource Lock State



Reflects the current locking state of the resource that is associated with the given session.

Range: No Lock (0), Exclusive Lock (1), Shared Lock (2)
 Default: No Lock (0)
 Access Privilege: Read Only Global

Resource Manufacturer Identification



A value corresponding to the VXI manufacturer ID of the manufacturer that created the implementation.

Range: 0h to 3FFFh
 Default: N/A
 Access Privilege: Read Only Global

Resource Manufacturer Name



A string that corresponds to the VXI manufacturer name of the manufacturer that created the implementation.

Range: N/A
 Default: N/A
 Access Privilege: Read Only Global

Resource Name



Unique identifier for a resource.

Range: N/A
 Default: N/A
 Access Privilege: Read Only Global

The address structure is shown in the following table. Optional parameters are shown in square brackets:

Interface	Grammar
GPIOB	<i>GPIOB[board]::primary address[::secondary address][::INSTR]</i>
GPIOB-VXI	<i>GPIOB-VXI[board]::VXI logical address[::INSTR]</i>
VXI	<i>VXI[board]::VXI logical address[::INSTR]</i>
Serial	<i>ASRL[board][::INSTR]</i>

Address String	Description
<code>GPIOB::1::0::INSTR</code>	A GPIOB device at primary address 1 and secondary address 0 in GPIOB interface 0.
<code>GPIOB-VXI::9::INSTR</code>	A VXI device at logical address 9 in a GPIOB-VXI controlled VXI system.

Address String	Description
VXI0::1::INSTR	A VXI device at logical address 1 in VXI interface VXI0.
ASRL0::INSTR	A serial device located on port 0.

Send End Enable



Specifies whether to assert END during the transfer of the last byte of the buffer.

Range: True/False
 Default: True
 Access Privilege: Read/Write Local

Slot



Specifies the physical slot location of the VXIbus device. If the slot number is not known, UNKNOWN SLOT is returned.

Range: 0-12; UNKNOWN SLOT (-1)
 Default: N/A
 Access Privilege: Read Only Global

Suppress End Enable



Specifies whether to suppress the END bit termination. If this attribute is set to TRUE, the END bit does not terminate read operations. If this attribute is set to FALSE, the END bit terminates read operations.

Range: True/False
 Default: False
 Access Privilege: Read/Write Local

Termination Character



The termination character. When the termination character is read and Termination Character Enable is enabled during a read operation, the read operation terminates. See the description for Termination Character Enable listed below.

Range: 0 to FFh
 Default: 0Ah (line feed)
 Access Privilege: Read/Write Local

Termination Character Enable



Determines whether the read operation should terminate when a termination character is received.

Range: True/False
 Default: False
 Access Privilege: Read/Write Local

Timeout Value



Specifies the timeout value to use (in milliseconds) when accessing the device associated with the given session. A timeout value of TMO_IMMEDIATE means that

operations should never wait for the device to respond. A timeout value of TMO INFINITE disables the timeout mechanism.

Range: IMMEDIATE (0); INFINITE (FFFFFFFFh)
Entire range 1 to FFFFFFFEh
Default: 2000
Access Privilege: Read/Write Local

Trigger Identifier



Identifier for the current triggering mechanism.

Range: GPIB: VI_TRIG_SW (-1)
GPIB-VXI, VXI: VI_TRIG_SW (-1);
VI_TRIG_TTL0 to VI_TRIG_TTL7 (0-7);
VI_TRIG_ECL0 to VI_TRIG_ECL1 (8-9)
Default: VI_TRIG_SW (-1)
Access Privilege: Read/Write Local



Note: *Trigger ID is Read/Write when the corresponding session is not enabled to receive trigger events. When the session is enabled to receive trigger events, the attribute is Read Only.*

User Data



Used privately by the application for a particular session. This data is not used by VISA for any purposes. It is provided to the application for its own use.

Range: 0h to FFFFFFFFh
Default: N/A
Access Privilege: Read/Write Local

Version of Implementation



Uniquely identifies each of the different revisions or implementations of a resource. This attribute value is defined by the individual manufacturer and increments with each new revision. The format of the value has the upper 12 bits as the major number of the version, the next lower 12 bits as the minor number of the version, and the lowest 8 bits as the sub-minor number of the version.

Range: 0h to FFFFFFFFh
Default: N/A
Access Privilege: Read Only Global

Version of Specification



Uniquely identifies the version of the VISA specification to which the implementation is compliant. The format of the value has the upper 12 bits as the major number of the version, the next lower 12 bits as the minor number of the version, and the lowest 8 bits as the sub-minor number of the version.

Range: 0h to FFFFFFFFh
Default: 00100000h
Access Privilege: Read Only Global

VXI Commander Logical Address



The logical address of the commander of the VXI device.

Range: 0 to 255
Default: N/A
Access Privilege: Read Only Global

VXI Logical Address



Specifies the logical address of the VXI device used by the given session.

Range: 0 to 255
 Default: N/A
 Access Privilege: Read Only Global

VXI Memory Address Space



Specifies the VXIbus address space used by the device. The three types are A16 only, A16/A24, or A16/A32 memory address space.

Range: A16 (1), A24 (2), A32 (3)
 Default: A16 (1)
 Access Privilege: Read Only Global

VXI Memory Base Address



Specifies the base address of the device in VXIbus memory address space. This base address is applicable to A24 or A32 address space.

Range: 0h to FFFFFFFFh
 Default: N/A
 Access Privilege: Read Only Global

VXI Memory Size



Specifies the size of memory requested by the device in VXIbus address space.

Range: 0h to FFFFFFFFh
 Default: N/A
 Access Privilege: Read Only Global

Window Access



Specifies the modes in which the current window may be accessed.

Range: NMAPPED (1), USE_OPERS (2)

Default: NMAPPED (1)

Access Privilege: Read Only Local

Window Base Address



Specifies the base address of the interface bus to which this window is mapped.

Range: 0h to FFFFFFFFh

Default: N/A

Access Privilege: Read Only Local

Window Size




Specifies the size of the region mapped to this window.

Range: 0h to FFFFFFFFh

Default: N/A


Access Privilege: Read Only Local

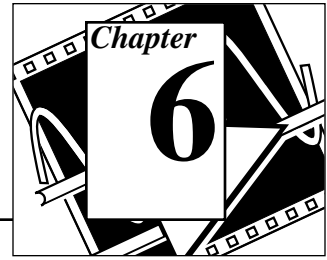

Traditional GPIB Functions

This chapter describes the traditional GPIB functions. You use these functions to handle situations outside the scope of the IEEE 488.2 standard.

The following figure shows the **GPIB Traditional Functions** palette which you access by selecting **Functions»Instrument I/O»GPIB Traditional**:

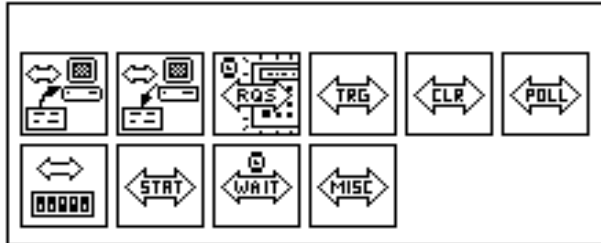

For examples of how to use the traditional GPIB functions, see `examples\instr\smplgpib.llb`.

You can find helpful information about individual VIs online by using the LabVIEW Help window (**Help»Show Help**). When you place the cursor on a VI icon, the wiring diagram and parameter names for that VI appear in the Help window. You can also find information for front panel controls or indicators by placing the cursor over the control or indicator with the Help window open. For more information on the LabVIEW Help window, refer to the *Getting Help* section in Chapter 2, *Creating VIs*, of the *LabVIEW User Manual*.



In addition to the Help window, LabVIEW has more extensive online information available. To access this information, select **Help»Online Reference**. For most block diagram objects, you can select **Online Reference** from the object's pop-up menu to access the online description. You can also access this information by pressing the button shown to the left, which is located at the bottom of LabVIEW's

Help window. For information on creating your own online reference files, see the *Creating Your Own Help Files* section in Chapter 25, *Managing Your Applications* of the *LabVIEW User Manual*.

Traditional GPIB Function Parameters

Most of the traditional GPIB functions use the following parameters:

- **address string** contains the address of the GPIB device with which the function communicates. You can input both the primary and secondary addresses in **address string** by using the form *primary+secondary*. Both *primary* and *secondary* are decimal values, so if *primary* is 2 and *secondary* is 3, **address string** is 2+3.

If you do not specify an address, the functions do not perform addressing before they attempt to read and write the string. They assume you have either sent these commands another way or that another Controller is in charge and therefore responsible for the addressing. If the Controller is supposed to address the device but does not do so before the time limit expires, the functions terminate with GPIB error 6 (timeout) and set bit 14 in **status**. If the GPIB is not the Controller-In-Charge, you must not specify an **address string**.

When there are multiple GPIB Controllers that LabVIEW can use, a prefix to the **address string** in the form ID:address (or ID: if no address is necessary) determines the Controller that a specific function uses. If a Controller ID is not present, the functions assume Controller (or bus) 0.


- **status** is a 16-bit Boolean array in which each bit describes a state of the GPIB Controller. If an error occurs, bit 15 is set. The error code field of the **error out** cluster is a GPIB error code only if bit 15 of **status** is set. Refer to GPIB Status in the *GPIB Function Descriptions* section of this chapter for status bit error codes.
- **error in** and **error out** terminals comprise the error clusters in each traditional GPIB function. The error cluster contains three fields. The **status field** is a Boolean which is TRUE when an error occurs, FALSE when no error occurs. The **code field** will be a GPIB error code value if an error occurs during a GPIB function. Table 6-3 lists the GPIB error codes. The **source field** is a string which describes where the error has occurred. If the **status field** of the **error in** parameter to a function is set, the function is not executed and the same error cluster is passed out. By wiring the **error out** of each function to the **error in** of the next function, the first error

condition is recorded and propagated to the end of the diagram where it is reported in only one place.

Traditional GPIB Function Behavior

The GPIB Read and GPIB Write functions leave the device in the addressed state when they finish executing. If your device cannot tolerate being left in the addressed state, use the GPIB Misc function to send the appropriate unaddress message or configure the NI-488.2 software to unaddress automatically for all devices on the GPIB.


The traditional GPIB Read and Write functions can execute asynchronously. This means other LabVIEW activity can continue while these GPIB functions are operating. When set to execute asynchronously, a small wristwatch icon appears as part of the function icons. A popup item on the Traditional GPIB Read and GPIB Write functions allows for switching their behavior to and from asynchronous operation.

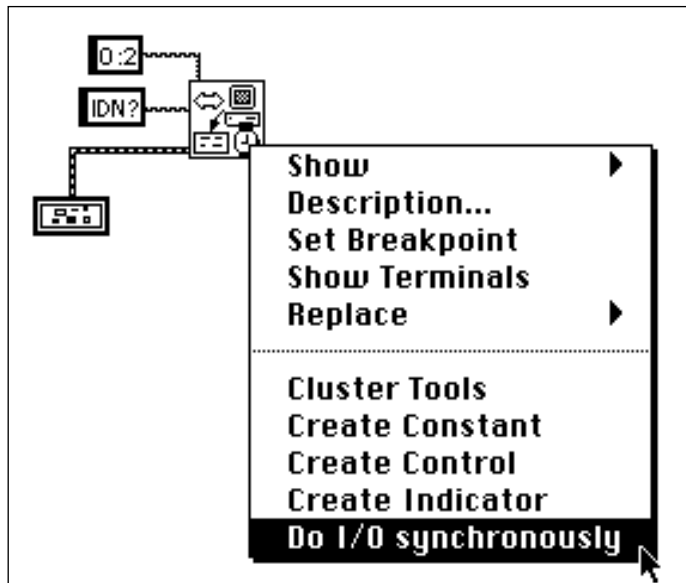

Traditional GPIB Function Descriptions

GPIB Clear

Sends either SDC (Selected Device Clear) or DCL (Device Clear).

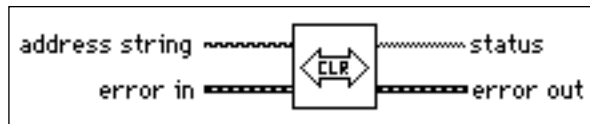

address string. If you specify an address in **address string**, the function sends SDC to that address. If **address string** is unwired, the function sends DCL. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.




status. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Initialization

Configures the GPIB interface at **address string**.

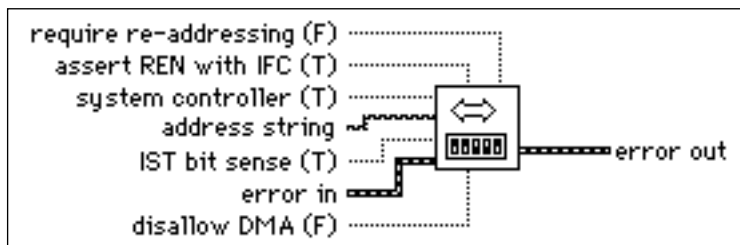



require re-addressing. If **require re-addressing** is TRUE, the function addresses the device before every read or write. If FALSE, the device must be able to retain addressing from one read or write to the next.



assert REN with IFC. If **assert REN with IFC** is TRUE, and if this Controller (specified by the ID in **address string**) is the System Controller, the function asserts the Remote Enable line.



system controller. If **system controller** is TRUE, this Controller acts as the System Controller.



address string. The **address string** sets the GPIB address of the GPIB controller itself. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



IST bit sense. If **IST bit sense** is TRUE, the Individual Status bit of the device responds TRUE to a parallel poll; if **IST bit sense** is FALSE, the Individual Status bit of the device responds FALSE to a parallel poll.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.




disallow DMA. If **disallow DMA** is TRUE, this device uses programmed I/O for data transfers.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Misc

Performs the GPIB operation indicated by **command string**. Use this low-level function when the previously described high-level functions are not suitable.

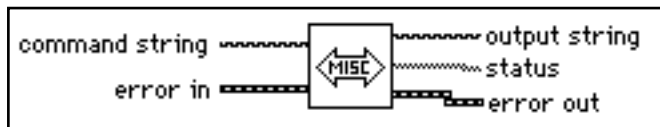



command string. The following table shows the functions you can specify in **command string**. See the *GPIB Device and Controller Functions* section for more information.

Table 6-1. Command String Functions

Device Functions	Description
<i>loc address</i>	Go to local.
<i>off address</i>	Take device offline.
<i>pct address</i>	Pass control.
<i>ppc byte address</i>	Parallel poll configure (enable or disable).
GPIB Controller Functions	Description
<i>cac 0/1</i>	Become active Controller.
<i>cmd string</i>	Send IEEE 488 commands.
<i>dma 0/1</i>	Set DMA mode or programmed I/O mode.
<i>gts 0/1</i>	Go from active Controller to standby.
<i>ist 0/1</i>	Set individual status bit.
<i>llo</i>	Local lockout.
<i>loc</i>	Place Controller in local state.
<i>off</i>	Take controller offline.
<i>ppc byte</i>	Parallel poll configure (enable or disable).
<i>ppu</i>	Parallel poll unconfigure all devices.
<i>rpp</i>	Conduct parallel poll.
<i>rsc 0/1</i>	Request or release system control.
<i>rsv byte</i>	Request service and/or set the serial poll status byte.

Table 6-1. Command String Functions (Continued)

Device Functions	Description
sic	Send interface clear.
sre 0/1	Set or clear remote enable.

To specify the GPIB Controller used by this function, use a *command string* in the form ID: *xxx*, where ID is the GPIB Controller (bus number) and *xxx* is the three-character command and its corresponding arguments, if any. If you do not specify a Controller ID, LabVIEW assumes 0.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



output string is the response the function returns when you conduct a parallel poll with the **rpp** command string.




status describes the state of the GPIB Controller. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Read

Reads **byte count** number of bytes from the GPIB device at **address string**.

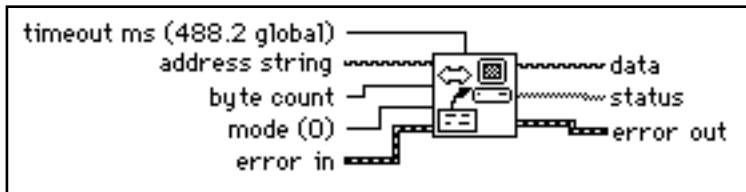

timeout ms. The operation aborts if it does not complete within **timeout ms**. If a timeout occurs, bit 14 of **status** is set. To disable timeouts, set **timeout ms** to 0. To use the 488.2 global timeout, leave this input unwired.

You use the SetTimeout function to change the default value (the 488.2 global timeout) of **timeout ms**. Initially, **timeout ms** defaults to 25,000. See the description of the SetTimeout function in Chapter 7, *GPIB 488.2 Functions*, for more information.



address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



byte count specifies the number of bytes the function reads from the GPIB device.



mode specifies conditions other than reaching **byte count** for terminating the read.

- 0: No EOS character. The EOS termination mode is disabled.
- 1: EOS character is CR. Read terminated on EOI, **byte count**, or CR.
- 2: EOS character is LF. Read terminated on EOI, **byte count**, or LF.
- x: Any other mode indicates the number (decimal) of the desired EOS character. See Appendix A, *Multiline Interface Messages*.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



data string. The function returns the data read in **data string**.



status. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

The GPIB Read Function terminates when it:

- Reads the number of bytes requested.
- Detects an error.
- Exceeds the time limit.
- Detects the END message (EOI asserted).
- Detects the EOS character (if this option is enabled by the value supplied to **mode**).




Note: *The function compares all eight bits when it checks for the EOS character.*



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Serial Poll

Performs a serial poll of the device indicated by **address string**.

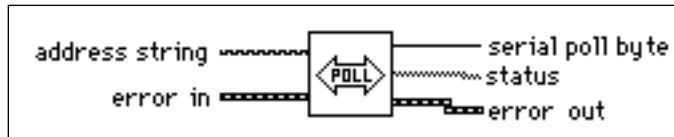

address string is the address of the device the function polls. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



serial poll byte is the response from the device. If the addressed device does not respond within the timeout limit, **serial poll byte** is -1.




status. If the addressed device does not respond within the timeout limit, the function returns GPIB error 6 (bit 14 of **status** is set). See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Status

Shows the status of the GPIB Controller indicated by **address string** after the previous GPIB operation.

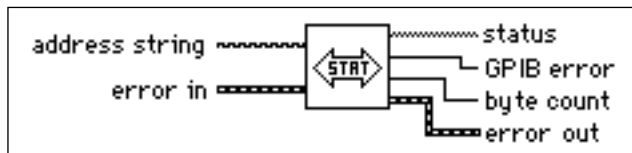



address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



status. The following table shows the numeric value and symbolic status of each bit in **status**. This table also includes a description of each bit.

Table 6-2. Status Bits

Status Bit	Numeric Value	Symbolic Status	Description
0	1	DCAS	Device Clear state
1	2	DTAS	Device Trigger State
2	4	LACS	Listener Active
3	8	TACS	Talker Active
4	16	ATN	Attention Asserted
5	32	CIC	Controller-In-Charge
6	64	REM	Remote State
7	128	LOK	Lockout State
8	256	CMPL	Operation Completed
12	4096	SRQI	SRQ Detected while CIC
13	8192	END	EOI or EOS Detected
14	16384	TIMO	Timeout
15	-32768	ERR	Error Detected



GPIB error contains the most recent error code reported by any of the GPIB functions. The following table shows the possible values for **GPIB error** if bit 15 of **status** is set.

Table 6-3. GPIB Error Codes

GPIB Error	Symbolic Status	Description
0	EDVR	Error Connecting to Driver
1	ECIC	Command Requires GPIB Controller to be CIC
2	ENOL	Write Detected No Listeners
3	EADR	GPIB Controller Not Addressed Correctly
4	EARG	Invalid Argument or Arguments
5	ESAC	Command Requires GPIB to be System Controller
6	EABO	I/O Operation Aborted
7	ENEB	Non-existent Board
8	EDMA	DMA Hardware Not Detected
9	EBTO	DMA Hardware uP Bus Timeout
11	ECAP	No Capability
12	EFSO	File System Operation Error
13	EOWN	Shareable Board Exclusively Owned
14	EBUS	GPIB Bus Error
15	ESTB	Serial Poll Byte Queue Overflow
16	ESRQ	SRQ Stuck On
17	ECMD	Unrecognized Command
19	EBNP	Board Not Present
20	ETAB	Table Error
30	NADDR	No GPIB Address Input

Table 6-3. GPIB Error Codes (Continued)

GPIB Error	Symbolic Status	Description
31	NSTRG	No String Input (Write)
32	NCNT	No Count Input (Read)
61	EPAR	Serial Port Parity Error
62	EORN	Serial Port Overrun
63	EOFL	Port Receive Buffer Overflow
64	EFRM	Serial Port Framing Error
65	SPTMO	Serial Port Timeout, Bytes Not Received at Serial Port




byte count is the number of bytes the previous GPIB operation sent.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Trigger

Sends GET (Group Execute Trigger) to the device indicated by **address string**.

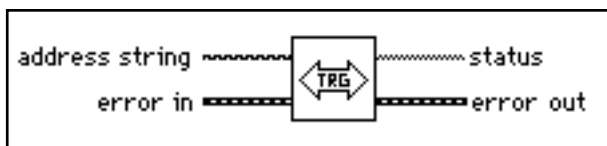

address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.




status. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Wait

Waits for the state(s) indicated by **wait state vector** at the device indicated by **address string**.

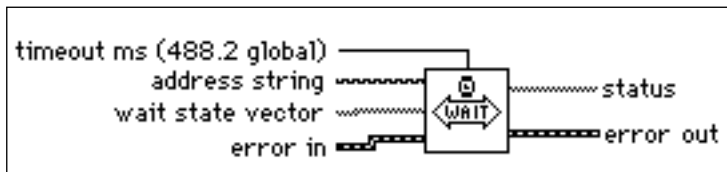

timeout ms. The operation aborts if it does not complete within **timeout ms**. If a timeout occurs, the function sets bit 14 of **status**. To disable timeouts, set **timeout ms** to 0. To use the 488.2 global timeout, leave this input unwired.

You use the SetTimeout function to change the default value (the 488.2 global timeout) of **timeout ms**. Initially, **timeout ms** defaults to 25,000. See the description of the SetTimeout function in Chapter 7, *GPIB 488.2 Functions*, for more information.



address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



wait state vector. Bits set to TRUE in **wait state vector** indicate states for which the function waits. If more than one bit is set, the function terminates when any one of the desired states exists.

The following table defines the bits that you can set in **wait state vector**. This table also lists the numeric value and description of each bit. These bits are the same as the ones that other GPIB functions return; however, only the bits listed are valid for this function.

Table 6-4. Wait State Vector Bits

Wait State Vector Bit	Numeric Value	Symbolic Status	Description
0	1	DCAS	Device Clear state
1	2	DTAS	Device Trigger State
2	4	LACS	Listener Active
3	8	TACS	Talker Active
4	16	ATN	Attention Asserted
5	32	CIC	Controller-In-Charge
6	64	REM	Remote State
7	128	LOK	Lockout State
12	4096	SRQI	SRQ Detected while CIC
13	8192	END	EOI or EOS Detected
14	16384	TIMO	Timeout



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



status. If you are waiting for multiple states, check **status** to see which state caused the function to terminate. Refer to the GPIB Status function for other **status** bit descriptions. See the *Traditional GPIB Function Parameters* section of this chapter for more information.


This function can run in parallel with other functions because LabVIEW alternately checks for status and executes other functions. In addition, multiple calls to this function can execute in parallel, so you can wait for different states on different Controllers at the same time or for multiple states to exist.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

Wait for GPIB RQS

Waits for the device indicated by **address string** to assert SRQ.

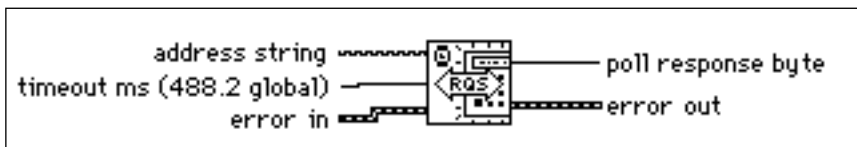

U32

timeout ms. The operation aborts if it does not complete within **timeout ms**. If a timeout occurs, **poll response byte** is -1. To disable timeouts, set **timeout ms** to 0. To use the 488.2 global timeout, leave this input unwired.

You use the SetTimeout function to change the default value (the 488.2 global timeout) of **timeout ms**. Initially, **timeout ms** defaults to 25,000. See the description of the SetTimeout function in Chapter 7, *GPIB 488.2 Functions*, for more information.

abc

address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

ERR

error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

I16

poll response byte. If SRQI is set, the function polls the device at the specified address to see if it requested service. When the specified device requests service (bit 6 is set in **poll response byte**), the function returns the serial poll response.


If the device indicated by **address string** does not respond within the timeout limit, **poll response byte** is -1.

ERR

error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Write

Writes **data** to the GPIB device identified by **address string**.

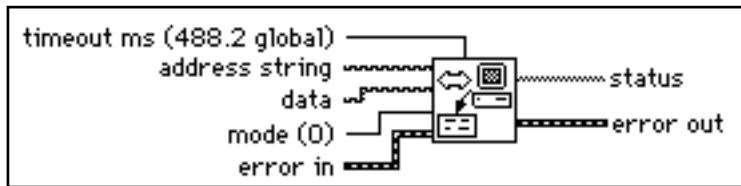

timeout ms. The operation aborts if it does not complete within **timeout ms**. If a timeout occurs, bit 14 of **status** is set. To disable timeouts, set **timeout ms** to 0. To use the 488.2 global timeout, leave this input unwired.

You use the SetTimeout function to change the default value (the 488.2 global timeout) of **timeout ms**. Initially, **timeout ms** defaults to 25,000. See the description of the SetTimeout function in Chapter 7, *GPIB 488.2 Functions*, for more information.



address string. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



data is the data the function writes to the GPIB device.



mode indicates how to terminate the GPIB Write.

- 0: Send EOI with the last character of the string.
- 1: Append CR to the string and send EOI with CR.
- 2: Append LF to the string and send EOI with LF.
- 3: Append CR LF to the string and send EOI with LF.
- 4: Append CR to the string but do not send EOI.
- 5: Append LF to the string but do not send EOI.
- 6: Append CR LF to the string but do not send EOI.
- 7: Do not send EOI.



error in. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



status. See the *Traditional GPIB Function Parameters* section of this chapter for more information.



error out. See the *Traditional GPIB Function Parameters* section of this chapter for more information.

GPIB Device and Controller Functions

This section describes the functions listed in the GPIB Misc function description. The device functions send configuration information to a specific instrument (device). The GPIB Controller functions configure the Controller or send IEEE 488 commands to which all instruments respond. Notice that there are both device and Controller versions of the **ppc** and **loc** commands. The syntax and use of the commands are slightly different for each version.

You can use these functions with all GPIB Controllers accessible by LabVIEW, unless stated otherwise in the function description below. An ECMD error (17) results when you execute a function for a GPIB Controller without the specified capability. The function syntax is strict. Each function recognizes only lowercase characters and allows only one space between the function name and the arguments.

Device Functions

loc – Go to local

syntax

loc *address*

loc temporarily moves devices from a remote program mode to a local mode.

address is the GPIB address of the device. This argument indicates both primary and secondary addresses if you use the form *primary+secondary*, where *primary* and *secondary* are the decimal values of the primary and secondary addresses. For example, if *primary* is 2 and *secondary* is 3, then *address* is 2+3.

loc sends the GTL (Go To Local) message to the GPIB device.

off– Take device offline

syntax **off** *address*

off takes the device at the specified GPIB address offline. This is only needed when sharing a device with another application which is using the NI 488 GPIB Library.

address is the GPIB address of the device. This argument indicates both primary and secondary addresses if you use the form *primary+secondary*, where *primary* and *secondary* are the decimal values of the primary and secondary addresses. For example, if *primary* is 2 and *secondary* is 3, then *address* is 2+3.

pct – Pass control

syntax **pct** *address*

pct passes Controller-in-Charge (CIC) authority to the device at the specified address. The GPIB Controller automatically goes into an idle state. The function assumes that the device to which **pct** passes control has Controller capability.

address is the GPIB address of the device. This argument indicates both primary and secondary addresses if you use the form *primary+secondary*, where *primary* and *secondary* are the decimal values of the primary and secondary addresses. For example, if *primary* is 2 and *secondary* is 3, then *address* is 2+3.

pct sends the following command sequence:

1. Talk address of the device
2. Secondary address of the device, if applicable
3. Take Control (TCT)

ppc – Parallel poll configure

syntax **ppc** *byte address*

ppc enables the instrument to respond to parallel polls.

byte is 0 or a valid parallel poll enable (PPE) command. If *byte* is 0, the parallel poll disable (PPD) byte 0x70 is sent to disable the device from responding to a parallel poll. Each of the 16 PPE messages selects a GPIB data line (DIO1 through DIO8) and sense (1 or 0) that the device must use when it responds to the Identify (IDY) message during a parallel poll. The device compares the **ist** sense and drives the indicated DIO line TRUE or FALSE.

address is the GPIB address of the device. This argument indicates both primary and secondary addresses if you use the form *primary+secondary*, where *primary* and *secondary* are the decimal values of the primary and secondary addresses. For example, if *primary* is 2 and *secondary* is 3, then *address* is 2+3.

Controller Functions

cac – Become active Controller

syntax

cac 0 (take control synchronously)

cac 1 (take control immediately)

cac takes control either synchronously or immediately (and in some cases asynchronously). You generally do not need to use the **cac** function because other functions, such as **cmd** and **rpp**, take control automatically.

If you try to take control synchronously when a data handshake is in progress, the function postpones the take control action until the handshake is complete. If a handshake is not in progress, the function executes the take control action immediately. Synchronous take control is not guaranteed if a read or write operation completes with a timeout or other error.

You should take control asynchronously when it is impossible to gain control synchronously (for example, after a timeout error).

The ECIC error results if the GPIB Controller is not CIC.

cmd – Send IEEE 488 commands

syntax

cmd *string*

cmd sends GPIB command messages. These command messages include device talk and listen addresses, secondary addresses, serial and parallel poll configuration messages, and device clear and trigger messages.

You do *not* use **cmd** to transmit programming instructions to devices. The GPIB Read and GPIB Write functions transmit programming instructions and other device-dependent information.

string contains the command bytes the Controller sends. ASCII characters represent these bytes in **cmd** *string*. If you must send nondisplayable characters, you can enable backslash codes on the string control or string constant or you can use a format function to list the commands in hexadecimal.

dma – Set DMA mode or programmed I/O mode

syntax

dma 0 (use programmed I/O)

dma 1 (use DMA)

dma indicates whether data transfers use DMA.

Some GPIB boards do not have DMA capability. If you try to execute **dma** 1, the function returns GPIB error 11 to indicate no capability.

gts – Go from active Controller to standby

syntax

gts 0 (no shadow handshaking)

gts 1 (shadow handshaking)

Description:

gts sets the GPIB Controller to the Controller Standby state and unasserts the ATN signal if it is the active Controller. Normally, the GPIB Controller is involved in the data transfer. **gts** permits GPIB devices to transfer data without involving the GPIB Controller.

If shadow handshaking is active, the GPIB Controller participates in the GPIB transfer as a Listener, but does not accept any data. When it detects the END message, the GPIB Controller asserts the Not Ready For Data (NRFD) to create a handshake holdoff state.

If shadow handshaking is not active, the GPIB Controller performs neither shadow handshaking nor a handshake holdoff.

If you activate the shadow handshake option, the GPIB Controller participates in a data handshake as a Listener without actually reading the data. It monitors the transfer for the END message and stops subsequent transfers. This mechanism allows the GPIB Controller to take control synchronously on subsequent operations such as **cmd** or **rpp**.

After sending the **gts** command, you should always wait for END before you initiate another GPIB command. You can do this with the GPIB Wait function.

The ECIC error results if the GPIB Controller is not CIC.

ist – Set individual status bit

syntax

ist 0 (individual status bit is cleared)

ist 1 (individual status bit is set)

ist sets the sense of the individual status (**ist**) bit.

You use **ist** when the GPIB Controller is not the CIC but participates in a parallel poll conducted by a device that is the active Controller. The CIC conducts a parallel poll by asserting the EOI and ATN signals, which send the Identify (IDY) message. While this message is active, each device that you configured to participate in the poll responds by asserting a predetermined GPIB data line either TRUE or FALSE, depending on the value of its local **ist** bit. For example, you can assign the GPIB Controller to drive the DIO3 data line TRUE if **ist** is 1 and FALSE if **ist** is 0. Conversely, you can assign it to drive DIO3 TRUE if **ist** is 0 and FALSE if **ist** is 1.

The Parallel Poll Enable (PPE) message in effect for each device determines the relationship among the value of **ist**, the line that is driven, and the sense at which the line is driven. The Controller is

capable of receiving this message either locally via **ppc** or remotely via a command from the CIC. Once the PPE message executes, **ist** changes the sense at which the GPIB Controller drives the line during the parallel poll, and the GPIB Controller can convey a one-bit, device-dependent message to the Controller.

llo – Local lockout

syntax **llo**

llo places all devices in local lockout state. This action usually inhibits recognition of inputs from the front panel of the device.

llo sends the Local Lockout (LLO) command.

loc – Place Controller in local state

syntax **loc**

loc places the GPIB Controller in a local state by sending the local message Return To Local (RTL) if it is not locked in remote mode (indicated by the LOK bit of status). You use **loc** to simulate a front panel RTL switch when you use a computer to simulate an instrument.

off– Take controller offline

syntax **off**

off takes the controller offline. This is only needed when sharing the controller with another application which is using the NI 488 Library.

ppc – Parallel poll configure (enable and disable)

syntax **ppc** *byte*

ppc configures the GPIB Controller to participate in a parallel poll by setting its Local Poll Enable (LPE) message to the value of *byte*. If the value of *byte* is 0, the GPIB Controller unconfigures itself.

Each of the 16 Parallel Poll Enable (PPE) messages selects the GPIB data line (DIO1 through DIO8) and sense (1 or 0) that the device must use when responding to the Identify (IDY) message during a parallel poll. The device interprets the assigned message and the current value of the individual status (**ist**) bit to determine if the selected line is driven TRUE or FALSE. For example, if PPE=0x64, DIO5 is driven TRUE if **ist** is 0 and FALSE if **ist** is 1. If PPE=0x68, DIO1 PPE message is in effect. You must know which PPE and PPD messages are sent and determine what the responses indicate.

ppu – Parallel poll unconfigure

syntax

ppu

ppu disables all devices from responding to parallel polls.

ppu sends the Parallel Poll Unconfigure (PPU) command.

rpp – Conduct parallel poll

syntax

rpp

rpp conducts a parallel poll of previously configured devices by asserting the ATN and EOI signals, which sends the IDY message.

rpp places the parallel poll response in the output string as ASCII characters.

rsc – Release or request system control

syntax

rsc 0 (release system control)

rsc 1 (request system control)

rsc releases or requests the capability of the GPIB Controller to send the Interface Clear (IFC) and Remote Enable (REN) messages to GPIB devices using the **sic** and **sre** functions. For the GPIB Controller to respond to IFC sent by another Controller, the GPIB Controller must not be the System Controller.

In most applications, the GPIB Controller is always the System Controller. You use **rsc** only if the computer is not the System Controller for the duration of the program execution.

rsv – Request service and/or set the serial poll status byte

syntax

rsv *byte*

rsv sets the serial poll status byte of the GPIB Controller to *byte*. If the 0x40 bit is set in *byte*, the GPIB Controller also requests service from the Controller by asserting the GPIB SRQ line. For instance, if you want to assert the GPIB SRQ line, send the ASCII character @, in which the 0x40 bit is set.

You use **rsv** to request service from the Controller using the Service Request (SRQ) signal and to provide a system-dependent status byte when the Controller serial polls the GPIB port.

sic – Send interface clear

syntax

sic


sic causes the Controller to assert the IFC signal for at least 100 msec if the Controller has System Controller authority. This action initializes the GPIB and makes the Controller port CIC. You generally use **sic** when you want a device to become CIC or to clear a bus fault condition.

The IFC signal resets only the GPIB functions of bus devices; it does not reset internal device functions. The Device Clear (DCL) and Selected Device Clear (SDC) commands reset the device functions. Consult the instrument documentation to determine the effect of these messages.

sre – Unassert or assert remote enable*syntax***sre** 0 (unassert Remote Enable)**sre** 1 (assert Remote Enable)

sre unasserts or asserts the GPIB REN line. Devices monitor REN when they select between local and remote modes of operation. A device does not actually enter remote mode until it receives its listen address.


The ESAC error occurs if the Controller is not System Controller.

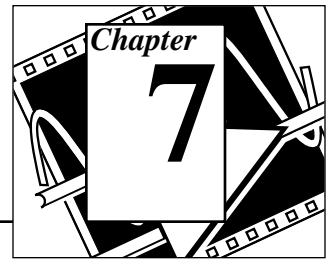

GPIB 488.2 Functions

This chapter describes the IEEE 488.2 (GPIB) Functions.

The following figure shows the **GPIB 488.2** palette which you access by selecting **Functions»Instrument I/O»GPIB 488.2**:

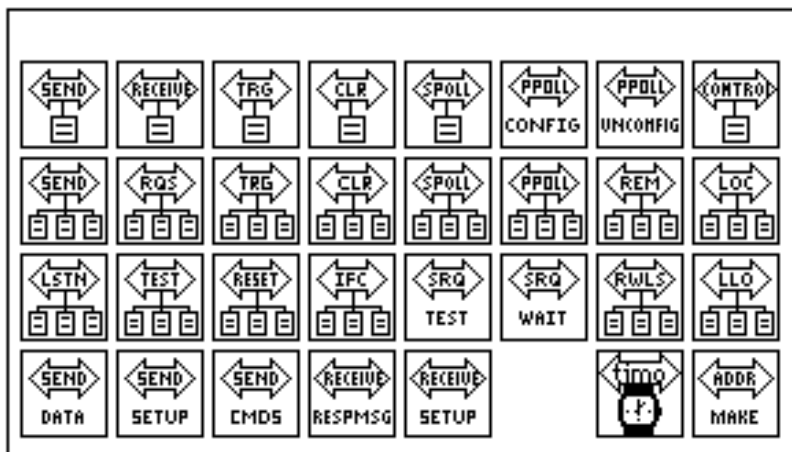

For examples of how to use the GPIB 488.2 Functions, see `examples\instr\smplgpib.llb`.

You can find helpful information about individual VIs online by using the LabVIEW Help window (**Help»Show Help**). When you place the cursor on a VI icon, the wiring diagram and parameter names for that VI appear in the Help window. You can also find information for front panel controls or indicators by placing the cursor over the control or indicator with the Help window open. For more information on the LabVIEW Help window, refer to the *Getting Help* section in Chapter 2, *Creating VIs*, of the *LabVIEW User Manual*.



In addition to the Help window, LabVIEW has more extensive online information available. To access this information, select **Help»Online Reference**. For most block diagram objects, you can select **Online Reference** from the object's pop-up menu to access the online

description. You can also access this information by pressing the button shown to the left, which is located at the bottom of LabVIEW's Help window. For information on creating your own online reference files, see the *Creating Your Own Help Files* section in Chapter 25, *Managing Your Applications* of the *LabVIEW User Manual*

GPIB 488.2 Common Function Parameters

Most of the GPIB 488.2 Functions use the following parameters:


- **address** contains the primary address of the GPIB device with which the function communicates. If a secondary address is required, use the `MakeAddr` function to put the primary and secondary addresses in the proper format. Unless specified otherwise, **address** and **address list** are data types integer and integer array, respectively.
- The default primary address of the GPIB board is 0, with no secondary address. It is designated as System Controller. The default timeout value for the functions is 10 seconds. If you want to change any of these parameters, use the configuration utility included with your GPIB board. You can also use the `GPIB Init` and `SetTimeOut` functions to set the primary address and to change the default timeout value at run time, but these functions affect the interface only when you use it with LabVIEW. For more information, see the documentation supplied with your hardware interface.
- **bus** refers to the GPIB bus number. If you have only one GPIB interface in your computer, the default bus number is 0. For additional GPIB interfaces, see the software installation instructions included with your GPIB board.
- **byte count** refers to the number of bytes that pass over the GPIB.
- **status** is a Boolean array in which each bit describes a state of the GPIB Controller. If an error occurs, the GPIB functions set bit 15. **GPIB error** is valid only if bit 15 of **status** is set. See the status bit and GPIB error tables in the GPIB Status function description in Chapter 6, *Traditional GPIB Functions*.
- **error in; error out**. See the *GPIB Traditional Function Parameters* section of Chapter 6 for a description.

GPIB 488.2 Function Descriptions

Single-Device Functions

DevClear

Clears a single device. To send the Selected Device Clear (SDC) message to several GPIB devices, use the DevClearList function.

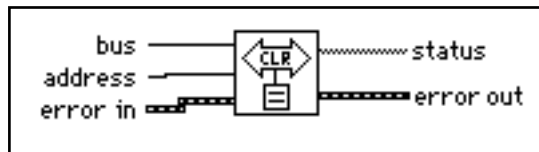

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function sends the GPIB SDC message to the device at **address**. If the **address** terminal is not wired, the function sends the Universal Device Clear message to all devices on the GPIB. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

PPollConfig

Configures a device for parallel polls.

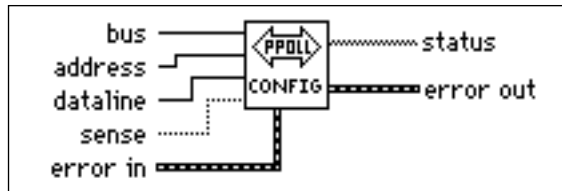

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function configures the GPIB device at **address** for parallel polls according to the **dataline** and **sense** parameters. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



dataline is the data line (1 to 8) on which the device responds to the parallel poll.



sense indicates the condition under which the data line is asserted or unasserted. The device compares this sense value (TRUE or FALSE) to its individual status bit and responds accordingly.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Devices have the option of configuring themselves for parallel polls, in which case they ignore attempts by the Controller to configure them. You should determine whether the device is locally or remotely configured before you use PPollConfig or PPollUnconfig.

PassControl

Passes control to another device with Controller capability.

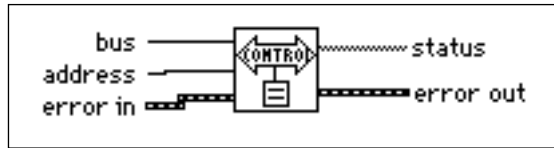

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function sends the GPIB Device Take Control message to the device specified by **address**. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

ReadStatus

Serial polls a single device to get its status byte.

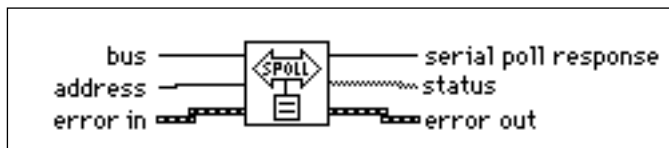

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function serial polls the device at **address**. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



serial poll response is the status byte of the device at **address**. You can learn how to interpret this byte from the documentation supplied with your device.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Receive

Reads data bytes from a GPIB device.

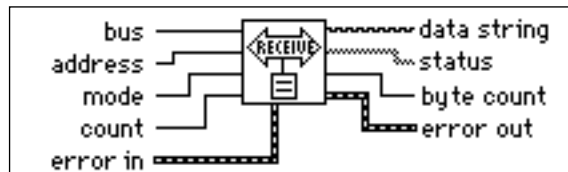

Receive terminates when the function does one of the following:

- reads the number of bytes requested
- detects an error
- exceeds the time limit
- detects the END message (EOI asserted)
- detects the EOS character (assuming the value supplied to **mode** has enabled this option)



bus. The function addresses the interface indicated by **bus** as a Listener. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function addresses the device indicated by **address** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

U16

mode selects the method that signals the end of the data. If **mode** is a value from decimal 0 through 255, the ASCII character that corresponds to it is the termination character, and the function stops the read when it detects the character. If **mode** is not wired, or is decimal 256, the function stops the read when it detects END.

I32

count. The function returns up to **count** data bytes from the device in **data string**.

FF

error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

abc

data string contains **count** data bytes from the GPIB device.

TF

status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

I32


byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

FF

error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Send

Sends data bytes to a single GPIB device.

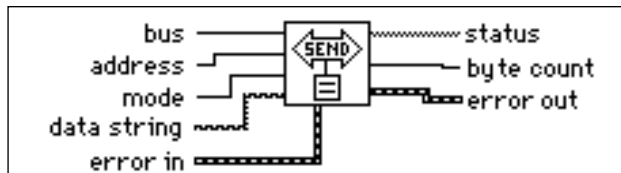**I16**

bus. The function addresses the interface indicated by **bus** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

I16

address. The function addresses the device indicated by **address** as a Listener. If you do not wire **address**, the function sends **data string** with no addressing, thereby sending it to all previously addressed Listeners. To send **data string** to several devices, use the SendList function. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



mode describes how to signal the end of the data to the Listener.

- 0: Do nothing to mark the end of the transfer.
- 1: Send NL (linefeed) with EOI after the data bytes.
- 2: Send EOI with the last data byte in the string.



data string contains the data bytes that the function sends to the GPIB device.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Trigger

Triggers a single device. To send a single message that triggers several GPIB devices, use the TriggerList function.

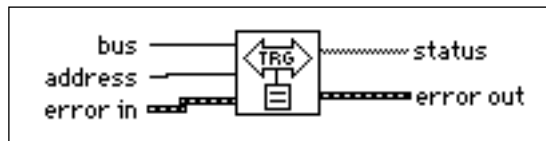

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function sends the GPIB Group Execute Trigger message to the device at **address**. If the **address** terminal is not wired, the function sends the Group Execute Trigger message with no addressing, and this triggers all previously addressed Listeners. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Multiple-Device Functions

AllSPoll

Serial polls all devices.

Although the AllSPoll function is general enough to serial poll any number of GPIB devices, you should use the ReadStatus function when you serial poll only one GPIB device.

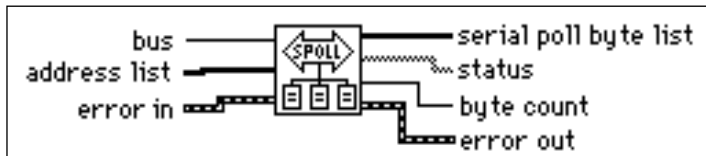

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function serial polls GPIB devices whose addresses appear in **address list**. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



serial poll byte list contains the responses of the devices whose addresses are contained in the corresponding elements of **address list**.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




byte count contains the index of any device that times out instead of responding to the poll.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

DevClearList

Clears multiple devices simultaneously.

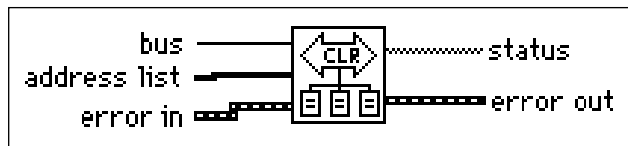

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function sends the GPIB Selected Device Clear (SDC) message to the devices listed in **address list**. If the **address list** terminal is not wired, the function sends the Universal Device Clear message to all devices on the GPIB. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

EnableLocal

Enables local mode for multiple devices.

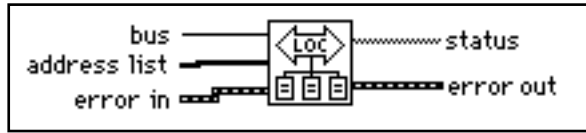

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function sends the GPIB Go To Local (GTL) command to the devices listed in **address list**. If the **address list** terminal is not wired, Remote Enable (REN) unasserts, which places all GPIB devices in local mode. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

EnableRemote

Enables remote programming of multiple GPIB devices.

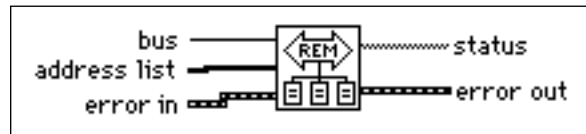

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function sends the GPIB Remote Enable (REN) command to devices listed in **address list**. If the **address list** terminal is not wired, the function

performs no addressing, and the Remote Enable (REN) line asserts. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

FindRQS

Determines which device is requesting service.

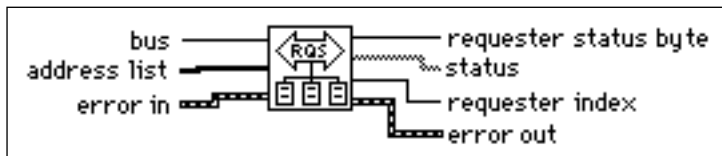

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. Starting at the beginning of the **address list** array, the function serial polls the GPIB devices listed in **address list** until the function finds the device asserting SRQ. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



requester status byte is the status byte for the device asserting SRQ.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




requester index is the index of the address in **address list** of the device asserting SRQ. Arrays are zero-indexed. A requester index of -1 indicates that none of the devices specified in the address list are requesting service.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

PPoll

Performs a parallel poll.

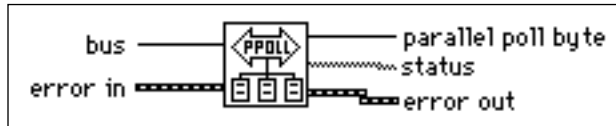

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



parallel poll byte is the result of the parallel poll. Each bit of the poll result returns one bit of status information from each device configured for parallel polls. The state of each bit (0 or 1) and the interpretation of these states are based on the latest parallel poll configuration sent to the devices and on the individual status of the devices.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

PPollUnconfig

Unconfigures devices for parallel polls. The function unconfigures the GPIB devices whose addresses are contained in the **address list** array for parallel polls; that is, they no longer participate in polls.

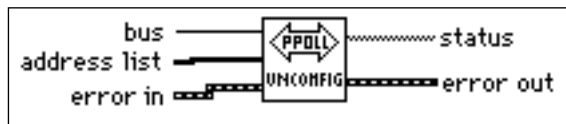

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. If the **address list** terminal is not wired, the function sends the GPIB Parallel Poll Unconfigure (PPU) message, unconfiguring all devices. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendList

Sends data bytes to multiple GPIB devices. This function is similar to Send, except that SendList sends data to multiple Listeners with only one transmission.

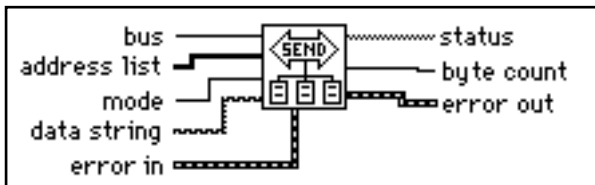

bus. The function addresses the interface indicated by **bus** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[I16]

address list. The **address list** array contains a list of GPIB addresses for devices, which the function addresses as Listeners. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[U8]

mode describes how to signal the end of the data to the Listener.

- 0: Do nothing to mark the end of the transfer.
- 1: Send NL (linefeed) with EOI after the data bytes.
- 2: Send EOI with the last data byte in the string.

[abc]

data string contains the data bytes that the function sends to the addressed devices.

[err]

error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[TF]

status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[I32]


byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[err]

error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

TriggerList

Triggers multiple devices simultaneously.

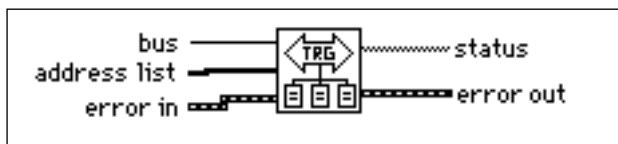**[I16]**

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

[I16]

address list. The function simultaneously triggers GPIB devices whose addresses appear in **address list**. If the **address list** terminal is not wired, the function sends the Group Execute Trigger message without addressing, which triggers all previously

addressed Listeners. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Bus Management Functions

FindLstn

Finds all Listeners on the GPIB. You normally use this function to detect the presence of devices at particular addresses because most GPIB devices have the ability to listen. When you detect them, you can usually interrogate the devices with to determine their identity.

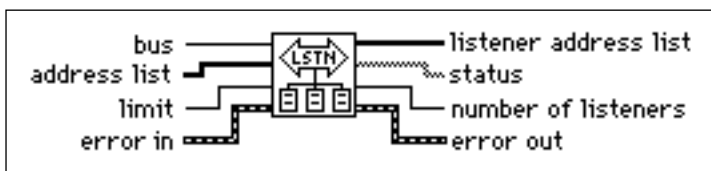

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function tests for the presence of a Listener for all devices listed in **address list**. If the function does not detect a listening device at the particular primary address, it tests all the secondary addresses associated with that primary address. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



limit. If more Listeners are present on the bus than **limit** specifies, the function truncates **listener address list** after **limit** entries.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.


- [I16]** **listener address list** contains the addresses of all the Listeners the function finds.
- [TF]** **status.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.
- [I32]** **number of listeners** contains the number of addresses placed in **listener address list**.
- [F6]** **error out.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

ResetSys

Performs bus initialization, message exchange initialization, and device initialization. First, the function asserts Remote Enable (REN), followed by Interface Clear (IFC), unaddressing all devices and making the GPIB board (the System Controller) the Controller-in-Charge.

Second, the function sends the Device Clear (DCL) message to all connected devices. This ensures that all IEEE 488.2-compatible devices can receive the Reset (RST) message that follows.

Third, the function sends the *RST message to all devices whose addresses are contained in the **address list** array. This message initializes device-specific functions within each device.

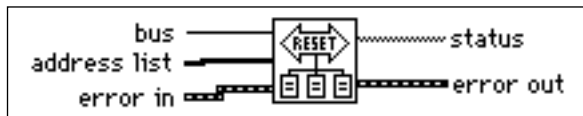

- [I16]** **bus.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.
- [I16]** **address list.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.
- [F6]** **error in.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.
- [TF]** **status.** See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendIFC

Clears the GPIB functions with Interface Clear (IFC). When you issue the GPIB Device IFC message, the interface functions of all connected devices return to their cleared states.

You should use this function as part of a GPIB initialization. It forces the GPIB board to be Controller of the GPIB and ensures that the connected devices are all unaddressed and that the interface functions of the devices are in their idle states.

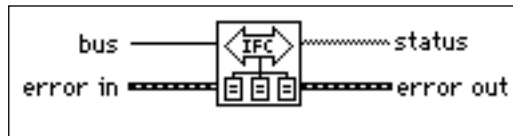

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendLLO

Sends the Local Lockout (LLO) message to all devices. When the function sends the GPIB Local Lockout message, a device cannot independently choose the local or remote state. While Local Lockout is in effect, only the Controller can alter the local or remote state of the devices by sending the appropriate GPIB messages. You should use SendLLO only in unusual local/remote situations, particularly those in which you must

lock all devices into local programming state. Use the SetRWLS Function when you want to place devices in Remote Mode With Lockout State.

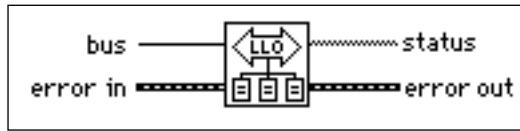

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SetRWLS

Places particular devices in the Remote With Lockout State. The function sends Remote Enable (REN) to the GPIB devices listed in **address list**. It also places all devices in Lockout State, which prevents them from independently returning to local programming mode without intervention by the Controller.

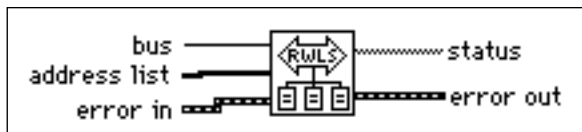

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

TestSRQ

Determines the current state of the SRQ line. This function is similar in format to the WaitSRQ function, except that WaitSRQ suspends itself while it waits for an occurrence of SRQ, and TestSRQ immediately returns the current SRQ state.

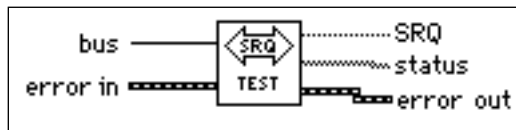

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



SRQ is TRUE if the SRQ line is asserted and FALSE if it is not.




status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

TestSys

Directs multiple devices to conduct IEEE 488.2 self-tests.

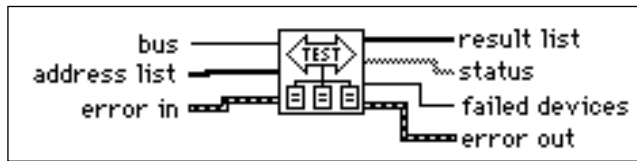

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function simultaneously sends to the GPIB devices whose addresses appear in **address list** a message that instructs them to conduct their self-test procedures. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



result list. Each device returns an integer code signifying the results of its tests, and the function returns these codes in the corresponding elements of **result list**. A result code of 0 indicates that the device passed its tests, and any other value indicates that the tests resulted in an error.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



failed devices contains the number of devices that failed their tests.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

WaitSRQ

Waits until a device asserts Service Request. The function suspends execution until a GPIB device connected on the GPIB asserts the Service Request (SRQ) line.

This function is similar in format to TestSRQ, except that TestSRQ returns the SRQ status immediately, whereas WaitSRQ suspends the program for the duration of the timeout period (but no longer) waiting for an SRQ to occur.

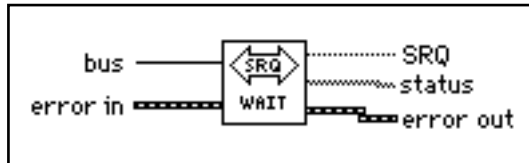

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



SRQ. If the SRQ occurs within the timeout period specified in the GPIB configuration, the function returns TRUE in **SRQ**. Otherwise, the function returns FALSE.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

Low-Level I/O Functions

RcvRespMsg

Reads data bytes from a previously addressed device. This function assumes that another function, such as ReceiveSetup, Receive, or SendCmds, has already addressed the GPIB Talkers and Listeners. You use RcvRespMsg specifically to skip the addressing step of GPIB management. You normally use the Receive function to perform the entire sequence of addressing and then to receive the data bytes.

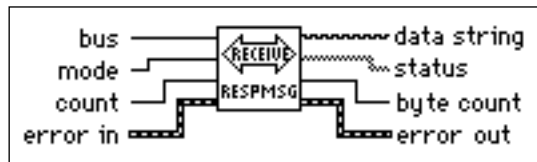

I16

bus. The function addresses the interface indicated by **bus** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

U16

mode determines the data termination method. If **mode** is from 0 through decimal 255, the ASCII character that corresponds to it is the termination character, and the function stops the read when it detects the character. If **mode** is not wired or is decimal 256, the read stops when the function detects the END message (EOI).

I32

count. The function reads up to **count** data bytes from the GPIB.

FF

error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

abc

data string contains the bytes read from the GPIB.

[TF]

status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

I32


byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

FF

error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

ReceiveSetup

Prepares a device to send data bytes and prepares the GPIB board to read data bytes. After you call this function, you can use a function such as RcvRespMsg to transfer the data from the Talker. In this way, you eliminate the need to re-address the devices between blocks of reads.

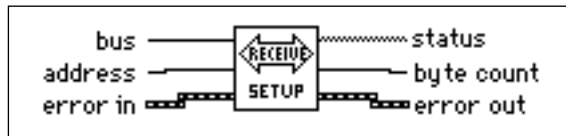

bus. The function addresses the interface indicated by **bus** as a Listener. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address. The function addresses the device specified by **address** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendCmds

Sends GPIB command bytes.

You normally do not need to use SendCmds for GPIB operation. You use it when specialized command sequences, not provided for in other functions, must be sent over the GPIB.

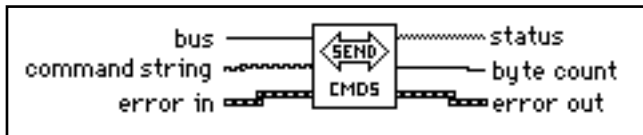

bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



command string contains the command bytes the function sends over the GPIB.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendDataBytes

Sends data bytes to previously addressed devices.

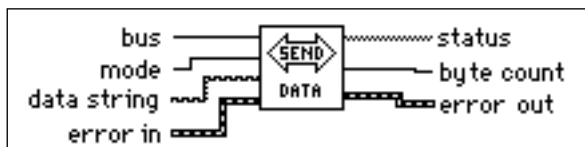



bus. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



mode indicates the method of signaling the end of the data to the Listeners.

- 0: Do nothing to mark the end of the transfer.
- 1: Send NL (linefeed) with EOI after the data bytes.
- 2: Send EOI with the last data byte in the string.



data string contains the data bytes the function sends over the GPIB.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SendSetup

Prepares particular devices to receive data bytes. You normally follow a call to this function with a call to a function such as `SendDataBytes` to actually transfer the data to the Listeners. This sequence eliminates the need to re-address the devices between blocks of sends.

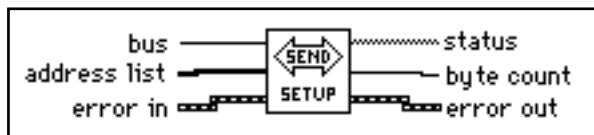

bus. The function addresses the interface indicated by **bus** as a Talker. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



address list. The function addresses GPIB devices whose addresses appear in **address list** as Listeners. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



status. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.



byte count. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

General Functions

MakeAddr

Combines **primary address** and **secondary address** in a specially formatted **packed address** for devices that require both a primary and secondary GPIB address.

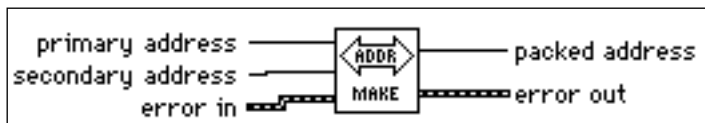

primary address. The function places **primary address** in the lower byte of the **packed address**.



secondary address. The function places **secondary address** in the upper byte of the **packed address**.



error in. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.




packed address contains **primary address** in its lower byte and **secondary address** in its upper byte. You can use **packed address** as the address input to the GPIB 488.2 Functions.



error out. See the *GPIB 488.2 Function Parameters* section of this chapter for more information.

SetTimeout

Changes the global timeout period for all GPIB 488.2 Functions. This function also sets the default timeout period for all GPIB functions.

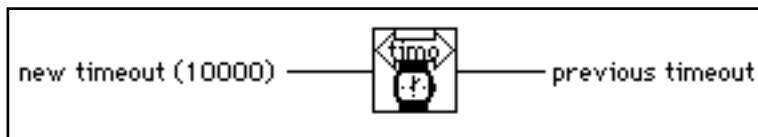


new timeout defaults to 10,000 ms.




previous timeout is the previous timeout period for all GPIB 448.2 Functions.

Serial Port VIs

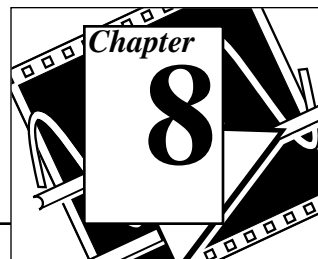

This chapter describes the VIs for serial port operations.

The following figure shows the **Serial** palette which you access by selecting **Functions»Instrument I/O»Serial**:

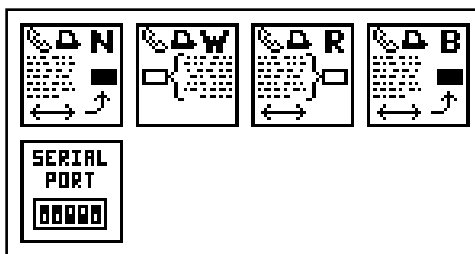

For examples of how to use the Serial Port VIs, see `examples\instr\smp1ser1.llb`.

You can find helpful information about individual VIs online by using the LabVIEW Help window (**Help»Show Help**). When you place the cursor on a VI icon, the wiring diagram and parameter names for that VI appear in the Help window. You can also find information for front panel controls or indicators by placing the cursor over the control or indicator with the Help window open. For more information on the LabVIEW Help window, refer to the *Getting Help* section in Chapter 2, *Creating VIs*, of the *LabVIEW User Manual*.



In addition to the Help window, LabVIEW has more extensive online information available. To access this information, select **Help»Online Reference**. For most block diagram objects, you can select **Online Reference** from the object's pop-up menu to access the online description. You can also access this information by pressing the button shown to the left, which is located at the bottom of LabVIEW's Help window. For information on creating your own online reference files, see the *Creating Your Own Help Files* section in Chapter 25, *Managing Your Applications* of the *LabVIEW User Manual*.

Common Serial Port VI Parameters

Port Number

When you use the serial port VIs under Windows 3.x, the **port number** parameter can have the following values:

0: COM1	5: COM6	10: LPT1
1: COM2	6: COM7	11: LPT2
2: COM3	7: COM8	12: LPT3
3: COM4	8: COM9	13: LPT4
4: COM5		

When you use the serial port VIs under Windows 95 or Windows NT, the **port number** parameter is 0 for COM1, 1 for COM2, and so on.

On a Sun SPARCstation under Solaris 1, the port number parameter for the serial port VIs is 0 for /dev/ttya, 1 for /dev/ttyb, and so on. Under Solaris 2, port 0 refers to /dev/cua/a, 1 to /dev/cua/b, and so on. Under HP-UX port number 0 refers to /dev/tty0p0, 1 to /dev/tty1p0, and so on.

On the Macintosh, port 0 is the modem, using the drivers .ain and .aout. Port 1 is the printer, using the drivers .bin and .bout. To get more ports on a Macintosh, you must install other boards, with the accompanying drivers.

For more information about configuring your serial ports, see the *Serial I/O* section of Appendix C, *Common Questions*.

Handshaking Modes

A common problem in serial communications is ensuring that both sender and receiver keep up with data transmission. The serial port driver can buffer incoming/outgoing information, but that buffer is of a finite size. When it becomes full, the computer ignores new data until you have read enough data out of the buffer to make room for new information.

Handshaking helps prevent this buffer from overflowing. With handshaking, the sender and the receiver notify each other when their buffers fill up. The sender can then stop sending new information until the other end of the serial communication is ready for new data.

You can perform two kinds of handshaking in LabVIEW—software handshaking and hardware handshaking. You can turn both of these forms of handshaking on or off using the Serial Port Init VI. By default, the VIs do not use handshaking.

Software Handshaking–XON/XOFF

XON/XOFF is a software handshaking protocol you can use to avoid overflowing serial port buffers. When the receive buffer is nearly full, the receiver sends XOFF (<control-S> [decimal 19]) to tell the other device to stop sending data. When the receive buffer is sufficiently empty, the receiver sends XON (<control-Q> [decimal 17]) to indicate that transmission can begin again. When you enable XON/XOFF, the devices always interpret <control-Q> and <control-S> as XON and XOFF characters, never as data. When you disable XON/XOFF, you can send <control-Q> and <control-S> as data. Do not use XON/XOFF with binary data transfers because <control-Q> or <control-S> may be embedded in the data, and the devices will interpret them as XON and XOFF instead of as data.

Error Codes


You can connect the **error code** parameter to one of the error handler VIs. These VIs can describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to the *Error Handling* section in Chapter 8, *Time & Dialog Functions*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Serial Port VI Descriptions

Bytes at Serial Port

Returns in **byte count** the number of bytes in the input buffer of the serial port indicated in **port number**.

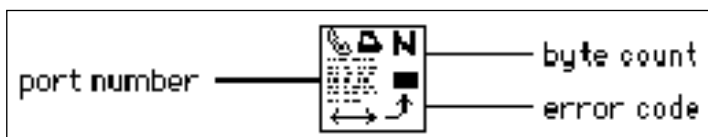



port number. See the *Common Serial Port VI Parameters* section of this chapter for a list of valid port numbers.



byte count is the number of bytes currently queued up in the serial port buffer.




error code. If **error code** is non-zero, an error occurred. Refer to Appendix D, *Error Codes*, for a list of error codes.

You can connect **error code** to one of the error handler VIs. These VIs describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to Chapter 12, *Error Handler VIs*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Serial Port Break

Sends a break on the output port specified by **port number** for a period of time at least as long as the **delay** input requests.

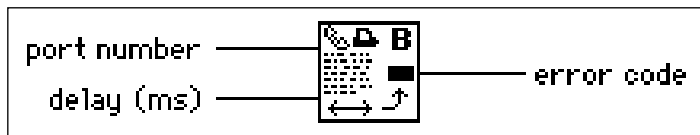

port number. See the *Common Serial Port VI Parameters* section of this chapter for a list of valid port numbers.



delay is the time interval in msec for which the break must last. The actual period of the break can be longer than this interval. On the Sun, **delay** is ignored; the break occurs for at least 0.25 s and no more than 0.5 s.




error code. If **error code** is non-zero, an error occurred. Refer to Appendix D, *Error Codes*, for a list of error codes.

You can connect **error code** to one of the error handler VIs. These VIs describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to Chapter 12, *Error Handler VIs*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Serial Port Init

Initializes the selected serial port to the specified settings.

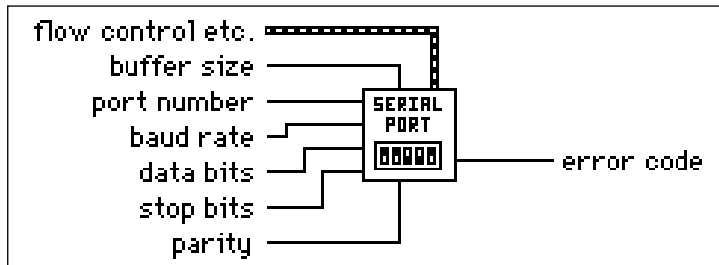

flow control etc. contains the following parameters.



Input XON/XOFF. See the *Common Serial Port VI Parameters* section of this chapter for more information.



Input HW Handshake. On the PC and SPARCstation, this parameter corresponds to Request To Send (RTS) handshaking.



Input alt HW Handshake. On the PC, this parameter corresponds to Data Terminal Ready (DTR) handshaking. On the SPARCstation, this parameter is ignored.



Output XON/XOFF. See the *Common Serial Port VI Parameters* section of this chapter for more information.



Output HW Handshake. On the PC and SPARCstation, this parameter corresponds to Clear to Send (CTS) handshaking.



Output alt HW Handshake. On the PC, this parameter corresponds to Data Set Ready (DSR) handshaking. On the SPARCstation, this parameter is ignored.



XOFF byte is the byte used for XOFF (^S).



XON byte is the byte used for XON (^Q).



Parity Error Byte. If the high byte is non-zero, the low byte is the character that is used to replace any parity errors found when **parity** is enabled.



buffer size indicates the size of the input and output buffers the VI allocates for communication through the specified port. If **buffer size** is less than or equal to 1 K, the VI uses 1 K as the **buffer size**. You may need to use larger buffers for large data transfers. The buffer size is in bytes.



port number. See the *Common Serial Port VI Parameters* section of this chapter for more information.



baud rate is the rate of transmission.



data bits is the number of bits in the incoming data. The value of **data bits** is between five and eight.



stop bits is 0 for one stop bit, 1 for one-and-a-half stop bits, or 2 for two stop bits.



parity is 0 for no parity, 1 for odd parity, 2 for even parity, 3 for mark parity, or 4 for space parity.




error code is -1 if **baud rate**, **data bits**, **stop bits**, **parity**, or **port number** are out of range, or if the serial port could not be initialized. Check the values of **baud rate**, **data bits**, **stop bits**, **parity**, and **port number**. If these values are valid, verify that the serial port has been initialized. Refer to Appendix D, *Error Codes*, for a list of error codes.

You can connect **error code** to one of the error handler VIs. These VIs can describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to Chapter 12, *Error Handler VIs*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Serial Port Read

Reads the number of characters specified by **requested byte count** from the serial port indicated in **port number**.

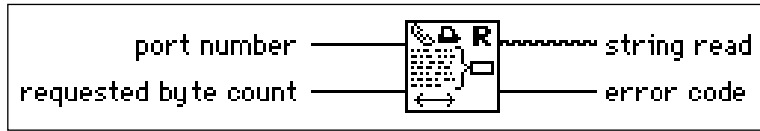

port number. See the *Common Serial Port VI Parameters* section of this chapter for a list of valid port numbers.



requested byte count specifies the number of characters to be read. If you want to read all of the characters currently at the serial port, first execute the Bytes at Serial Port VI to determine the exact number of bytes ready to be read. Then use the **byte count** output of that VI as the **requested byte count** input to the Serial Port Read VI.



string read. The VI returns the bytes read in **string read**.




error code. If **error code** is non-zero, an error occurred. Refer to Appendix D, *Error Codes*, for a list of error codes.

You can connect **error code** to one of the error handler VIs. These VIs describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to Chapter 12, *Error Handler VIs*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Serial Port Write

Writes the data in **string to write** to the serial port indicated in **port number**.

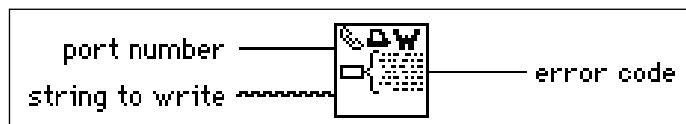



port number. See the *Common Serial Port VI Parameters* section of this chapter for a list of valid port numbers.



string to write is the data to be written to the serial port.



error code. If **error code** is non-zero, an error occurred. Refer to Appendix D, *Error Codes*, for a list of error codes.

You can connect **error code** to one of the error handler VIs. These VIs describe the error and give you options on how to proceed when an error occurs. For more information on using the error handler VIs, refer to Chapter 12, *Error Handler VIs*, of the *LabVIEW Function Reference Manual*.

Some error codes returned by the serial port VIs are platform-specific. Please refer to your system documentation for a list of error codes.

Multiline Interface Messages



This appendix lists multiline interface messages, which are commands that IEEE 488 defines. Multiline interface messages manage the GPIB—they perform tasks such as initializing the bus, addressing and unaddressing devices, and setting device modes for local or remote programming. These multiline interface messages are sent and received with ATN TRUE. The following list includes the mnemonics and messages that correspond to the interface functions.

For more information on these messages, refer to the ANSI/IEEE Standard 488.1-1987, *IEEE Standard Digital Interface for Programmable Instrumentation*.

Multiline Interface Messages

Hex	Oct	Dec	ASCII	Msg
00	000	0	NUL	
01	001	1	SOH	GTL
02	002	2	STX	
03	003	3	ETX	
04	004	4	EOT	SDC
05	005	5	ENQ	PPC
06	006	6	ACK	
07	007	7	BEL	
08	010	8	BS	GET
09	011	9	HT	TCT
0A	012	10	LF	
0B	013	11	VT	
0C	014	12	FF	
0D	015	13	CR	
0E	016	14	SO	
0F	017	15	SI	
10	020	16	DLE	
11	021	17	DC1	LLO
12	022	18	DC2	
13	023	19	DC3	

Hex	Oct	Dec	ASCII	Msg
14	024	20	DC4	DCL
15	025	21	NAK	PPU
16	026	22	SYN	
17	027	23	ETB	
18	030	24	CAN	SPE
19	031	25	EM	SPD
1A	032	26	SUB	
1B	033	27	ESC	
1C	034	28	FS	
1D	035	29	GS	
1E	036	30	RS	
1F	037	31	US	
20	040	32	SP	MLA0
21	041	33	!	MLA1
22	042	34	“	MLA2
23	043	35	#	MLA3
24	044	36	\$	MLA4
25	045	37	%	MLA5
26	046	38	&	MLA6
27	047	39	,	MLA7
28	050	40	(MLA8
29	051	41)	MLA9

Hex	Oct	Dec	ASCII	Msg
2A	052	42	*	MLA10
2B	053	43	+	MLA11
2C	054	44	,	MLA12
2D	055	45	-	MLA13
2E	056	46	.	MLA14
2F	057	47	/	MLA15
30	060	48	0	MLA16
31	061	49	1	MLA17
32	062	50	2	MLA18
33	063	51	3	MLA19
34	064	52	4	MLA20
35	065	53	5	MLA21
36	066	54	6	MLA22
37	067	55	7	MLA23
38	070	56	8	MLA24
39	071	57	9	MLA25
3A	072	58	:	MLA26
3B	073	59	;	MLA27
3C	074	60	<	MLA28
3D	075	61	=	MLA29
3E	076	62	>	MLA30
3F	077	63	?	UNL

Hex	Oct	Dec	ASCII	Msg
40	100	64	@	MTA0
41	101	65	A	MTA1
42	102	66	B	MTA2
43	103	67	C	MTA3
44	104	68	D	MTA4
45	105	69	E	MTA5
46	106	70	F	MTA6
47	107	71	G	MTA7
48	110	72	H	MTA8
49	111	73	I	MTA9
4A	112	74	J	MTA10
4B	113	75	K	MTA11
4C	114	76	L	MTA12
4D	115	77	M	MTA13
4E	116	78	N	MTA14
4F	117	79	O	MTA15
50	120	80	P	MTA16
51	121	81	Q	MTA17
52	122	82	R	MTA18
53	123	83	S	MTA19
54	124	84	T	MTA20
55	125	85	U	MTA21

Hex	Oct	Dec	ASCII	Msg
56	126	86	V	MTA22
57	127	87	W	MTA23
58	130	88	X	MTA24
59	131	89	Y	MTA25
5A	132	90	Z	MTA26
5B	133	91	[MTA27
5C	134	92	\	MTA28
5D	135	93]	MTA29
5E	136	94	^	MTA30
5F	137	95	_	UNT
60	140	96	`	MSA0,PPE
61	141	97	a	MSA1,PPE
62	142	98	b	MSA2,PPE
63	143	99	c	MSA3,PPE
64	144	100	d	MSA4,PPE
65	145	101	e	MSA5,PPE
66	146	102	f	MSA6,PPE
67	147	103	g	MSA7,PPE
68	150	104	h	MSA8,PPE
69	151	105	i	MSA9,PPE
6A	152	106	j	MSA10,PPE
6B	153	107	k	MSA11,PPE

Hex	Oct	Dec	ASCII	Msg
6C	154	108	l	MSA12,PPE
6D	155	109	m	MSA13,PPE
6E	156	110	n	MSA14,PPE
6F	157	111	o	MSA15,PPE
70	160	112	p	MSA16,PPD
71	161	113	q	MSA17,PPD
72	162	114	r	MSA18,PPD
73	163	115	s	MSA19,PPD
74	164	116	t	MSA20,PPD
75	165	117	u	MSA21,PPD
76	166	118	v	MSA22,PPD
77	167	119	w	MSA23,PPD
78	170	120	x	MSA24,PPD
79	171	121	y	MSA25,PPD
7A	172	122	z	MSA26,PPD
7B	173	123	{	MSA27,PPD
7C	174	124		MSA28,PPD
7D	175	125	}	MSA29,PPD
7E	176	126	~	MSA30,PPD
7F	177	127	DEL	

Message Definitions

Mnemonics	Definition
DCL	Device Clear
GET	Group Execute Trigger
GTL	Go To Local
LLO	Local Lockout
MLA	My Listen Address
MSA	My Secondary Address
MTA	My Talk Address
PPC	Parallel Poll Configure
PPD	Parallel Poll Disable
PPE	Parallel Poll Enable
PPU	Parallel Poll Unconfigure
SDC	Selected Device Clear
SPD	Serial Poll Disable
SPE	Serial Poll Enable
TCT	Take Control
UNL	Unlisten
UNT	Untalk



Operation of the GPIB

This appendix describes basic concepts you need to understand to operate the GPIB. It also contains a description of the physical and electrical characteristics of the GPIB and configuration requirements of the GPIB.

Types of Messages

The GPIB carries device-dependent messages and interface messages.

- Device-dependent messages, often called *data* or *data messages*, contain device-specific information such as programming instructions, measurement results, machine status, and data files.
- Interface messages manage the bus itself. They are usually called *commands* or *command messages*. Interface messages perform such tasks as initializing the bus, addressing and unaddressing devices, and setting device modes for remote or local programming.

Do not confuse the term *command* as used here with some device instructions, which can also be called commands. These device-specific instructions are actually data messages.

Talkers, Listeners, and Controllers

GPIB devices can be Talkers, Listeners, and/or Controllers. A digital voltmeter, for example, is a Talker and may be a Listener as well. A Talker sends data messages to one or more Listeners. The Controller manages the flow of information on the GPIB by sending commands to all devices.

The GPIB is like an ordinary computer bus, except that the computer has its circuit cards interconnected via a backplane bus, whereas the GPIB has stand-alone devices interconnected via a cable bus.

The role of the GPIB Controller is similar to the role of the CPU of a computer, but a better analogy is to the switching center of a city

telephone system. The switching center (Controller) monitors the communications network (GPIB). When the center (Controller) notices that a party (device) wants to make a call (send a data message), it connects the caller (Talker) to the receiver (Listener).

The Controller addresses a Talker and a Listener before the Talker can send its message to the Listener. After the Talker transmits the message, the Controller may unaddress both devices.

Some bus configurations do not require a Controller. For example, one device may always be a Talker (called a Talk-only device) and there may be one or more Listen-only devices.

A Controller is necessary when you must change the active or addressed Talker or Listener. A computer usually handles the Controller function.

With the GPIB board and its software, your personal computer plays all three roles:

- Controller—to manage the GPIB
- Talker—to send data
- Listener—to receive data

The Controller-In-Charge and System Controller

Although there can be multiple Controllers on the GPIB, only one Controller at a time is active or Controller-In-Charge (CIC). You can pass active control from the current CIC to an idle Controller. Only one device on the bus—the System Controller—can make itself the CIC. The GPIB board is usually the System Controller.

GPIB Signals and Lines

The interface system consists of 16 signal lines and 8 ground-return or shield-drain lines.

The 16 signal lines are divided into three groups:

- Eight data lines
- Three handshake lines
- Five interface management lines

Data Lines

The eight data lines, DIO1 through DIO8, carry both data and command messages. All commands and most data use the 7-bit ASCII or International Standards Organization (ISO) code set, in which case the eighth bit, DIO8, is unused or is used for parity.

Handshake Lines

Three lines asynchronously control the transfer of message bytes among devices. This process is called a three-wire interlocked handshake, and it guarantees that message bytes on the data lines are sent and received without transmission error.

NRFD (not ready for data)

NRFD indicates whether a device is ready to receive a message byte. All devices drive NRFD when they receive commands, and Listeners drive it when they receive data messages.

NDAC (not data accepted)

NDAC indicates whether a device has accepted a message byte. All devices drive NDAC when they receive commands, and Listeners drive it when they receive data messages.

DAV (data valid)

DAV tells whether the signals on the data lines are stable (valid) and whether devices can accept them safely. The Controller drives DAV when sending commands, and the Talker drives it when sending data messages.

Interface Management Lines

Five lines manage the flow of information across the interface.

ATN (attention)

The Controller drives ATN true when it uses the data lines to send commands and drives ATN false when a Talker can send data messages.

IFC (interface clear)

The System Controller drives the IFC line to initialize the bus and become CIC.

REN (remote enable)

The System Controller drives the REN line, which places devices in remote or local program mode.

SRQ (service request)

Any device can drive the SRQ line to asynchronously request service from the Controller.

EOI (end or identify)

The EOI line has two purposes. The Talker uses the EOI line to mark the end of a message string. The Controller uses the EOI line to tell devices to respond in a parallel poll.

Physical and Electrical Characteristics

You usually connect devices with a cable assembly consisting of a shielded 24-conductor cable which has both a plug and a receptacle connector at each end. With this design, you can link devices in either a linear or a star configuration, or a combination of the two.

The standard connector is the Amphenol or Cinch Series 57 *Microribbon* or *Amp Champ* type. You can use an adapter cable with a non-standard cable and/or connector for special interconnection applications.

The GPIB uses negative logic with standard transistor-transistor logic (TTL) level. When DAV is true, for example, it is a TTL low level (≤ 0.8 V), and when DAV is false, it is a TTL high level (≥ 2.0 V).

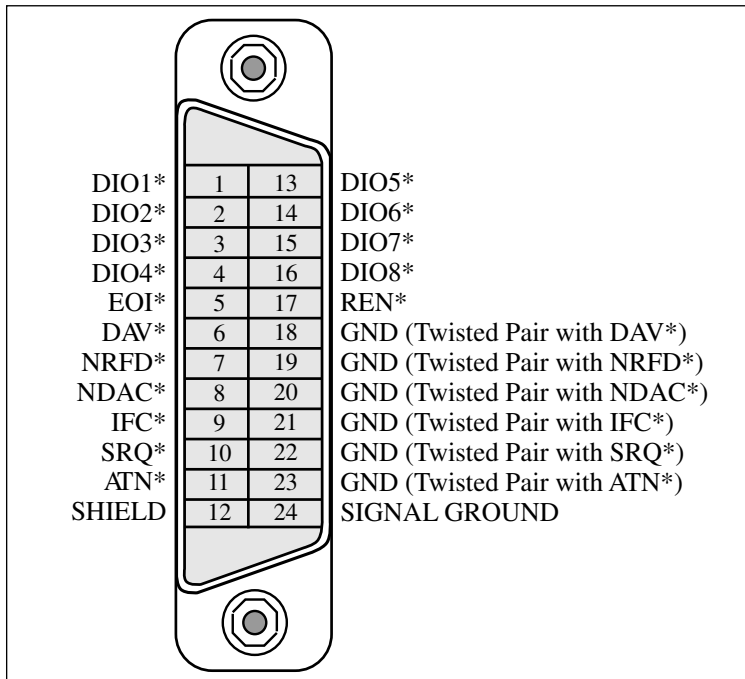


Figure B-1. GPIB Connector Showing Signal Assignment

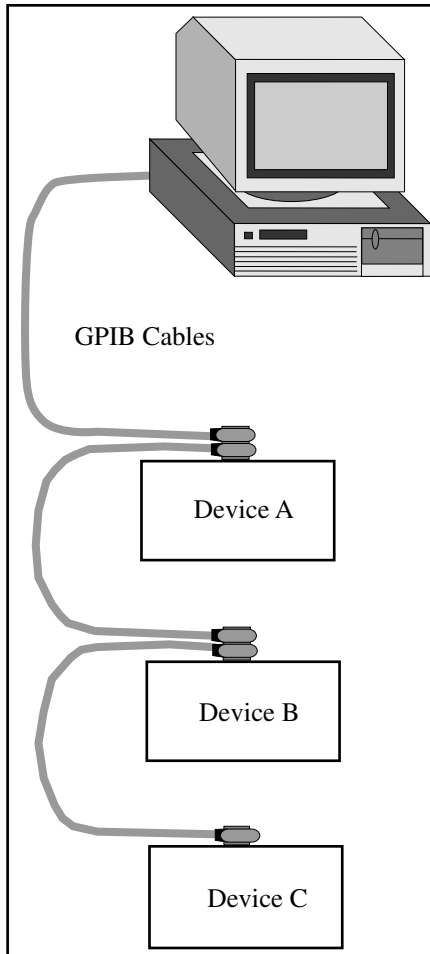


Figure B-2. Linear Configuration

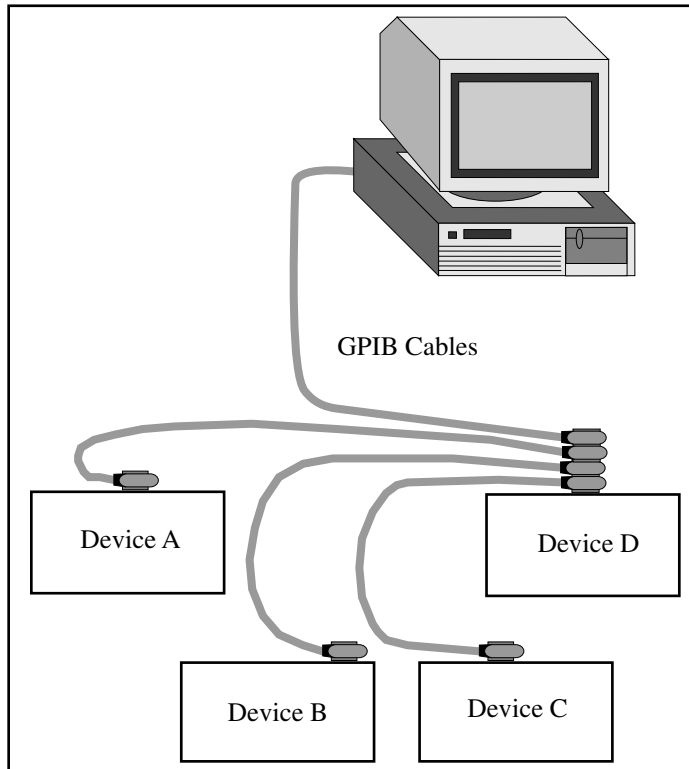


Figure B-3. Star Configuration

Configuration Requirements

To achieve the high data transfer rate for which the GPIB was designed, the physical distance between devices and the number of devices on the bus must be limited. The following restrictions are typical:

- A maximum separation of 4 m between any two devices and an average separation of 2 m over the entire bus.
- A maximum total cable length of 20 m.
- No more than 15 devices connected to each bus, with at least two-thirds powered on.

Contact National Instruments for bus extenders if your requirements exceed these limits.

Instrument I/O VI

Common Questions



This appendix answers common questions about LabVIEW and Instrument I/O. Questions are divided into GPIB and Serial I/O sections according to the relevant platform: All, Windows, Macintosh, Sun and HP-UX. Please contact National Instruments if you have further questions or suggestions regarding LabVIEW.

GPIB

Questions for All Platforms

When using a LabVIEW instrument driver, I have trouble talking with my instrument.

Ensure that your GPIB interface is working properly. Use the example `LabVIEW<->GPIB.vi` located in `examples\instr\smpIgpib.llb`. Try a simple command to the instrument; for example, if the instrument is 488.2, the string `*IDN?` will ask the instrument to send its identification string (about 100 characters).

Once communication is established, you should have no problems with the instrument driver.

I get timeout errors with GPIB Read/Write in LabVIEW.

Try running a simple program to establish communication between LabVIEW and GPIB. Use the example `LabVIEW<->GPIB.vi`, located in `examples\instr\smpIgpib.llb`. Try a simple command to the instrument; for example, if the instrument is 488.2, the string `*IDN?` will ask the instrument to send its identification string (about 100 characters).

If you still receive errors with GPIB, you may have a configuration problem. Open the GPIB configuration utility (Macintosh: `NI-488 Config`; Windows: `wibconf`; Sun: `ibconf`; HP-UX: `ibconf`). Verify that the settings match your hardware settings. Exit the

configuration utility and run the `ibic` (Interface Bus Interactive Control) utility for your platform (Macintosh: `ibic`, which ships with your GPIB board; Windows: `wibic`; Sun: `ibic`; HP-UX: `ibic`). Try the following sequence:

<code>: ibfind gpib0</code>	Find the GPIB interface
<code>id = 32000</code>	
<code>gpib0: ibsic</code>	Clear the GPIB bus (Send Interface Clear)
<code>[0130] [cmpl cic atn]</code>	Operation completed successfully
<code>gpib0: ibfind dev1</code>	Find device 1. Use appropriate instrument address.
<code>id = 32xxx</code>	
<code>dev1: ibwrt "*IDN?"</code>	Write string to instrument. 488.2 instruments recognize this command and return their identification string.
<code>count = 5</code>	Five bytes were sent.
<code>dev1: ibrd 100</code>	Read up to 100 bytes from the instrument.
<code>Fluke xxx Multimeter...</code>	Instrument returns identification string.

If you have any configuration errors, you will get an error message in one of these steps. The NI 488.2 Software Reference Manual that came with your GPIB board has detailed descriptions of the error messages.

Why can I communicate with my GPIB instrument with a LabVIEW VI running in execution highlighting mode but not when it is running full-speed?

This sounds like a timing problem. VIs run much slower with execution highlighting enabled than they will otherwise. Your instrument may need more time to prepare the data to send. Add a delay function or use service requests before the `GPIB Read.vi` to give the instrument enough time to generate the data it needs to send back to the computer.

Why can I write successfully to my GPIB instrument but can't read back from it?

When `GPIB Write.vi` executes, the computer is in talk mode and the instrument is in listen mode. When `GPIB Read.vi` executes, the device is supposed to switch to talk mode and the computer to listen mode. The device is prompted to switch modes by a termination signal which may be a character (End Of String) or a GPIB bus line (End Or Identify). So, if `GPIB Read.vi` times out or returns an EABO (Operation aborted) error, it means that the device is not receiving the right termination signal. To determine the termination mode for a given instrument, refer to its manual. As a rule of thumb, all IEEE 488.2 devices terminate on <CR><LF> and assertion of the EOI (End Or Identify) line in the GPIB bus.

Use the configuration utility for your platform (Macintosh: NI-488 Config; Windows: `wibconf`; Sun: `ibconf`; HP-UX: `ibconf`) to change the termination character.

A VI which talks with a GPIB instrument works fine on one platform and not on another.

Make sure that the instrument is configured properly in `wibconf`, `NI-488 Config`, or `ibconf`. Some older, 488.1 instruments do not automatically go to remote state. Go to the GPIB configuration utility for your platform and set the field "Assert REN when SC." This will ensure that the instrument goes to remote state (as opposed to local state) when it is addressed.

Windows Only

I can communicate with my instrument using a Quick Basic program but not from LabVIEW.

The GPIB board has separate handlers for DOS and Windows. Quick Basic accesses the DOS handler, but LabVIEW accesses the Windows handler. Make sure that the board and device are configured properly through `wibconf.exe`.

Version 3.0.1 does not surrender time to other applications in this situation and thus greatly improves the performance of GPIB.

Serial I/O

Questions for All Platforms

Why doesn't my instrument respond to commands sent with Serial Port Write.vi?

Many instruments expect a carriage return or line feed to terminate the command string. The `Serial Port Write.vi` in LabVIEW sends only those characters included in the string input; no termination character is appended to the string. Many terminal emulation (for example, Windows Terminal) packages automatically append a carriage return to the end of all transmissions. With LabVIEW, you will need to include the proper termination character in your string input to `Serial Port Write.vi` if your instrument requires it.

Some instruments require a carriage return (`\r`); others require a line feed (`\n`). When you enter a return on the keyboard (on PC keyboards, this is the Enter key on the main alphanumeric keypad), LabVIEW inserts a `\n` (line feed). To insert a carriage return, use `Concatenate Strings` and append a Carriage Return constant to the string, or manually enter (`\r`) after selecting **'\ Codes Display** from the string pop-up menu. To learn more about **'\ Codes Display**, see the *'\ Codes Display Option* section in Chapter 12, *String Controls and Indicators* of the *LabVIEW User Manual*.

Make sure that your cable works. Many of our technical support questions are related to bad cables. In computer to computer communication with serial I/O, use a null-modem to reverse the receive and transmit signals.

See the example `LabVIEW<->Serial.vi` to establish communication with your instrument. It is located in `examples\instr\smp1ser1.llb`. The VI also demonstrates the use of `Bytes at Serial Port.vi` before reading data back from the serial port.

How do I close a serial port so that other applications can use it?

You may wish to close the serial port after use. For example, on Windows, a VI may write information using `Serial Port`

`Write.vi` to `lpt1`, connected to a printer. After the operation is complete, LabVIEW still has control over the serial port. Other applications cannot use this port until LabVIEW has released control.

LabVIEW contains a `Close Serial Driver.vi` on all platforms. This VI tells LabVIEW to release control over the specified port. The `Close Serial Driver.vi` is not in the **Serial** palette; it is located in `vi.lib\Instr_sersup.llb`. To access the VI, use the **Functions»VI...** or **File»Open...** commands.

How do I reset or clear a serial port?

Use `Serial Port Init.vi` to reinitialize the port. This VI will automatically clear the serial port buffers associated with the port number.

How can I add additional serial ports to my computer?

You can add additional serial ports to your IBM-compatible PC using AT bus serial interface boards from National Instruments. (See Windows below). If you are using another platform, you can use a third party board to add serial ports to your computer. Some third party boards require a special language interface that does not conform to the standard API for serial ports on the platform. In this case, you will need to write your own interface, probably through a CIN or DLL, to the driver. The following sections explain how to use the serial port VIs included in LabVIEW to address boards which use the standard serial port interface on the different platforms.

Windows 3.x

LabVIEW for Windows uses the standard Microsoft Windows interface to the serial port. Thus, limitations in the use of serial ports in LabVIEW are due to limitations in Windows. For example, because Windows can only address ports COM1 through COM9, LabVIEW can only address these ports. Further, Windows allows access to only eight serial ports at any time; thus, LabVIEW can control a maximum of eight serial ports at any one time.

National Instruments now sells Plug and Play AT serial boards for a variety of solutions for serial communications. The AT-485 and AT-232 asynchronous serial interface boards are available in either 2 or 4-port configurations. Full Plug and Play compatibility gives you the benefits of switchless configuration for easier installation and maintenance. The AT-485 and AT-232 include the following software

components for use with Windows: device driver, hardware diagnostic test, and configuration utility. These boards have been tested with LabVIEW for Windows. Contact National Instruments for more information:

Manufacturer:	National Instruments
Product Name:	AT-485 and AT-232
U.S. Office:	Tel: 512-794-0100 Fax: 512-794-8411
FaxBack:	800-329-7177 , order number 1430

A third party board which uses some work-around to overcome Windows limitations will in all probability not work well, if at all, with LabVIEW. It is possible to write a CIN or DLL to access such boards, but this is not recommended. In general, boards that use the standard Windows interface should be completely compatible with LabVIEW.

Windows 95 and Windows NT

Unlike Windows 3.x, Windows 95 and Windows NT are not limited to only eight serial ports. Under Windows 95 and Windows NT, LabVIEW can address as many as 256 serial ports if you want. The default port number parameter is 0 for COM1, 1 for COM2, 2 for COM3 and so on.

The file `labview.ini` contains the LabVIEW configuration options. To set the devices which will be used by the serial port VIs, set the configuration option `labview.serialDevices` to the list of devices to be used. For example, to set up your devices in a way similar to Windows 3.x:

```
labview.serialDevices="COM1; COM2; COM3; COM4;  
COM5; COM6; COM7; COM8; COM9; COM10; LPT1; LPT2;  
LPT3; LPT4;"
```

The above should appear as a single line in your configuration .

Macintosh

LabVIEW uses standard system INITs to talk with the serial ports. By default, LabVIEW uses the drivers `.aIn` and `.aOut`, the modem port, as port 0, and drivers `.bIn` and `.bOut`, the printer port, for port 1. To

access additional ports, you must install additional boards with the accompanying INITs.

After the board and INIT(s) are installed, tell LabVIEW how to use the additional ports. The method used in LabVIEW 4.0 is slightly different to that used in previous versions.

Modify the global `serpOpen.vi` (located in `vi.lib\Instr_serpup.llb`) to accommodate the additional ports. All serial port VIs call `Open Serial Driver.vi`, which reads the appropriate input and/or output driver names from the `serpOpen.vi`. The front panel of `serpOpen.vi` contains two string arrays, named `input driver names` and `out driver names`.

When `Open Serial Driver.vi` is called, it uses the port number input as the index into the string arrays. Therefore, to allow LabVIEW to recognize additional serial ports, add additional string elements into the input and output driver names, then select **Make Current Values Default** and save the changes to `serpOpen.vi`.

Each serial port on a plug-in board has two names, one for input and one for output. The exact names and instructions for installing the drivers comes with the documentation for the board. Contact the board manufacturer if the instructions are missing or unclear.

The following boards work with LabVIEW for Macintosh:

Manufacturer:	Creative Solutions, Inc.
Product Name:	Hurdler HQS (4 ports) or HDS (2 ports)
Phone:	301-984-0262
Fax:	301-770-1675

Manufacturer:	Greensprings
Product Name:	RM 1280 (4 ports)
Phone:	415-327-1200
Fax:	415-327-3808

Sun and HP-UX

When you use the serial port VIs on a Sun SPARCstation under Solaris 1, the port number parameter is 0 for `/dev/ttya`, 1 for `/dev/ttyb`, 2 for `/dev/ttyc`, and so on. Under Solaris 2, port 0 refers to `/dev/cua/a`, 1 to `/dev/cua/b`, and so on.. Under HP-UX port number 0 refers to `/dev/tty00`, 1 to `/dev/tty01`, and so on.

Because other vendor's serial port boards can have arbitrary device names, LabVIEW has developed an easy interface to keep the numbering of ports simple. In LabVIEW for Sun version 3.x, a configuration option exists to tell LabVIEW how to address the serial ports. LabVIEW will support any board which uses standard UNIX devices. Some manufacturers suggest using `cua` rather than `tty` device nodes with their boards. LabVIEW for Sun can address both types of nodes.

The file `.labviewrc` (`.Xdefaults` in previous versions) contains the LabVIEW configuration options. To set the devices which will be used by the serial port VIs, set the configuration option `labview.serialDevices` to the list of devices to be used. For example, the default is:

```
labview.serialDevices:/dev/ttya:/dev/ttyb:/dev/
ttyc:...:/dev/ttyz
```

**Note:**

This requires that any third party serial board installation include a method of creating a standard `/dev` file (node) and that the user knows the name of that file.

The following boards should work with LabVIEW for Sun:

Manufacturer:	Sun
Product Name:	SBus Serial Parallel Controller (8 serial, 1 parallel)
Manufacturer:	Artecon, Inc.
Product Name:	SUNX-SB-300P ArtePort SBus Card with 3 Ser / 1 Par Ports SUNX-SB-400P ArtePort SBus Card with 4 Ser / 1 Par Ports SUNX-SB-1600 ArtePort SBus Card with 16 Serial Ports

For any of these products contact SunExpress at **800-873-7869**.

How can I control the DTR and RTS serial lines?

`Serial Port Init.vi` can be used to configure the serial port for hardware handshaking; however, some applications may require manual toggling of the DTR and RTS lines. Because the interface to the serial ports is platform-dependent, each platform has a separate mechanism to control the lines.

Windows

The LabVIEW for Windows distribution contains a VI which you can use to drive the DTR and RTS serial lines. The VI `serial_line_ctrl.vi`, located in `vi.lib\Instr_sersup.llb`, can be used to control these lines. The VI will toggle these lines according to the function input. Valid codes for the input are:

- 0 - noop
- 1 - clear DTR
- 2 - set DTR
- 3 - clear RTS
- 4 - set RTS
- 5 - set DTR protocol
- 6 - clr DTR protocol
- 7 - noop2

Macintosh

On the Macintosh, you can use the `Device Control/Status` function to control the serial port. *Inside Macintosh* (see Volume II, pages 245-259 and Volume. IV, pages 225-228), contains specific information on what `csCodes` can be sent to the serial port. A summary of the codes and their functions is listed here:

<u>code</u>	<u>param</u>	<u>Effect</u>
13	baudRate	Set baud rate (actual rate, as an integer)
14	serShk	Set handshake parameters
16	byte	Set miscellaneous control options
17		Asserts DTR
18		Negates DTR
19	char	Replace parity errors
20	2 chars	Replace parity errors with alternate character
21		Unconditionally set XOff for output flow control

22	Unconditionally clear XOff for output flow control
23	Send XOn for input flow control if XOff was sent last
24	Unconditionally send XOn for input flow control
25	Send XOff for input flow control if XOn was sent last
26	Unconditionally send XOff for input flow control
27	Reset SCC channel

Sun

LabVIEW for Sun contains no specific support to toggle the hardware handshaking lines of the serial ports. To manually control these lines, you must write a CIN. See the section *Steps for Creating a CIN* in Chapter 1, *CIN Overview*, of the *LabVIEW Code Interface Reference Manual* for further details on how to write a CIN.

Why can't I allocate a serial buffer larger than 32kbytes?

You can not use a buffer size on `Serial Port Init.vi` that is larger than 32kbytes because Windows and Macintosh limit the serial port buffer to 32k; thus, if you allocate a buffer larger than this, LabVIEW truncates the buffer size to 32k. This is not a problem on the Sun.

Windows Only

How do I access the parallel port?

In LabVIEW for Windows, port 10 is LPT1, port 11 is LPT2, and so on. To send data to a printer connected to a parallel port, use the `Serial Port Write.vi`. See Chapter 8, *Serial Port VIs*, for more details.

What do the error numbers received from the serial port VIs mean?

The Serial Port VIs in LabVIEW for Windows return the errors reported by the Windows `GetCommError` function. Error numbers returned by the Serial Port VIs are 0x4000 (16,384) 'Or'-ed with the error numbers in the following table. Note that the error returned reflects the status of the serial port; the error may have been generated as the result of a previous serial port function. The return values can be a combination of the following errors:

Hex Value	Error Name	Meaning
0x0001	CE_RXOVER	Receiving queue overflowed. There was either no room in the input queue or a character was received after the end-of-file character was received.
0x0002	CE_OVERRUN	Character was not read from the hardware before the next character arrived. The character was lost.
0x0004	CE_RXPARITY	Hardware detected a parity error.
0x0008	CE_FRAME	Hardware detected a framing error.
0x0010	CE_BREAK	Hardware detected a break condition.
0x0020	CE_CTSTO	CTS (clear-to-send) timeout. While a character was being transmitted, CTS was low for the duration specified by the <code>fCtsHold</code> member of <code>COMSTAT</code> .
0x0040	CE_DSRTTO	DSR (data-set-ready) timeout. While a character was being transmitted, DSR was low for the duration specified by the <code>fDsrHold</code> member of <code>COMSTAT</code> .
0x0080	CE_RLSDTO	RLSD (receive-line-signal- detect) timeout. While a character was being transmitted, RLSD was low for the duration specified by the <code>fRlsdHold</code> member of <code>COMSTAT</code> .

Hex Value	Error Name	Meaning
0x0100	CE_TXFULL	Transmission queue was full when a function attempted to queue a character.
0x0200	CE_PTO	Timeout occurred during an attempt to communicate with a parallel device.
0x0400	CE_IOE	I/O error occurred during an attempt to communicate with a parallel device.
0x0800	CE_DNS	Parallel device was not selected.
0x1000	CE_OOP	Parallel device signaled that it is out of paper.
0x8000	CE_MODE	Requested mode is not supported, or the <i>idComDev</i> parameter is invalid. If set, CE_MODE is the only valid error.

To use this table, take the error number and dissect it into its error components. For example, if `Serial Port Write.vi` returns the error 16,408, then the errors returned are CE_BREAK and CE_FRAME (16,408 = 16,384 + 16 + 8 = 0x4000 + 0x0010 + 0x0008).

Sun Only

I receive an error -37 when performing serial I/O.

Error -37 means that LabVIEW cannot find the appropriate serial device. This indicates that either a) the `/dev/tty?` files do not exist on your machine, or b) LabVIEW cannot find the file `serpdrv`.

By default, LabVIEW addresses `/dev/ttya` as port 0, `/dev/ttyb` as port 1, and so on. These devices must exist and the user must have read and write permissions to access the devices. You can change the devices LabVIEW accesses with the serial port VIs by adding the `serialDevices` configuration option to your `.Xdefaults` file. See the *LabVIEW Configuration Options* section in the *LabVIEW for Sun Release Notes* on how to use this option.

The `serpdrv` file is shipped with LabVIEW and serves as the interface between LabVIEW and the Sun serial ports. This file needs to be in the location specified by the `libdir` configuration option, set

to the LabVIEW directory by default. This means that `serpdrv` needs to be in the same directory as `gpibdrv` and `vi.lib`.

Serial I/O hangs on a Solaris 1.x machine.

LabVIEW for Sun uses asynchronous I/O calls when performing serial port operations. In the `Generic_Small` kernel, asynchronous I/O has been commented out. To access the serial ports from LabVIEW for Sun, the user must use the standard `Generic` kernel (not `Generic_Small`), or rebuild the `Generic_Small` kernel and reboot the SPARC.



Error Codes

This appendix lists the error codes returned by LabVIEW functions, including the name and description of each code.

VISA Error Codes

The following table lists the VISA error codes and their descriptions:

Table D-1. VISA Error Codes

Code	Name	Description
-1073807360	VI_ERROR_SYSTEM_ERROR	Unknown system error (miscellaneous error).
-1073807346	VI_ERROR_INV_OBJECT VI_ERROR_INV_SESSION	The given session or object reference is invalid.
-1073807344	VI_ERROR_INV_EXPR	Invalid expression specified for search.
-1073807343	VI_ERROR_RSRC_NFOUND	Insufficient location information or the device or resource is not present in the system.
-1073807342	VI_ERROR_INV_RSRC_NAME	Invalid resource reference specified. Parsing error.
-1073807341	VI_ERROR_INV_ACC_MODE	Invalid access mode.
-1073807339	VI_ERROR_TMO	Timeout expired before operation completed.
-1073807338	VI_ERROR_CLOSING_FAILED	Unable to deallocate the previously allocated data structures corresponding to this session or object reference.

Table D-1. VISA Error Codes (Continued)

Code	Name	Description
-1073807331	VI_ERROR_NSUP_ATTR	The specified attribute is not defined or supported by the referenced resource.
-1073807330	VI_ERROR_NSUP_ATTR_STATE	The specified state of the attribute is not valid, or is not supported as defined by the resource.
-1073807329	VI_ERROR_ATTR_READONLY	The specified attribute is read-only.
-1073807322	VI_ERROR_INV_EVENT	Specified event type is not supported by the resource.
-1073807321	VI_ERROR_INV_MECH	Invalid mechanism specified.
-1073807320	VI_ERROR_HNDLR_NINSTALLED	A handler was not installed.
-1073807319	VI_ERROR_INV_HNDLR_REF	The given handler reference is invalid.
-1073807318	VI_ERROR_INV_CONTEXT	Specified event context is invalid.
-1073807308	VI_ERROR_RAW_WR_PROT_VIOL	Violation of raw write protocol occurred during transfer.
-1073807307	VI_ERROR_RAW_RD_PROT_VIOL	Violation of raw read protocol occurred during transfer.
-1073807306	VI_ERROR_OUTP_PROT_VIOL	Device reported an output protocol error during transfer.
-1073807305	VI_ERROR_INP_PROT_VIOL	Device reported an input protocol error during transfer.
-1073807304	VI_ERROR_BERR	Bus error occurred during transfer.
-1073807302	VI_ERROR_INV_SETUP	Unable to start operation because setup is invalid (due to attributes being set to an inconsistent state).

Table D-1. VISA Error Codes (Continued)

Code	Name	Description
-1073807300	VI_ERROR_ALLOC	Insufficient system resources to perform necessary memory allocation.
-1073807299	VI_ERROR_INV_MASK	Invalid buffer mask specified.
-1073807298	VI_ERROR_IO	Could not perform read/write operation because of I/O error.
-1073807297	VI_ERROR_INV_FMT	A format specifier in the format string is invalid.
-1073807295	VI_ERROR_NSUP_FMT	A format specifier in the format string is not supported.
-1073807294	VI_ERROR_LINE_IN_USE	The specified trigger line is currently in use.
-1073807286	VI_ERROR_SRQ_NOCCURRED	Service request has not been received for the session.
-1073807282	VI_ERROR_INV_SPACE	Invalid address space specified.
-1073807279	VI_ERROR_INV_OFFSET	Invalid offset specified.
-1073807276	VI_ERROR_NSUP_OFFSET	Specified offset is not accessible from this hardware.
-1073807273	VI_ERROR_WINDOW_NMAPPED	The specified session is not currently mapped.
-1073807265	VI_ERROR_NLISTENERS	No Listeners condition is detected (both NRFD and NDAC are deasserted).
-1073807264	VI_ERROR_NCIC	The interface associated with this session is not currently the controller in charge.
-1073807257	VI_ERROR_NSUP_OPER	The given session or object reference does not support this operation.

Table D-1. VISA Error Codes (Continued)

Code	Name	Description
-1073807242	VI_ERROR_NSUP_WIDTH	Specified width is not supported by this hardware.
-1073807239	VI_ERROR_INV_PROT	The protocol specified is invalid.
-1073807237	VI_ERROR_INV_SIZE	Invalid size of window specified.
-1073807232	VI_ERROR_WINDOW_MAPPED	The specified session already contains a mapped window.
-1073807231	VI_ERROR_NIMPL_OPER	The given operation is not implemented.

Instrument Driver Error Codes

The following table lists the error codes reserved for LabVIEW instrument drivers. Error codes -1200 to -1299 are reserved for instrument driver specific I/O error. -1300 and above are reserved for instrument specific error codes. All instrument specific error codes should be represented with a corresponding error message in the PREFIX Error Message VI.

Table D-2. Instrument Driver Error Codes

Status	Description	Set By
0	No error: the call was successful	
-1200	Invalid syntax string	VISA Transition Library
-1201	Error finding instruments	VISA Transition Library
-1202	Unable to initialize interface or instrument	VISA Transition Library
-1205	Invalid Instrument Handle	VISA Transition Library
-1210	Parameter out of range	Instrument Driver VI
-1215	Error closing instrument	VISA Transition Library
-1218	Error getting attribute	VISA Transition Library

Table D-2. Instrument Driver Error Codes (Continued)

Status	Description	Set By
-1219	Error setting attribute	VISA Transition Library
-1223	Instrument identification query failed	Instrument Driver Initialize VI
-1224	Error clearing instrument	VISA Transition Library
-1225	Error triggering instrument	VISA Transition Library
-1226	Error polling instrument	VISA Transition Library
-1230	Error writing to instrument	VISA Transition Library
-1231	Error reading from instrument	VISA Transition Library
-1236	Error interpreting instrument response	Instrument Driver VI
-1240	Instrument timed out	VISA Transition Library
-1250	Error mapping VXI address	VISA Transition Library
-1251	Error unmapping VXI address	VISA Transition Library
-1252	Error accessing VXI address	VISA Transition Library
-1260	Function not supported by controller	VISA Transition Library
-1300	Instrument-specific error	
-13xx	Developer specified error codes	

GPIB Error Codes

The following table lists the GPIB error codes and their descriptions:

Table D-3. GPIB Error Codes

Code	Name	Description
0	EDVR	Error connecting to driver.
1	ECIC	Command requires GPIB Controller to be CIC.

Table D-3. GPIB Error Codes (Continued)

Code	Name	Description
2	ENOL	Write detected no Listeners.
3	EADR	GPIB Controller not addressed correctly.
4	EARG	Invalid argument or arguments.
5	ESAC	Command requires GPIB Controller to be SC.
6	EABO	I/O operation aborted.
7	ENEB	Non-existent board.
8	EDMA	DMA Hardware error detected.
9	EBTO	DMA hardware μ P bus timeout.
11	ECAP	No capability.
12	EFSO	File system operation error.
13	EOWN	Shareable board exclusively owned.
14	EBUS	GPIB bus error.
15	ESTB	Serial poll byte queue overflow.
16	ESRQ	SRQ stuck on.
17	ECMD	Unrecognized command.
19	EBNP	Board not present.
20	ETAB	Table error.
30	NADDR	No GPIB address input.
31	NSTRG	No string input (write).
32	NCNT	No count input (read).


Serial Port Error Codes

The following table lists the Serial Port error codes and their descriptions:

Table D-4. Serial Port Error Codes

Code	Name	Description
61	EPAR	Serial port parity error.
62	EORN	Serial port overrun error.
63	EOFL	Serial port receive buffer overflow.
64	EFRM	Serial port framing error.
65	SPTMO	Serial port timeout, bytes not received at serial port.

Customer Communication

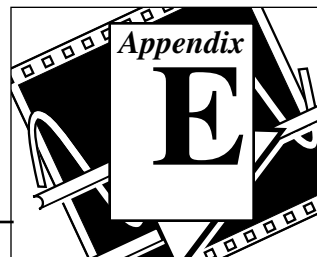

For your convenience, this appendix contains forms to help you gather the information necessary to help us solve your technical problems and a form you can use to comment on the product documentation. When you contact us, we need the information on the Technical Support Form and the configuration form, if your manual contains one, about your system configuration to answer your questions as quickly as possible.

National Instruments has technical assistance through electronic, fax, and telephone systems to quickly provide the information you need. Our electronic services include a bulletin board service, an FTP site, a FaxBack system, and e-mail support. If you have a hardware or software problem, first try the electronic support systems. If the information available on these systems does not answer your questions, we offer fax and telephone support through our technical support centers, which are staffed by applications engineers.

Electronic Services



Bulletin Board Support

National Instruments has BBS and FTP sites dedicated for 24-hour support with a collection of files and documents to answer most common customer questions. From these sites, you can also download the latest instrument drivers, updates, and example programs. For recorded instructions on how to use the bulletin board and FTP services and for BBS automated information, call (512) 795-6990. You can access these services at:

United States: (512) 794-5422 or (800) 327-3077

Up to 14,400 baud, 8 data bits, 1 stop bit, no parity

United Kingdom: 01635 551422

Up to 9,600 baud, 8 data bits, 1 stop bit, no parity

France: 1 48 65 15 59

Up to 9,600 baud, 8 data bits, 1 stop bit, no parity



FTP Support

To access our FTP site, log on to our Internet host, [ftp.natinst.com](ftp://ftp.natinst.com), as anonymous and use your Internet address, such as joesmith@anywhere.com, as your password. The support files and documents are located in the /support directories.



FaxBack Support

FaxBack is a 24-hour information retrieval system containing a library of documents on a wide range of technical information. You can access FaxBack from a touch-tone telephone at the following numbers:

(512) 418-1111 or (800) 329-7177



E-Mail Support (currently U.S. only)

You can submit technical support questions to the appropriate applications engineering team through e-mail at the Internet addresses listed below. Remember to include your name, address, and phone number so we can contact you with solutions and suggestions.

GPIB: gpib.support@natinst.com

DAQ: daq.support@natinst.com

VXI: vxi.support@natinst.com

LabWindows: lw.support@natinst.com

LabVIEW: lv.support@natinst.com

HiQ: hiq.support@natinst.com

VISA: visa.support@natinst.com

Fax and Telephone Support

National Instruments has branch offices all over the world. Use the list below to find the technical support number for your country. If there is no National Instruments office in your country, contact the source from which you purchased your software to obtain support.



Fax



Telephone

Australia	03 9 879 9422	03 9 879 9179
Austria	0662 45 79 90 0	0662 45 79 90 19
Belgium	02 757 00 20	02 757 03 11
Canada (Ontario)	519 622 9310	519 622 9311
Canada (Quebec)	514 694 8521	514 694 4399
Denmark	45 76 26 00	45 76 71 11
Finland	90 527 2321	90 502 2930
France	1 48 14 24 24	1 48 14 24 14
Germany	089 741 31 30	089 714 60 35
Hong Kong	2645 3186	2686 8505
Italy	02 48301892	02 48301915
Japan	03 5472 2970	03 5472 2977
Korea	02 596 7456	02 596 7455
Mexico	95 800 010 0793	5 520 3282
Netherlands	0348 433466	0348 430673
Norway	32 84 84 00	32 84 86 00
Singapore	2265886	2265887
Spain	91 640 0085	91 640 0533
Sweden	08 730 49 70	08 730 43 70
Switzerland	056 200 51 51	056 200 51 55
Taiwan	02 377 1200	02 737 4644
U.K.	01635 523545	01635 523154

Technical Support Form

Photocopy this form and update it each time you make changes to your software or hardware, and use the completed copy of this form as a reference for your current configuration. Completing this form accurately before contacting National Instruments for technical support helps our applications engineers answer your questions more efficiently.

If you are using any National Instruments hardware or software products related to this problem, include the configuration forms from their user manuals. Include additional pages if necessary.

Name _____

Company _____

Address _____

Fax (____) _____ Phone (____) _____

Computer brand _____ Model _____ Processor _____

Operating system (include version number) _____

Clock speed _____MHz RAM _____MB Display adapter _____

Mouse ____yes ____no Other adapters installed _____

Hard disk capacity _____MB Brand _____

Instruments used _____

National Instruments hardware product model _____ Revision _____

Configuration _____

National Instruments software product _____ Version _____

Configuration _____

The problem is: _____

List any error messages: _____

The following steps reproduce the problem: _____

Documentation Comment Form

National Instruments encourages you to comment on the documentation supplied with our products. This information helps us provide quality products to meet your needs.

Title: *LabVIEW Instrument I/O VI Reference Manual*

Edition Date: November 1995

Part Number: 320537C-01

Please comment on the completeness, clarity, and organization of the manual.

If you find errors in the manual, please record the page numbers and describe the errors.

Thank you for your help.

Name

Title


Company

Address

Phone ()

Mail to: Technical Publications
National Instruments Corporation
6504 Bridge Point Parkway, MS 53-02
Austin, TX 78730-5039

Fax to: Technical Publications
National Instruments Corporation
(512) 794-5678

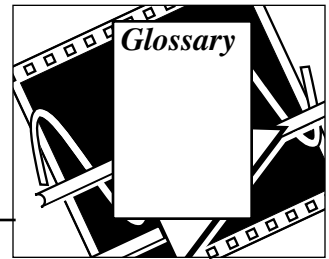

Prefix	Meaning	Value
n-	nano-	10^{-9}
μ -	micro-	10^{-6}
m-	milli-	10^{-3}
k-	kilo-	10^3
M-	mega-	10^6

A

ANSI	American National Standards Institute.
array	Ordered, indexed set of data elements of the same type.
ASCII	American Standard Code for Information Interchange.
asynchronous	Mode in which multiple processes share processor time; one executes
execution	while the others wait for interrupts, such as while performing device I/O or waiting for a clock tick.
ATN	The GPIB attention line that distinguishes between commands and data messages. When ATN is asserted, bytes on the GPIB DIO lines are commands.

C

CIC	<i>See</i> Controller-In-Charge.
-----	----------------------------------

commands or command messages	GPIB message that performs tasks such as initializing the bus, addressing and unaddressing devices, and setting device modes.
Controller or Controller-In-Charge	The device that manages the GPIB by sending interface messages to other devices. The GPIB itself can also act as Controller-In-Charge.
CPU	Central processing unit.
CR	Carriage return.
CTS	Clear To Send (line).

D

data or data messages	Device-dependent messages that contain programming instructions, measurement results, machine status, and so on.
DAV	Data valid. One of the three GPIB handshake lines. <i>See</i> handshake.
device	An instrument, peripheral, computer, or other electronics equipment that can be programmed over the GPIB.
device function	A function that operates on or otherwise pertains to a GPIB device rather than to the GPIB board in the computer.
DIO	Digital Input/Output (line).
DIO1 through DIO8	The GPIB lines that transmit command or data bytes from one device to another.
DMA	Direct memory access. High-speed data transfer between the GPIB and memory that is not handled directly by the CPU. DMA is not available on some systems.
driver	Common term for software used to manipulate a device or interface board.
DSR	Data Set Ready (line).
DTR	Data Terminal Ready (line).

E

EABO	An error code that means the I/O operation aborted.
ECIC	An error code that means the function requires the GPIB board to be Controller-In-Charge.
EOI	End or identify. A GPIB line that is used to signal either the last byte of a data message (END) or carry the parallel poll identify (IDY) message.
EOS or EOS byte	A 7- or 8-bit End-of-String character that is sent as the last byte of a data message.
ESAC	An error code which means the GPIB board is not System Controller as required.
ETAB	An error code which means there is a table problem.

G

General Purpose Interface Bus *See* GPIB.

GPIB	General Purpose Interface Bus is the common name for the communications interface system defined in ANSI/IEEE Standard 488.1-1987 and ANSI/IEEE Standard 488.2-1987. Hewlett-Packard, the inventor of the bus, calls it the HP-IB.
GPIB address	The address of a device on the GPIB, composed of a primary address (MLA and MTA) and perhaps a secondary address (MSA). The GPIB board has both a GPIB address and an I/O address.
GPIB board	Reference to the National Instruments family of GPIB boards.
GTL	Go to local. The GPIB command used to place an addressed Listener in local (front panel) control mode.

H

handshake	The mechanism used to transfer bytes from the source handshake function of one device to the acceptor handshake function of another device. The three GPIB lines DAV, NRFD, and NDAC
-----------	--

are used in an interlocked fashion to signal the phases of the transfer, so that bytes can be sent asynchronously (for example, without a clock) at the speed of the slowest device.

hex Hexadecimal.

high-level function A device function that combines several rudimentary board operations into one function so that the user does not have to be concerned with bus management or other GPIB protocol matters. *See* low-level function.

I

IDY Mnemonic for *identify*, a remote GPIB message.

IEEE Institute of Electrical and Electronic Engineers.

IFC Mnemonic for *interface clear*, a remote GPIB message.

interface message A broadcast message sent from the Controller to all devices and used to manage the GPIB. Common interface messages include Interface Clear, listen addresses, talk addresses, and Serial Poll Enable/Disable. *See* data or data messages.

I/O Input/output. In the context of this manual, the transmission of commands or messages between the computer via the GPIB board and other devices on the GPIB.

ISO International Standards Organization.

ist An individual status bit of the status byte used in the parallel poll configure function.

L

LabVIEW Laboratory Virtual Instrument Engineering Workbench.

LF Linefeed.

Listener A GPIB device that receives data messages from a Talker.

LLO Local lockout. The GPIB command used to tell all devices that they can or should ignore remote (GPIB) data messages or local (front panel) controls, depending on whether the device is in local or remote program mode.

LOK Lockout bit.

low-level function A rudimentary board or device function that performs a single operation.

M

m Meters.

MB Megabytes of memory.

MLA My Listen Address.

MSA My Secondary Address.

MTA My Talk Address.

N

NaN Digital display value for a floating-point representation of *not a number*, typically the result of an undefined operation, such as $\log(-1)$.

NDAC Not data accepted. One of the three GPIB handshake lines. *See* handshake.

NL New line or linefeed.

NRFD Not ready for data. One of the three GPIB handshake lines. *See* handshake.

O

oct Octal.

P

parallel poll configure	The process of polling all configured devices at once and reading a composite poll response. <i>See</i> serial poll.
parallel poll configure	<i>See</i> PPC.
parallel poll disable	<i>See</i> PPD.
parallel poll enable	<i>See</i> PPE.
parallel poll unconfigure	<i>See</i> PPU.
PPC	The GPIB command that configures an addressed Listener to participate in polls.
PPD	The GPIB command that disables a configured device from participating in polls. There are 16 PPD commands.
PPE	The GPIB command that enables a configured device to participate in polls and to assign a DIO response line. There are 16 PPE commands.
PPU	The GPIB command that disables any device from participating in polls.

R

REN	Remote enable. A GPIB line controlled by the System Controller but used by the CIC to place devices in remote program mode.
RTS	Request To Send (line).

S

SCPI	Standard Commands for Programmable Instruments.
SCSI	Small Computer System Interface (bus).
sec	Seconds.
serial poll	The process of polling and reading the status byte of one device at a time. <i>See</i> parallel poll.

SRQ Service request. The GPIB line that a device asserts to notify the CIC that the device needs servicing.

status byte The data byte sent by a device when it is serially polled.

System Controller The single designated Controller that can assert control (become CIC of the GPIB) by sending the Interface Clear (IFC) message. Other devices can become CIC only by having control passed to them.

T

Talker A GPIB device that sends data messages to Listeners.

timeout A feature of LabVIEW that prevents I/O functions from hanging indefinitely when there is a problem on the GPIB.

TTL Transistor-transistor logic.

U

UNL Unlisten. The GPIB command that unaddresses any active Listeners.

UNT Untalk. The GPIB command that unaddresses an active Talker.


V

VI Virtual instrument. LabVIEW program; so called because it models the appearance and function of a physical instrument.

VISA Virtual Instrument Software Architecture. A single interface library for controlling VXI, GPIB, RS-232, and other types of instruments.

W

wire Data path between nodes.

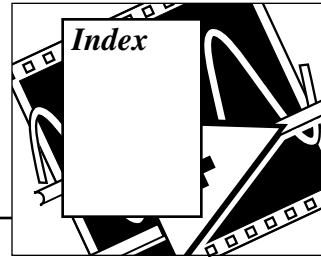

A

action VIs, for instrument drivers, 2-7
address parameter, GPIB 488.2 functions
 specifying, 1-13
address parameter, GPIB functions
 changing primary and secondary
 addresses, 7-2
 default primary address, 7-2
 primary address, 7-2
 secondary address, 7-2
address string parameter, traditional GPIB
 functions, 6-2
AllSPoll function, 7-9 to 7-10
ANSI/IEEE Standard 488.1-1987, 1-9
application instrument driver VIs
 definition, 2-6
 example, 2-6
 instrument driver development
 considerations, 3-24
 purpose and use, 2-6
Assert or unassert remote enable. *See* sre
 (Unassert or assert remote enable function).
ATN (attention) line, B-3

B

block diagram
 for instrument drivers, 3-19
bulletin board access for instrument driver, 1-2
bus extenders for GPIB devices, B-7
bus management functions
 FindLstn, 1-13, 7-16 to 7-17
 list of functions, 1-15 to 1-16
 overview, 1-15
 ResetSys, 7-17 to 7-18
 SendIFC, 7-18

SendLLO, 7-18 to 7-19
SetRWLS, 7-19 to 7-20
TestSRQ, 7-20
TestSys, 7-21
WaitSRQ, 7-22

bus parameters, GPIB 488.2 functions, 7-2
byte count parameters, GPIB 488.2
 functions, 7-2
Bytes at Serial Port VI, 8-7

C

cac (Become active Controller) function,
 6-19
CD-ROM distribution of instrument
 drivers, 1-2
CIC (Controller-In-Charge), B-2
Close VI
 PREFIX Close VI, 4-6 to 4-7
 required for instrument drivers, 2-7
cmd (Send IEEE 488 commands)
 function, 6-20
command messages, GPIB, B-1
command strings
 functions (table), 6-6 to 6-7
compatible hardware. *See* hardware
 compatible with GPIB.
component instrument drivers VIs
 action/status, 2-7
 close VI, 2-7
 configuration VIs, 2-7
 data VIs, 2-7
 developer-specified VIs, 2-8
 initialize VI, 2-7
 template VIs, 2-8
 utility VIs, 2-7

- Conduct parallel poll function. *See* rpp (Conduct parallel poll) function.
- configuration utility, for changing GPIB function parameters, 7-2
- configuration, GPIB
 - linear configuration (illustration), B-6
 - requirements, B-7
 - star configuration (illustration), B-7
- connector pane, for instrument drivers
 - rules for creating, 3-20 to 3-21
 - selecting pattern with extra terminals (note), 3-21
- connector, GPIB, B-5
- Controller functions
 - cac (Become active Controller), 6-19
 - cmd (Send IEEE 488 commands), 6-20
 - dma (Set DMA mode or programmed I/O mode), 6-20
 - gts (Go from active Controller to standby), 6-20 to 6-21
 - ist (Set individual status bit), 6-21
 - llo (Local lockout), 6-22
 - loc (Place Controller in local state), 6-22
 - off (Take controller offline), 6-22
 - ppc (Parallel poll configure, enable and disable), 6-22 to 6-23
 - ppu (Parallel poll unconfigure), 6-23
 - rpp (Conduct parallel poll), 6-23
 - rsc (Release or request system control), 6-23
 - rsv (Request service and/or set the serial poll status byte), 6-24
 - sic (Send interface clear), 6-24
 - sre (Unassert or assert remote enable), 6-25
- Controller-In-Charge (CIC), B-2
- Controllers
 - address string when using multiple Controllers, 6-2
 - compatibility of National Instruments GPIB boards, 1-10 to 1-12
 - controlling IEEE 488.1 and IEEE 488.2 devices, 1-14

- forcing GPIB board to be Controller, 7-18
- GPIB error codes (table), 6-11 to 6-12
- IEEE Standard 488.2, 1-9
- LLO (Local Lockout)
 - message, 7-18 to 7-19
- operation of GPIB Controller, B-2
- PPE (Parallel Poll Enable)
 - message, 6-21
- required capabilities, 1-10
- status bits (table), 6-10
- controls
 - descriptions for online help, 3-22
 - front panel of instrument drivers
 - guidelines for placing, 3-17 to 3-18
 - required controls, 3-16 to 3-17
- customer communication, xix, E-1

D

- data acquisition boards
 - compatibility with LabVIEW software
 - LabVIEW for HP-UX, 1-11
 - LabVIEW for
 - Macintosh, 1-11 to 1-12
 - LabVIEW for Sun, 1-11
 - LabVIEW for Windows 3.x, 1-10
 - LabVIEW for Windows NT, 1-11
- data dependency
 - in instrument drivers, 3-13
- data lines, GPIB, B-3
- data messages, GPIB, B-1
- data types
 - icons for, *xvii* to *xviii*
- data VIs, for instrument drivers, 2-7
- DAV (data valid) line, B-3
- DCL (Device Clear) message, 6-4, 6-24, 7-17
- Design, 3-4
- DevClear function, 7-3
- DevClearList function, 7-10
- developer-specified VIs, for instrument drivers, 2-8
- Device Clear message. *See* DCL (Device Clear) message.

device functions

- command strings functions
 - (table), 6-6 to 6-7
- loc (Go to local), 6-17
- off (Take device offline), 6-18
- pc (Pass control), 6-18
- ppc (Parallel poll configure), 6-18 to 6-19

dma (Set DMA mode or programmed I/O mode) function, 6-20

documentation

- conventions, *xvi* to *xviii*
- organization of manual, *xv* to *xvi*
- related documentation, *xviii*

DTR serial lines

- Macintosh, C-9 to C-10
- Sun, C-10
- Windows, C-9

E

electrical characteristics of GPIB, B-4 to B-7

EnableLocal function, 7-11

EnableRemote function, 7-11 to 7-12

EOI (end or identify) line, B-4

error codes

- GPIB (table), D-5 to D-6
- instrument driver (table), D-4 to D-5
- serial port VIs, 8-3
- serial ports (table), D-7
- VISA (table), D-1 to D-4

error handler VIs

- connecting error codes with, 8-7

error message VI

- PREFIX Error Message VI, 4-7 to 4-8

error parameters

- traditional GPIB functions, 6-2

error query VI

- PREFIX Error Query VI, 4-9 to 4-10

error reporting

- instrument driver front panel, 3-17
- instrument driver template VIs, 4-5
- instrument drivers, 2-9

errors

- ESAC error, 6-25

- GPIB error codes (table), 6-11 to 6-12
- obtaining error information, 1-13

ESAC error, 6-25

F

Fast Data Channel Mode, 5-27

Fast Data Channel Number, 5-27

Fast Data Channel Pair, 5-27

Fast Data Channel Signal Enable, 5-28

FindLstn function

- description, 7-16 to 7-17
- implementation of Find All Listeners protocol, 1-13

FindRQS function, 7-12 to 7-13

Format & Append function, 3-12

front panel

- for instrument drivers
 - control guidelines, 3-17 to 3-18
 - general requirements, 3-16
 - required controls, 3-16 to 3-17

functional body of instrument drivers

- application VIs, 2-5
- component VIs, 2-5
- external interface model, 2-2
- internal interface model, 2-4 to 2-5

Functions palettepalette

- GPIB 488.2, 7-1
- GPIB Traditional option, 6-1
- Instrument I/O option, 6-1, 7-1, 8-1

functions. *See* Controller functions; device functions.

G

general functions

- list of functions, 1-17
- overview, 1-17
- MakeAddr, 7-27
- SetTimeout, 7-28

General Purpose Interface Bus. *See* GPIB.

GET (Group Execute Trigger) message, 6-12

Go from active Controller to standby function. *See* gts (Go from active Controller

- to standby)
- Go To Local function. *See* loc (Go to local) device function.
- Go To Local message. *See* GTL (Go To Local) message.
- GPIO
 - ANSI/IEEE Standard 488.1-1987, 1-9
 - command messages, B-1
 - common questions, C-1 to C-3
 - compatible GPIO hardware
 - LabVIEW for
 - Macintosh, 1-11 to 1-12
 - LabVIEW for Sun, 1-11
 - LabVIEW for Windows 3.x, 1-10
 - configuration
 - linear configuration
 - (illustration), B-6
 - requirements, B-7
 - star configuration (illustration), B-7
 - Controller-In-Charge, B-2
 - Controllers, B-2 to B-4
 - data lines, B-3
 - data messages, B-1
 - definition, 1-9
 - electrical characteristics, B-4 to B-7
 - error codes, D-5 to D-6
 - GPIO connector signal assignment, B-5
 - handshake lines, B-3
 - history, 1-9
 - IEEE 488.1 standard, 1-9
 - IEEE 488.2 standard, 1-9
 - IEEE Standard 488-1975, 1-9
 - interface management lines, B-3 to B-4
 - Listeners, B-1 to B-2
 - physical characteristics, B-4 to B-7
 - questions and answers about
 - instrument I/O, C-1 to C-3
 - System Controller, B-2
 - Talkers, B-1 to B-2
 - types of messages, B-1
- GPIO 488.2 functions
 - bus management functions
 - FindLstn, 7-16 to 7-17
 - list of functions, 1-15 to 1-16
 - ResetSys, 7-17 to 7-18
 - SendIFC, 7-18
 - SendLLO, 7-18 to 7-19
 - SetRWLS, 7-19 to 7-20
 - TestSRQ, 7-20
 - TestSys, 7-21
 - WaitSRQ, 7-22
 - compared with traditional GPIO functions, 1-12
 - enhancements and additions, 1-12
 - example in `smplgpib.llb`, 7-1
 - general functions
 - list of functions, 1-17
 - MakeAddr, 7-27
 - SetTimeout, 7-28
 - getting help, 7-1
 - low-level I/O functions
 - list of functions, 1-16
 - RcvRespMsg, 7-23
 - ReceiveSetup, 7-24
 - SendCmds, 7-25
 - SendDataBytes, 7-25 to 7-26
 - SendSetup, 7-26 to 7-27
 - multiple-device functions
 - AllSPoll, 7-9 to 7-10
 - DevClearList, 7-10
 - EnableLocal, 7-11
 - EnableRemote, 7-11 to 7-12
 - FindRQS, 7-12 to 7-13
 - list of functions, 1-15
 - PPoll, 7-13
 - PPollUnconfig, 7-14
 - SendList, 7-14 to 7-15
 - TriggerList, 7-15 to 7-16
 - overview, 1-9 to 1-10
 - parameter descriptions, 7-2
 - single-device functions
 - DevClear, 7-3
 - PassControl, 7-5
 - PPollConfig, 7-4
 - ReadStatus, 7-5 to 7-6
 - Receive, 7-6 to 7-7

- Send, 7-7 to 7-8
- Trigger, 7-8 to 7-9
- GPIO 488.2 option, Functions palette, 7-1
- GPIO 488.2 palette, 7-1
- GPIO 488.2 VIs
 - bus management function
 - FindLstn, 1-13
- GPIO Clear function, 6-4
- GPIO error parameter, GPIO 488.2
 - functions, 7-2
- GPIO Initialization function
 - changing primary address, 7-2
 - description, 6-4 to 6-5
- GPIO Misc Function
 - changes from previous versions, 1-13
- GPIO Misc function
 - description, 6-5 to 6-7
- GPIO option, Function palette
 - GPIO 488.2 palette, 7-1
- GPIO Primary Address, 5-28
- GPIO Read function
 - description, 6-7 to 6-8
- GPIO Secondary Address, 5-28
- GPIO Serial Poll function, 6-9
- GPIO Status function
 - description, 6-9 to 6-12
 - no longer used, 1-13
- GPIO Traditional Functions palette, 6-1
- GPIO Trigger function, 6-12 to 6-13
- GPIO Wait function
 - description, 6-13 to 6-14
- GPIO Write function
 - description, 6-16
 - leaving devices in addressed state, 6-3
- Group Executive Trigger message. *See* GET (Group Executive Trigger) message.
- GTL (Go To Local) message, 7-11
- gts (Go from active Controller to standby)
 - function, 6-20 to 6-21

H

- handshake lines, GPIO
 - DAV (data valid), B-3

- NDAC (not data accepted), B-3
- NRFD (not ready for data), B-3
- handshaking
 - handshaking modes, serial port
 - VIs, 8-2 to 8-3
 - software handshaking
 - (XON/XOFF), 8-3
- hardware, 1-9
- hardware compatible with GPIO
 - LabVIEW for HP-UX, 1-11
 - LabVIEW for Macintosh, 1-11 to 1-12
 - LabVIEW for Sun, 1-11
 - LabVIEW for Windows NT, 1-11
- Hewlett-Packard's HP-IB, 1-9

I

- I/O interfaces of instrument drivers, 2-3
- icon, for instrument drivers, 3-20
- IDY (Identify) message, 6-19
- IFC (interface clear) line, B-3
- IFC (Interface Clear) message, 6-23, 6-24, 7-17, 7-18
- Immediate Servant, 5-29
- Increment Destination Count, 5-29
- Increment Source Count, 5-29
- indicators
 - descriptions for online help, 3-22
- initialize VI
 - PREFIX Initialize VI, 4-10 to 4-12
 - required for instrument drivers, 2-7
- instrument driver development
 - adding component VI VIs, 3-9 to 3-10
 - application VIs, 3-23
 - assembling command
 - strings, 3-11 to 3-12
 - block diagram guidelines, 3-19
 - checklist, 3-24 to 3-27
 - connector pane guidelines, 3-20 to 3-21
 - data dependency considerations, 3-13
 - design of structure, 3-1 to 3-7
 - example, 3-6 to 3-7
 - guidelines and recommendations, 3-4 to 3-6

- organization, 3-3 to 3-4
 - structure and VI hierarchy, 3-2
 - error reporting, 3-21
 - front panel guidelines
 - control guidelines, 3-17 to 3-18
 - general requirements, 3-16
 - required controls, 3-16 to 3-17
 - icon guidelines, 3-20
 - loop termination conditions, 3-10 to 3-11
 - modifying templates, 3-8 to 3-9
 - online help information, 3-21 to 3-23
 - control and indicator descriptions, 3-22
 - VI descriptions, 3-22
 - programmatically checking parameter ranges, 3-13 to 3-15
 - instrument driver template VIs
 - advantages, 1-3
 - close VI, 2-7
 - initialize VI, 2-7
 - lists of available templates, 1-3
 - modifying, 3-8 to 3-9
 - overview, 1-3
 - PREFIX Close, 4-6 to 4-7
 - PREFIX Error Message, 4-7 to 4-8
 - PREFIX Error Query, 4-9 to 4-10
 - PREFIX Initialize, 4-10 to 4-12
 - PREFIX Message-Based Template, 4-13 to 4-14
 - PREFIX Register-Based Template, 4-14 to 4-16
 - PREFIX Reset, 4-16 to 4-17
 - PREFIX Revision Query, 4-17 to 4-19
 - PREFIX Self-Test, 4-19 to 4-20
 - purpose and use, 2-7, 4-1
 - instrument drivers
 - definition, 1-1, 2-1
 - distribution methods, 1-1 to 1-2
 - bulletin board access, 1-2
 - CD-ROM distribution, 1-2
 - Internet access, 1-2
 - error codes (table), D-4 to D-5
 - external interface model
 - functional body, 2-2
 - I/O interface, 2-3
 - illustration, 2-2
 - Interactive Developer Interface, 2-3
 - programmatically interface, 2-3
 - subroutine interface, 2-4
 - internal design model
 - application VIs, 2-6
 - component VIs, 2-6 to 2-8
 - error reporting, 2-9
 - functional body, 2-4 to 2-5
 - illustration, 2-5
 - libraries of drivers, 1-1
 - overview, 1-1 to 1-3
 - Instrument I/O option, 6-1, 7-1, 8-1
 - instrument I/O, questions about
 - GPIO, C-1 to C-3
 - serial I/O, C-4 to C-13
 - Interactive Developer Interface for instrument drivers, 2-3
 - Interface Clear (IFC) message. *See* IFC (Interface Clear) message.
 - interface management lines, GPIO
 - ATN (attention), B-3
 - EOI (end or identify), B-4
 - IFC (interface clear), B-3
 - REN (remote enable), B-4
 - SRQ (service request), B-4
 - Interface Number, 5-30
 - Interface Type, 5-30
 - Internet access to instrument drivers, 1-2
 - IO Protocol, 5-28
 - ist (Set individual status bit) function, 6-21
- ## L
- LabVIEW GPIO 488.2 VIs. *See* GPIO 488.2 VIs.
 - LabVIEW software
 - compatibility with data acquisition boards
 - LabVIEW for HP-UX, 1-11
 - LabVIEW for Macintosh, 1-11
 - LabVIEW for Sun, 1-11

- LabVIEW for Windows 3.x, 1-10
- LabVIEW for Windows NT, 1-11
- LabVIEW-compatible GPIB hardware. *See* hardware compatible with GPIB.
- lines. *See* signals and lines, GPIB.
- Listeners, B-1 to B-2
- llo (Local lockout) function, 6-22
- LLO (Local Lockout) message, 7-18 to 7-19
- loc (Place Controller in local state) function, 6-22
- Local lockout function. *See* llo (Local lockout) function.
- Local Lockout message. *See* LLO (Local Lockout) message.
- Lockout State, setting, 7-19
- loop termination method, including in instrument drivers, 3-10 to 3-11
- low-level I/O functions
 - list of functions, 1-16
 - overview, 1-16
 - RcvRespMsg, 7-23
 - ReceiveSetup, 7-24
 - SendCmds, 7-25
 - SendDataBytes, 7-25 to 7-26
 - SendSetup, 7-26 to 7-27
- LPE (Local Poll Enable) message, 6-22

M

- Mainframe Logical Address, 5-30
- MakeAddr function
 - changing primary and secondary address formats, 7-2
 - description, 7-27
- manual. *See* documentation.
- Manufacturer ID, 5-30
- Maximum Queue Length, 5-31
- message-based template VI
 - PREFIX Message-Based Template VI, 4-13 to 4-14
- mnemonics for multiline interface messages, A-8
- Model Code, 5-31
- multiline interface messages, A-2 to A-7

- definitions, A-8
- mnemonics, A-8
- multiple-device functions
 - AllSPoll, 7-9 to 7-10
 - DevClearList, 7-10
 - EnableLocal, 7-11
 - EnableRemote, 7-11 to 7-12
 - FindRQS, 7-12 to 7-13
 - list of functions, 1-15
 - overview, 1-15
 - PPoll, 7-13
 - PPollUnconfig, 7-14
 - SendList, 7-14 to 7-15
 - TriggerList, 7-15 to 7-16

N

- NDAC (not data accepted) line, B-3
- NRFD (not ready for data) line, B-3

O

- off (Take controller offline) function, 6-22
- online help information for instrument drivers
 - control and indicator descriptions, 3-22
 - VI descriptions, 3-22

P

- Parallel Poll Enable (PPE)
 - message, 6-19, 6-21, 6-23
- parallel poll functions
 - ppc (Parallel poll configure) Controller function, 6-22 to 6-23
 - ppc (Parallel poll configure) device function, 6-18 to 6-19
 - ppu (Parallel poll unconfigure) function, 6-23
 - rpp (Conduct parallel poll) function, 6-23
- parameter range checking, for instrument drivers, 3-13 to 3-15
- Pass control function. *See* pct (Pass control)

function.
 PassControl function, 7-5
 Pick Line & Append function, 3-12
 Place Controller in local state function. *See* loc
 (Place Controller in local state) function.
 port number parameter, serial port VIs, 8-2
 ppc (Parallel poll configure, enable and
 disable) function, 6-22 to 6-23
 PPE (Parallel Poll Enable)
 message, 6-19, 6-21, 6-23
 PPoll function, 7-13
 PPollConfig function, 7-4
 PPollUnconfig function, 7-14
 ppu (Parallel poll unconfigure) function, 6-23
 PREFIX Close VI, 4-6 to 4-7
 PREFIX Error Message VI, 4-7 to 4-8
 PREFIX Error Query VI, 4-9 to 4-10,
 4-10 to 4-12
 PREFIX Message-Based Template VI,
 4-13 to 4-14
 PREFIX Register-Based Template VI,
 4-14 to 4-16
 PREFIX Reset VI, 4-16 to 4-17
 PREFIX Revision Query VI, 4-17 to 4-19
 PREFIX Self-Test VI, 4-19 to 4-20
 primary address, GPIB. *See* address
 parameter, GPIB 488.2 VIs.

Q

questions and answers about I/O
 GPIB, C-1 to C-3
 serial I/O, C-4 to C-13

R

RcvRespMsg function, 7-23
 ReadStatus function, 7-5 to 7-6
 Receive function, 7-6 to 7-7
 ReceiveSetup function, 7-24
 register-based template VI
 PREFIX Register-Based Template VI,
 4-14 to 4-16
 Release or request system control function.

See rsc (Release or request system control)
 function.
 Remote Enable (REN) message. *See* REN
 (Remote Enable) message.
 Remote With Lockout State, 7-19
 REN (remote enable) line, B-4
 REN (Remote Enable) message
 releasing or requesting with rsc
 function, 6-23
 sending with SetRWLS function, 7-19
 unasserting or asserting, 6-25, 7-11, 7-17
 Request service and/or set the serial poll
 status byte. *See* rsv (Request service and/or
 set the serial poll status byte) function.
 Reset (RST) message, 7-17
 reset VI
 PREFIX Reset VI, 4-16 to 4-17
 ResetSys function, 7-17 to 7-18
 Resource Lock State, 5-31
 Resource Manufacturer Identification, 5-31
 Resource Manufacturer Name, 5-32
 Resource Name, 5-32 to 5-33
 responsibility for addressing, 6-2
 revision query VI
 PREFIX Revision Query VI,
 4-17 to 4-19
 rpp (Conduct parallel poll) function, 6-23
 rsc (Release or request system control)
 function, 6-23
 rsv (Request service and/or set the serial poll
 status byte) function, 6-24
 RTS serial lines
 Macintosh, C-9 to C-10
 Sun, C-10
 Windows, C-9

S

SCPI (Standard Commands for
 Programmable Instruments), 1-10
 SDC (Selected Device Clear)
 message, 6-4, 6-24
 secondary address, GPIB. *See* address
 parameter, GPIB 488.2 VIs.

- Selected Device Clear (SDC) message. *See* SDC (Selected Device Clear) message.
- self-test VI
 - PREFIX Self-Test VI, 4-19 to 4-20
- Send End Enable, 5-33
- Send function, 7-7 to 7-8
- Send IEEE 488 commands functions. *See* cmd (Send IEEE 488 commands) function.
- Send interface clear function. *See* sic (Send interface clear) function.
- SendCmds function, 7-25
- SendDataBytes function, 7-25 to 7-26
- SendIFC function, 7-18
- SendList function, 7-14 to 7-15
- SendLLO function, 7-18 to 7-19
- SendSetup function, 7-26 to 7-27
- serial I/O, questions about, C-4 to C-13
- Serial option, 8-1
- Serial palette, 8-1
- serial poll functions
 - GPIOB Serial Poll VI, 6-9
- Serial Port
 - error codes, D-7
- serial port boards, accessing
 - HP-UX, C-7 to C-8
 - Macintosh, C-6 to C-7
 - Windows 3.x, C-5 to C-6
 - Windows 95 and Windows NT, C-6
- serial port buffer
 - handshaking modes, 8-2 to 8-3
 - software handshaking (XON/XOFF), 8-3
- serial port VIs
 - Bytes at Serial Port, 8-3 to 8-4, 8-7
 - common parameters, 8-2
 - examples in `smplsrl.llb`, 8-1
 - getting help, 8-1
 - handshaking modes, 8-2 to 8-3
 - list of VIs, 1-17
 - overview, 1-17
 - port number parameter, 8-2
 - Serial Port Break, 8-4 to 8-5
 - Serial Port Init, 8-5 to 8-7
 - Serial Port Write, 8-7 to 8-8
 - software handshaking (XON/XOFF), 8-3
 - Set DMA mode or programmed I/O mode function. *See* dma (Set DMA mode or programmed I/O mode) function.
 - Set individual status bit function. *See* ist (set individual status bit) function.
 - SetRWLS function, 7-19 to 7-20
 - SetTimeout function
 - changing timeout value, 7-2
 - description, 7-28
 - sic (Send interface clear) function, 6-24
 - signal and lines, GPIB
 - data lines, B-3
 - handshake lines
 - DAV (data valid), B-3
 - NDAC (not data accepted), B-3
 - NRFD (not ready for data), B-3
 - interface management lines
 - ATN (attention), B-3
 - EOI (end or identify), B-4
 - IFC (interface clear), B-3
 - REN (remote enable), B-4
 - SRQ (service request), B-4
 - single-device functions
 - list of functions, 1-14
 - overview, 1-14
 - DevClear, 7-3
 - PassControl, 7-5
 - PPollConfig, 7-4
 - ReadStatus, 7-5 to 7-6
 - Receive, 7-6 to 7-7
 - Send, 7-7 to 7-8
 - Trigger, 7-8 to 7-9
 - Slot, 5-33
 - software handshaking (XON/XOFF), 8-3
 - sre (Unassert or assert remote enable) function, 6-25
 - SRQ, B-4
 - asserting, 7-12
 - description, B-4
 - determining current state, 7-20
 - requesting with `rsv` function, 6-24

standards

- ANSI/IEEE Standard 488.1-1987, 1-9
- IEEE 488.1 standard, 1-9
- IEEE Standard 488.2, 1-9
- IEEE Standard 488-1975, 1-9
- Standard Commands for Programmable Instruments (SCPI), 1-10

status bits

- ist (Set individual status bit) function, 6-21
- table, 6-10

status parameters

- GPIO 488.2 functions, 7-2
- traditional functions, 6-2
- traditional GPIB functions, 6-2

status VIs, for instrument drivers, 2-7

subroutine interface for instrument drivers, 2-4

Suppress End Enable, 5-34

System Controller. *See also* Controllers.

- designation of, 7-2
- GPIO Controller as, 6-21
- making Controller-in-Charge, 7-17
- operation of, B-2

T

Talkers, B-1 to B-2

template VIs. *See* instrument driver

template VIs.

Termination Character, 5-34

Termination Character Enable, 5-34, 5-34

TestSRQ function, 7-20

TestSys function, 7-21

Timeout Value, 5-34

timeout value

- changing with SetTimeout function, 7-2, 7-28

traditional GPIB functions

- behavior, 6-3
- compared with GPIO 488.2 functions, 1-13
- examples in `smplgpiplib`, 6-1
- getting help, 6-1
- GPIO Clear function, 6-4

GPIO Initialization, 6-4 to 6-5

GPIO Misc, 1-13, 6-5 to 6-7

GPIO Read, 6-7 to 6-8

GPIO Serial Poll, 6-9

GPIO Status, 6-9 to 6-12

GPIO Trigger, 6-12 to 6-13

GPIO Wait, 6-13 to 6-14

GPIO Write, 6-16

lists of functions, 1-12 to 1-13

overview, 1-12

parameter description, 6-2

Wait for GPIO RQS, 6-15

Trigger function, 7-8 to 7-9

Trigger Identifier, 5-35

TriggerList function, 7-15 to 7-16

troubleshooting. *See* questions and answers about instrument I/O.**U**Unassert or assert remote enable. *See* `Unassert or assert remote enable` function.

User Data, 5-35

utility VIs, for instrument drivers, 2-7

V

Version of Implementation, 5-36

Version of Specification, 5-36

VI descriptions, 3-22

Virtual Instrument Software Architecture (VISA), 1-4

VISA Assert Trigger function, 5-3 to 5-4

VISA Attribute Node, 5-26 to 5-38

description, 1-6

list of attributes, 1-6 to 1-8

VISA Attribute Descriptions

Fast Data Channel Mode, 5-27

Fast Data Channel Number, 5-27

Fast Data Channel Pair, 5-27

Fast Data Channel Signal

Enable, 5-28

GPIO Primary Address, 5-28

GPIO Secondary Address, 5-28

- Immediate Servant, 5-29
- Increment Destination Count, 5-29
- Increment Source Count, 5-29
- Interface Number, 5-30
- Interface Type, 5-30
- IO Protocol, 5-28
- Mainframe Logical Address, 5-30
- Manufacturer ID, 5-30
- Maximum Queue Length, 5-31
- Model Code, 5-31
- Resource Lock State, 5-31
- Resource Manufacturer
 - Identification, 5-31
- Resource Manufacturer Name, 5-32
- Resource Name, 5-32 to 5-33
- Send End Enable, 5-33
- Slot, 5-33
- Suppress End Enable, 5-34
- Termination Character, 5-34
- Termination Character Enable, 5-34, 5-34
- Timeout Value, 5-34
- Trigger Identifier, 5-35
- User Data, 5-35
- Version of Implementation, 5-36
- Version of Specification, 5-36
- VXI Commander Logical
 - Address, 5-36
- VXI Logical Address, 5-37
- VXI Memory Address Space, 5-37
- VXI Memory Base Address, 5-37
- VXI Memory Size, 5-37
- Window Access, 5-38
- Window Base Address, 5-38
- Window Size, 5-38
- VISA Clear function, 5-4
- VISA Close function, 5-4 to 5-5
- VISA Disable Event function, 5-12 to 5-13
- VISA Discard Events function, 5-13 to 5-14
- VISA Enable Event function, 5-14
- VISA Find Resource function, 5-5 to 5-6
- VISA function
 - list of functions, 1-4 to 1-6
- VISA functions overview, 1-4
- VISA In8 / In16 / In32 function, 5-16 to 5-17
- VISA Library Function, 5-1
- VISA Library palette, 5-1
- VISA Library Reference
 - definition, 1-4
 - overview, 1-4
- VISA Lock function, 5-6
- VISA Map Address function, 5-23, 5-25
- VISA Memory Allocation function, 5-18, 5-23
- VISA Memory Free function, 5-18, 5-24
- VISA Move In8 / Move In16 / Move In32
 - function, 5-19 to 5-20
- VISA Move Out8 / Move Out16 / Move Out32
 - function, 5-20 to 5-21
- VISA Open function, 5-7 to 5-8
- VISA Operation Descriptions
 - Event Handling functions
 - VISA Disable Event function, 5-12 to 5-13
 - VISA Discard Events function, 5-13 to 5-14
 - VISA Enable Event function, 5-14
 - VISA WaitOnEvent function, 5-15 to 5-16
 - High Level Register Access functions
 - VISA In8 / In16 / In32 function, 5-16 to 5-17
 - VISA Memory Allocation function, 5-18
 - VISA Memory Free function, 5-18
 - VISA Move In8 / Move In16 / Move In32
 - function, 5-19 to 5-20
 - VISA Move Out8 / Move Out16 / Move Out32
 - function, 5-20 to 5-21
 - VISA Out8 / Out16 / Out32
 - function, 5-21 to 5-22
 - Low Level Register Access functions
 - VISA Map Address function, 5-23, 5-25

- VISA Memory Allocation
 - function, 5-23
- VISA Memory Free function, 5-24
- VISA Peek8 / Peek16 / Peek32
 - function, 5-24
- VISA Poke 8 / Poke 16 / Poke 32
 - function, 5-25
- VISA Assert Trigger function, 5-3 to 5-4
- VISA Clear function, 5-4
- VISA Close function, 5-4 to 5-5
- VISA Find Resource function, 5-5 to 5-6
- VISA Lock function, 5-6
- VISA Open function, 5-7 to 5-8
- VISA Read function, 5-8 to 5-9
- VISA Read STB function, 5-9 to 5-10
- VISA Status Description function, 5-10
- VISA Unlock function, 5-11
- VISA Write function, 5-11 to 5-12
- VISA Operation Parameters, 5-2 to 5-3
- VISA Out8 / Out16 / Out32 function,
 - 5-21 to 5-22
- VISA Peek8 / Peek16 / Peek32 function, 5-24
- VISA Poke8 / Poke16 / Poke32 function, 5-25
- VISA Read function, 5-8 to 5-9
- VISA Read STB function, 5-9 to 5-10
- VISA session input control, required for
 - instrument driver front panel, 3-17
- VISA session output control, required for
 - instrument driver front panel, 3-17
- VISA Status Description function, 5-10
- VISA Unlock function, 5-11
- VISA WaitOnEvent function, 5-15 to 5-16
- VISA Write function, 5-11 to 5-12
- VXI Commander Logical Address, 5-36
- VXI Logical Address, 5-37
- VXI Memory Address Space, 5-37
- VXI Memory Base Address, 5-37
- VXI Memory Size, 5-37
- Window Access, 5-38
- Window Base Address, 5-38
- Window Size, 5-38

W

- Wait for GPIB RQS, 6-15
- wait state vector bits (table), 6-14
- WaitSRQ function, 7-22