COMPREHENSIVE SERVICES

We offer competitive repair and calibration services, as well as easily accessible documentation and free downloadable resources.

SELL YOUR SURPLUS

We buy new, used, decommissioned, and surplus parts from every NI series. We work out the best solution to suit your individual needs.

Sell For Cash Get Credit Receive a Trade-In Deal

OBSOLETE NI HARDWARE IN STOCK & READY TO SHIP

We stock New, New Surplus, Refurbished, and Reconditioned NI Hardware.



Bridging the gap between the manufacturer and your legacy test system.

0

1-800-915-6216



www.apexwaves.com

sales@apexwaves.com

All trademarks, brands, and brand names are the property of their respective owners.

Request a Quote



PXIe-4309

SC Express

PXIe-4309 User Manual

32 Ch (8 ADC), 2 MS/s, 18 - 28 bit, Flexible Resolution PXI Analog Input Module



Worldwide Technical Support and Product Information

Worldwide Offices

Visit ni.com/niglobal to access the branch office websites, which provide up-to-date contact information, support phone numbers, email addresses, and current events.

National Instruments Corporate Headquarters

11500 North Mopac Expressway Austin, Texas 78759-3504 USA

Tel: 866 ASK MYNI (275 6964)

For further support information, refer to the *NI Services* appendix. To comment on National Instruments documentation, refer to the National Instruments website at ni.com/info and enter the Info Code feedback.

© 2017 National Instruments. All rights reserved.

Legal Information

Limited Warranty

This document is provided 'as is' and is subject to being changed, without notice, in future editions. For the latest version, refer to ni.com/manuals. NI reviews this document carefully for technical accuracy; however, NI MAKES NO EXPRESS OR IMPLIED WARRANTIES AS TO THE ACCURACY OF THE INFORMATION CONTAINED HEREIN AND SHALL NOT BE LIABLE FOR ANY ERRORS.

NI warrants that its hardware products will be free of defects in materials and workmanship that cause the product to fail to substantially conform to the applicable NI published specifications for one (1) year from the date of invoice.

For a period of ninety (90) days from the date of invoice, NI warrants that (i) its software products will perform substantially in accordance with the applicable documentation provided with the software and (ii) the software media will be free from defects in materials and workmanship.

If NI receives notice of a defect or non-conformance during the applicable warranty period, NI will, in its discretion: (i) repair or replace the affected product, or (ii) refund the fees paid for the affected product. Repaired or replaced Hardware will be warranted for the remainder of the original warranty period or ninety (90) days, whichever is longer. If NI elects to repair or replace the product, NI may use new or refurbished parts or products that are equivalent to new in performance and reliability and are at least functionally equivalent to the original part or product.

You must obtain an RMA number from NI before returning any product to NI. NI reserves the right to charge a fee for examining and testing Hardware not covered by the Limited Warranty.

This Limited Warranty does not apply if the defect of the product resulted from improper or inadequate maintenance, installation, repair, or calibration (performed by a party other than NI); unauthorized modification; improper environment; use of an improper hardware or software key; improper use or operation outside of the specification for the product; improper voltages; accident, abuse, or neglect; or a hazard such as lightning, flood, or other act of nature.

THE REMEDIES SET FORTH ABOVE ARE EXCLUSIVE AND THE CUSTOMER'S SOLE REMEDIES, AND SHALL APPLY EVEN IF SUCH REMEDIES FAIL OF THEIR ESSENTIAL PURPOSE.

EXCEPT AS EXPRESSLY SET FORTH HEREIN, PRODUCTS ARE PROVIDED "AS IS" WITHOUT WARRANTY OF ANY KIND AND NI DISCLAIMS ALL WARRANTIES, EXPRESSED OR IMPLIED, WITH RESPECT TO THE PRODUCTS, INCLUDING ANY IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, TITLE OR NON-INFRINGEMENT, AND ANY WARRANTIES THAT MAY ARISE FROM USAGE OF TRADE OR COURSE OF DEALING, NI DOES NOT WARRANT, GUARANTEE, OR MAKE ANY REPRESENTATIONS REGARDING THE USE OF OR THE RESULTS OF THE USE OF THE PRODUCTS IN TERMS OF CORRECTNESS, ACCURACY, RELIABILITY, OR OTHERWISE, NI DOES NOT WARRANT THAT THE OPERATION OF THE PRODUCTS WILL BE UNINTERRUPTED OR ERROR FREE.

In the event that you and NI have a separate signed written agreement with warranty terms covering the products, then the warranty terms in the separate agreement shall control.

Copyright

Under the copyright laws, this publication may not be reproduced or transmitted in any form, electronic or mechanical, including photocopying, recording, storing in an information retrieval system, or translating, in whole or in part, without the prior written consent of National Instruments Corporation.

National Instruments respects the intellectual property of others, and we ask our users to do the same. NI software is protected by copyright and other intellectual property laws. Where NI software may be used to reproduce software or other materials belonging to others, you may use NI software only to reproduce materials that you may reproduce in accordance with the terms of any applicable license or other legal restriction.

End-User License Agreements and Third-Party Legal Notices

You can find end-user license agreements (EULAs) and third-party legal notices in the following locations:

- Notices are located in the <National Instruments>_Legal Information and <National Instruments> directories.
- EULAs are located in the <National Instruments>\Shared\MDF\Legal\license directory.
- Review <National Instruments>_Legal Information.txt for information on including legal information in
 installers built with NI products.

U.S. Government Restricted Rights

If you are an agency, department, or other entity of the United States Government ("Government"), the use, duplication, reproduction, release, modification, disclosure or transfer of the technical data included in this manual is governed by the Restricted Rights provisions under Federal Acquisition Regulation 52.227-14 for civilian agencies and Defense Federal Acquisition Regulation Supplement Section 252.227-7014 and 252.227-7015 for military agencies.

Trademarks

Refer to the NI Trademarks and Logo Guidelines at ni.com/trademarks for more information on National Instruments trademarks.

ARM, Keil, and µVision are trademarks or registered of ARM Ltd or its subsidiaries.

LEGO, the LEGO logo, WEDO, and MINDSTORMS are trademarks of the LEGO Group.

TETRIX by Pitsco is a trademark of Pitsco, Inc.

FIELDBUS FOUNDATION™ and FOUNDATION™ are trademarks of the Fieldbus Foundation.

EtherCAT® is a registered trademark of and licensed by Beckhoff Automation GmbH.

CANopen® is a registered Community Trademark of CAN in Automation e.V.

DeviceNet[™] and EtherNet/IP[™] are trademarks of ODVA.

Go!, SensorDAQ, and Vernier are registered trademarks of Vernier Software & Technology. Vernier Software & Technology and vernier.com are trademarks or trade dress.

Xilinx is the registered trademark of Xilinx, Inc.

Taptite and Trilobular are registered trademarks of Research Engineering & Manufacturing Inc.

FireWire® is the registered trademark of Apple Inc.

Linux® is the registered trademark of Linus Torvalds in the U.S. and other countries.

Handle Graphics[®], MATLAB[®], Real-Time Workshop[®], Simulink[®], Stateflow[®], and xPC TargetBox[®] are registered trademarks, and TargetBox[™] and Target Language Compiler[™] are trademarks of The MathWorks, Inc.

Tektronix®, Tek, and Tektronix, Enabling Technology are registered trademarks of Tektronix, Inc.

The Bluetooth® word mark is a registered trademark owned by the Bluetooth SIG, Inc.

The ExpressCard™ word mark and logos are owned by PCMCIA and any use of such marks by National Instruments is under license

The mark LabWindows is used under a license from Microsoft Corporation. Windows is a registered trademark of Microsoft Corporation in the United States and other countries.

Other product and company names mentioned herein are trademarks or trade names of their respective companies.

Members of the National Instruments Alliance Partner Program are business entities independent from National Instruments and have no agency, partnership, or joint-venture relationship with National Instruments.

Patents

For patents covering National Instruments products/technology, refer to the appropriate location: **Help»Patents** in your software, the patents.txt file on your media, or the *National Instruments Patent Notice* at ni.com/patents.

Export Compliance Information

Refer to the Export Compliance Information at ni.com/legal/export-compliance for the National Instruments global trade compliance policy and how to obtain relevant HTS codes, ECCNs, and other import/export data.

WARNING REGARDING USE OF NATIONAL INSTRUMENTS PRODUCTS

YOU ARE ULTIMATELY RESPONSIBLE FOR VERIFYING AND VALIDATING THE SUITABILITY AND RELIABILITY OF THE PRODUCTS WHENEVER THE PRODUCTS ARE INCORPORATED IN YOUR SYSTEM OR APPLICATION, INCLUDING THE APPROPRIATE DESIGN, PROCESS, AND SAFETY LEVEL OF SUCH SYSTEM OR APPLICATION.

PRODUCTS ARE NOT DESIGNED, MANUFACTURED, OR TESTED FOR USE IN LIFE OR SAFETY CRITICAL SYSTEMS, HAZARDOUS ENVIRONMENTS OR ANY OTHER ENVIRONMENTS REQUIRING FAIL-SAFE PERFORMANCE, INCLUDING IN THE OPERATION OF NUCLEAR FACILITIES; AIRCRAFT NAVIGATION; AIR TRAFFIC CONTROL SYSTEMS; LIFE SAVING OR LIFE SUSTAINING SYSTEMS OR SUCH OTHER MEDICAL DEVICES; OR ANY OTHER APPLICATION IN WHICH THE FAILURE OF THE PRODUCT OR SERVICE COULD LEAD TO DEATH, PERSONAL INJURY, SEVERE PROPERTY DAMAGE OR ENVIRONMENTAL HARM (COLLECTIVELY, "HIGH-RISK USES"). FURTHER, PRUDENT STEPS MUST BE TAKEN TO PROTECT AGAINST FAILURES, INCLUDING PROVIDING BACK-UP AND SHUT-DOWN MECHANISMS. NI EXPRESSLY DISCLAIMS ANY EXPRESS OR IMPLIED WARRANTY OF FITNESS OF THE PRODUCTS OR SERVICES FOR HIGH-RISK USES.

Contents

Chapter 1	
Getting Started	
Installation	1_1
Module Specifications	
wioduic specifications	1-1
Chapter 2	
Using the PXIe-4309	
Connecting Signals	2.1
Connecting Signals	
Analog Input Circuitry	
Module Pinout	
I/O Connector Signal Descriptions	
PXIe-4309 Block Diagram	
Signal Acquisition Considerations	
ADC	
Single Channel per ADC	
Multiple Channels per ADC	
Sample Modes	
Averaging Acquisition	
Single Convert Acquisition	
Sample Timing Modes	
Software-Timed Acquisitions	
Hardware-Timed Acquisitions	
Offset Error Cancellation Techniques	
Auto Zero	
Chopping	
Timing and Triggering	
Sample Clock Timebase	
External Clock	
Digital Triggering	
PFI	
Using PFI Terminals as Timing Input Signals	
Exporting Timing Output Signals Using PFI Terminals	
Connecting PFI Input Signals	
PFI Filters	2-19
I/O Protection	2-20
Synchronization	2-20
Reference Clock Synchronization	2-20
Accessory Auto-Detection	2-21

Chapter 3		
SC Expre	ss Considerations	
	lock and Trigger Signals	3-1
	K100	
_	Κ10	
_	gers	
	AR Trigger	
	TAR <a.c></a.c>	
Appendix	A	
NI Service	98	
Figures		
Figure 2-1.	Analog Input Circuitry	2-2
Figure 2-2.	Signal Source with Common-Mode Connected to an Analog Input	2.2
Figure 2-3.	Channel	
Figure 2-3.	Recommended Wiring for a Grounded Source with Common-Mode.	
Figure 2-5.	Recommended Wiring for a Grounded Source	
Figure 2-6.	PXIe-4309 Block Diagram	
Figure 2-7.	Single Channel Averaging	
Figure 2-8.	Multiple Channel Averaging	
Figure 2-9.	Single Convert Acquisition of a Single Channel	
Figure 2-10.	Single Convert Acquisition of Two Channels on the Same ADC	
Figure 2-11.	Analog Input Auto Zero	
Figure 2-12.	Analog Input Chopping	
Figure 2-13.	PFI Circuitry	
Figure 2-14.	Low-to-High Transition for a Custom Filter with N Set to 5	
Tables		
Table 2-1.	PXIe-4309 Front Connector Signal Pin Assignments	2-5
Table 2-2.	I/O Connector Signal Descriptions	2-6
Table 2-3.	AI Channel to ADC Mapping	2-8
Table 2-4.	Maximum Sample Rate	2-10
Table 2-5.	Maximum Sample Rate (Hardware-Timed Single Point,	
	On-Demand, and External Sample Clock)	
Table 2-6.	Maximum Sample Rate (Auto Zero Every Sample)	
Table 2-7.	AI Channel Connection Pairing for Chopping Channels	
Table 2-8.	Maximum Sample Rate (Chopping Enabled)	
Table 2-9.	Debouncing Filters	2-19
Table 3-1	PXIe DSTAR Line Descriptions	3-2

Getting Started

The PXIe-4309 is a flexible resolution voltage input module that supports up to 32 channels. The module multiplexes its analog input channels across eight simultaneously sampling ADCs.

The PXIe-4309 can sample eight input channels up to 2 MS/s/ch (single channel per ADC), and 32 input channels up to 400 kS/s aggregate (up to 4 channels per ADC).

Installation

Refer to the PXIe-4309 and TB-4309 (ST)/TB-4309 (MT) Getting Started Guide and Terminal Block Specifications document for step-by-step installation instructions, accessory details, and a complete listing of supported accessories.



Note For a complete list of terminal blocks supported by a specific release of NI-DAQmx, refer to the NI-DAQmx Readme, available on the version-specific download page or installation media.

Module Specifications

Refer to the *PXIe-4309 Specifications* document for module specifications.

Using the PXIe-4309

This chapter describes how to connect voltage signals to the PXIe-4309. It also provides the I/O connector pin assignments of the module.

Driver support for the PXIe-4309 was first available in NI-DAQmx 17.1. For the list of devices supported by a specific release, refer to the NI-DAOmx Readme, available on the version-specific download page or installation media.

Connecting Signals

This section provides information regarding the connection of voltage signals.



Caution To ensure the specified EMC performance, operate this product only with shielded, twisted pair cables, and shielded accessories.



Caution To ensure the specified EMC performance, the length of all I/O cables must be no longer than 3 m (10 ft.).



This icon denotes a caution, which advises you of precautions to take to avoid injury, loss of data, or a system crash.


Connecting Voltage Signals



Note The TB-4309 (ST) or TB-4309 (MT) terminal block is available for measuring voltage signals.

Analog Input Circuitry


Figure 2-1. Analog Input Circuitry

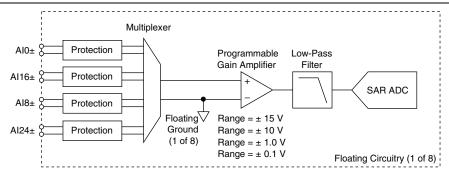

Input signals on each channel are conditioned, buffered, and then sampled by a SAR ADC.

The PXIe-4309 has eight independent signal conditioning paths and ADCs, enabling simultaneous sampling of up to eight input channels. Each ADC can scan through four independent input channels for a total of 32 input channels.

Figure 2-2. Signal Source with Common-Mode Connected to an Analog Input Channel

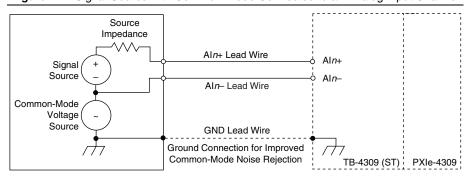

Source Impedance

Best measurement performance is achieved with source impedance of 50 Ω or lower; refer to the Source Impedance Error section of the PXIe-4309 Specifications for information on the error introduced by source impedances greater than 50 Ω .

Common-Mode Rejection Ratio (CMRR)

The analog input circuitry is floating to provide best in class CMRR of DC common-mode voltages.

Floating analog input circuits have an inherent input impedance mismatch. To minimize the effect of this mismatch on CMRR of AC common-mode voltages, keep the DC resistance and AC impedance of the AIn- lead wire to a minimum.

Common-Mode Noise

Enclosure

Floating power supplies inject common-mode current onto its floating ground. To minimize the effect of this current on measurement noise, keep the DC resistance and AC impedance of the AIn-lead wire to a minimum.

Shielding

For best noise immunity use double-shielded twisted pair cabling.

Aln+
AlnTB-4309 (ST)
PXIe-4309

Figure 2-3. Recommended Wiring for a Floating Source

Figure 2-4. Recommended Wiring for a Grounded Source with Common-Mode


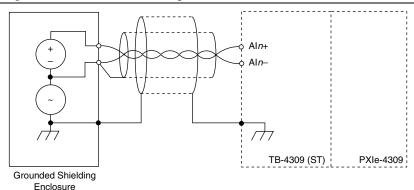


Figure 2-5. Recommended Wiring for a Grounded Source

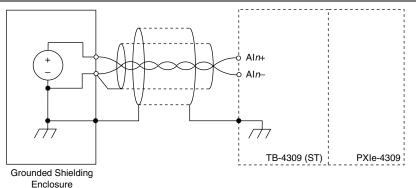

Module Pinout

Table 2-1 shows the pinout of the front connector of the PXIe-4309. Refer to the *I/O Connector* Signal Descriptions section for definitions of each signal. Refer to the PXIe-4309 and TB-4309 (ST)/TB-4309 (MT) Getting Started Guide and Terminal Block Specifications for signal locations on the terminal blocks.

Table 2-1. PXIe-4309 Front Connector Signal Pin Assignments

Front C	onne	ctor	Diagram	Row Number	Column A	Column B	Column C
				32	GND	AI16+	AI0+
	0	olum	ın.	31	AI8+	AI16-	AI0-
	Α	В	С	30	AI8-	AI24-	AI24+
32	0	0	0	29	GND	AI20+	AI4+
31	0	0	0	28	AI12+	AI20-	AI4-
30	0	0	0	27	AI12-	AI28-	AI28+
29	0	0	0	26	GND	AI21+	AI5+
28 27	0	0	0	25	AI13+	AI21-	AI5-
26	0	0	0	24	AI13-	AI29-	AI29+
25	0	0	0	23	GND	AI17+	AI1+
24	0	0	0	22	AI9+	AI17-	AI1-
23 22	0	0	0	21	AI9-	AI25-	AI25+
21	0	0	0	20	GND	NC	NC
20	0	0	0	19	REF+	NC NC	NC NC
19	0	0	0	-			
18	0	0	0	18	REF-	NC	NC
17	0	0	0	17	GND	AI18+	AI2+
16	0	0	0	16	AI10+	AI18-	AI2-
15 14	0	0	0	15	AI10-	AI26-	AI26+
13	0	0	0	14	GND	AI22+	AI6+
12	0	0	0	13	AI14+	AI22-	AI6-
11	0	0	0	12	AI14-	AI30-	AI30+
10 9	0	0	0	11	GND	AI23+	AI7+
8	0	0	0	10	AI15+	AI23-	AI7-
7	0	0	0	9	AI15-	AI31-	AI31+
6	0	0	0	8	GND	AI19+	AI3+
5	0	0	0	7	AI11+	AI19-	AI3-
4	0	0	0				
3	0	0	0	6	AI11-	AI27-	AI27+
2	0	0	0	5	PFI GND	PFI0	PFI1
'	Ľ			4	NC	GND	GND
				3	NC	NC	NC
RS	VD is	rese	rved	2	RSVD	RSVD	RSVD
			ection	1	RSVD	RSVD	RSVD

I/O Connector Signal Descriptions

Table 2-2 describes the signals found on the I/O connector.

Table 2-2. I/O Connector Signal Descriptions

Signal Names	Direction	Description
GND	_	Ground
AI<031>+	Input	Positive inputs of the differential analog input channels 0 to 31.
AI<031>-	Input	Negative inputs of the differential analog input channels 0 to 31.
RSVD	_	Reserved for communication with the accessory.
PFI<01>	Input/Output	Configurable 5.0 V tolerant digital input or 3.3 V digital output for sending or receiving trigger and synchronization signals. PFI lines are referenced to PFI GND.
PFI GND	_	Digital ground reference for PFI lines.
REF+	Output	Positive output of the internal reference used for self-calibration. Refer to the <i>PXIe-4309</i> Specifications and <i>PXIe-4309 Calibration</i> Procedure.
REF-	Output	Negative output of the internal reference used for self-calibration. Refer to the PXIe-4309 Specifications and PXIe-4309 Calibration Procedure.

PXIe-4309 Block Diagram

Input signals on each channel are conditioned, buffered, and then sampled by a SAR ADC. The PXIe-4309 supports averaging and auto zero in hardware.

The PXIe-4309 has 8 independent signal conditioning paths and ADCs, enabling simultaneous sampling of up to eight input channels. Each ADC can scan through four independent input channels for a total of 32 input channels. Figure 2-6 shows the block diagram of the PXIe-4309.

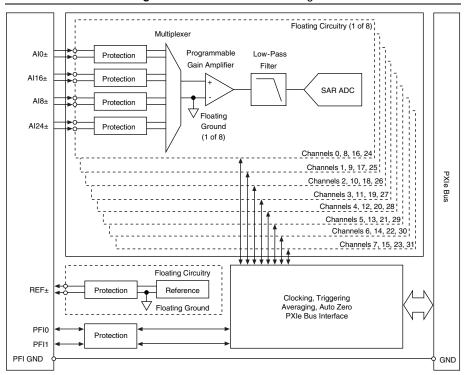


Figure 2-6. PXIe-4309 Block Diagram

Signal Acquisition Considerations

This section contains information about signal acquisition concepts, including operation modes, SAR converters, timing, triggering, and synchronization.

ADC

The PXIe-4309 has eight simultaneously sampled 18-bit 2 MS/s successive approximation (SAR) ADCs. SAR involves a binary search through all possible quantization levels before finally converging upon a digital output for each conversion.

Single Channel per ADC

The PXIe-4309 can sample eight simultaneous input channels up to 2 MS/s when only one channel per ADC is selected.

Multiple Channels per ADC

The PXIe-4309 can sample 32 input channels when multiple channel per ADC are selected. In this mode of operation, each ADC can scan through four independent input channels up to an aggregate sample rate of 400 kS/s. The mapping of input channels to ADCs is shown in Table 2-3

ADC	Al Channels
0	ai0, ai8, ai16, ai24
1	ai1, ai9, ai17, ai25
2	ai2, ai10, ai18, ai26
3	ai3, ai11, ai19, ai27
4	ai4, ai12, ai20, ai28
5	ai5, ai13, ai21, ai29
6	ai6, ai14, ai22, ai30
7	ai7, ai15, ai23, ai31

Table 2-3. Al Channel to ADC Mapping

The signal conditioning circuitry of each ADC has a non-zero settling time when it changes the input multiplexer to scan through multiple channels.

The settling time increases as the aggregate sample rate decreases to allow for more settling time. This reduces settling and ghosting errors.

Sample Modes

The PXIe-4309 has two sampling modes: averaging and single sample. The sampling mode is selected based on the timing mode. Refer to Sample Timing Modes section for more information.

Terminology for Averaging Acquisition and Single Convert Acquisition sections:


- Convert is a single ADC quantized digitization of the input signal.
- Sample is the data returned from hardware.

Averaging Acquisition

The PXIe-4309 supports averaging in hardware. Based on the sample rate, the PXIe-4309 will average as many converts per sample as possible to deliver the best resolution for a given sample rate. Refer to the PXIe-4309 Specifications for Noise versus Sample Rate information.

Figure 2-7 depicts an averaging acquisition of two samples of a single channel.

Figure 2-7. Single Channel Averaging

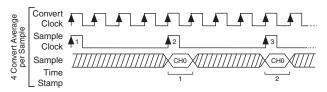



Samples are timestamped to the time when the sample clock takes place, but are acquired throughout the sample period. This introduces a phase shift of half a sample period when channel expanding to other types of data acquisition devices.

Figure 2-8 depicts an averaging acquisition of two converts per sample of two channels on the same ADC. The settling time increases as the aggregate sample rate decreases to allow for more settling time to reduce settling and ghosting errors.

Figure 2-8. Multiple Channel Averaging

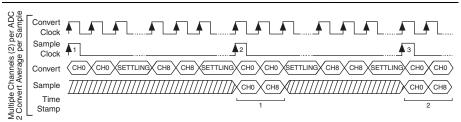


Table 2-4 shows the max sampling rate based on the maximum number of channels per ADC.

Table 2-4. Maximum Sample Rate

Maximum Channels per ADC	Maximum Sample Rate	Aggregate Rate
1	2 MHz	2 MHz
2	200 kHz	
3	133 kHz	400 kHz
4	100 kHz	

Single Convert Acquisition

The PXIe-4309 takes a single convert per channel and returns that as a sample. Figure 2-9 depicts a single convert acquisition of a single channel.

Figure 2-9. Single Convert Acquisition of a Single Channel

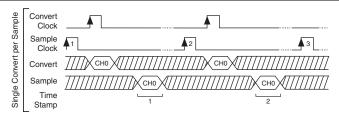


Figure 2-10 depicts a single convert acquisition of two channels on the same ADC. The settling time increases as the aggregate sample rate decreases to allow for more settling time to reduce settling and ghosting errors.

Figure 2-10. Single Convert Acquisition of Two Channels on the Same ADC

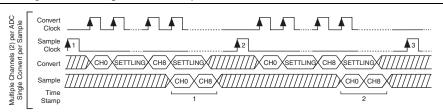


Table 2-5 shows the max sampling rate based on the maximum number of channels per ADC. Table 2-5. Maximum Sample Rate (Hardware-Timed Single Point, On-Demand, and External Sample Clock)

Maximum Channels per ADC	Maximum Sample Rate	Aggregate Rate
1	900 kHz	900 kHz
2	180 kHz	
3	120 kHz	360 kHz
4	90 kHz	

Sample Timing Modes

The PXIe-4309 supports four timing modes that fall into two categories: software-timed and hardware-timed.



Note Software-timed, hardware timed single point (HWTSP), and external sample clock use single sample acquisition to minimize latency. Refer to the Single Convert Acquisition section for more information.

Software-Timed Acquisitions

Software controls the rate of the acquisition. Software sends a separate command to the hardware to acquire each sample. In NI-DAOmx, software-timed acquisitions are referred to as having on-demand timing. Software-timed acquisitions are also referred to as immediate or static acquisitions and are typically used for reading a single sample of data.

Hardware-Timed Acquisitions

A digital hardware signal (AI Sample Clock) controls the rate of the acquisition. This signal can be generated internally by the PXIe-4309 or provided externally.

Hardware-timed acquisitions have several advantages over software-timed acquisitions:

- The time between samples can be much shorter.
- The timing between samples is deterministic.
- Hardware-timed acquisitions can use hardware triggering.

Hardware-timed operations can be buffered or HWTSP. A buffer is a temporary storage in computer memory for to-be-transferred samples.

- **Buffered**—In a buffered acquisition, data is moved from the onboard FIFO memory of the DAQ device to a PC buffer using DMA before it is transferred to application memory. Buffered acquisitions typically allow for much faster transfer rates than HWTSP acquisitions because data is moved in large blocks, rather than one point at a time.
 - One property of buffered I/O operations is the sample mode. The sample mode can be either finite or continuous:
 - Finite sample mode acquisition refers to the acquisition of a specific, predetermined number of data samples. Once the specified number of samples has been read, the acquisition stops. If you use a reference trigger, you must use finite sample mode.
 - Continuous acquisition refers to the acquisition of an unspecified number of samples. Instead of acquiring a set number of data samples and stopping, a continuous acquisition continues until you stop the operation. Continuous acquisition is also referred to as double-buffered or circular-buffered acquisition.
 - If data cannot be transferred across the bus fast enough, the FIFO becomes full. For more information refer to NI-DAQmx Help.

HWTSP—Typically used to read single samples at known time intervals. While buffered operations are optimized for high throughput, HWTSP operations are optimized for low latency and low jitter. In addition, HWTSP can notify software if it falls behind hardware. These features make HWTSP ideal for real time control applications. HWTSP operations, in conjunction with the wait for next sample clock function, provide tight synchronization between the software layer and the hardware layer. Refer to the NI-DAOmx Hardware-Timed Single Point Lateness Checking document, for more information. To access this document, go to ni.com/info and enter the Info Code daghwtsp.

Offset Error Cancellation Techniques

The PXIe-4309 supports two types of offset error cancellation techniques: auto zero and chopping.

Auto Zero

Each signal conditioning path of the PXIe-4309 has an internal channel to measure its offset error as shown in Figure 2-11. When Auto Zero is set to once or every sample, the offset error is measured and subtracted from the input measurement, effectively canceling the offset error.

Multiplexer AI0+ Protection AIO-Programmable Gain Amplifier AI16+ Protection AI16-AI8+ Protection Floating Range = \pm 15 V Ground Range = \pm 10 V AI24+ (1 of 8) Range = \pm 1.0 V Protection AI24-Range = $\pm 0.1 \text{ V}$ Short < Floating Circuitry (1 of 8)

Figure 2-11. Analog Input Auto Zero

Use the DAQmx Channel Property Node (Analog Input»General Properties»Advanced» High Accuracy Settings» Auto Zero Mode) to configure Auto Zero behavior as none, once, or every sample.



Note To read the offset error add _aignd_vs_aignd<0..7> to the acquisition task, where each channel index represents the ADC index. The offset error can only be read with Auto Zero set to none.



Note For a single PXIe-4309 module, the Auto Zero setting must be consistent across all channels for a given task; for channel expansion across multiple PXIe-4309, each module can have a different Auto Zero setting.



Note When Auto Zero is set to Once or Every Sample, the Input Voltage Range setting must be consistent across all channels per ADC for a given task.

Auto Zero Once

With Auto Zero set to once, the offset error is measured when the acquisition task starts. The offset error measurement introduces a delay from start to the first sample clock. The start delay may be gueried with the DAOmx Trigger Property Node (Start»More»Delay).



Note The start delay is in units of 100 MHz timebase ticks by default.

With auto zero set to once, the offset error temperature coefficient is the same as with auto zero none, making Auto Zero once an effective offset cancellation technique so long as the ambient temperature is stable throughout an acquisition.



Note Auto Zero Once has no impact on low frequency noise performance.

Refer to the *PXIe-4309 Specifications* for performance specifications of Auto Zero Once.

Auto Zero Every Sample

With Auto Zero set to every sample, the offset error is measured at the beginning of every sample period. Refer to the PXIe-4309 Specifications for performance specifications of Auto Zero Every Sample.

With Auto Zero Every Sample the offset error temperature coefficient is reduced, making auto zero every sample an effective offset cancellation technique in the presence of ambient temperature fluctuations.



Note Auto Zero Every Sample improves low frequency noise performance.

Refer to the PXIe-4309 Specifications for performance specifications of Auto Zero Every Sample.

Table 2-6 shows the maximum sample rate with Auto Zero set to every sample based on the maximum number of channels per ADC.

Maximum Channels per ADC	Maximum Sample Rate	Aggregate Rate
1	10 kHz	
2	5 kHz	10.1.11
3	3.3 kHz	10 kHz
4	2.5 kHz	

Table 2-6. Maximum Sample Rate (Auto Zero Every Sample)

Chopping

Figure 2-12 shows the connections required to use two input channels of the same ADC to enable chopping for the measurement of two voltage sources.

Multiplexer AI0+ Protection AI0-Signal V_{S1} Source Programmable AI16+ Gain Amplifier Protection AI16-AI8+ Protection Signal Floating Range = \pm 15 V AI8-Source Ground Range = \pm 10 V AI24+ (1 of 8)Range = \pm 1.0 V Protection Range = $\pm 0.1 \text{ V}$ AI24-Floating Circuitry (1 of 8)

Figure 2-12. Analog Input Chopping

The first channel measures the source voltage, while the second channel measures the negative source voltage; in both measurements the offset error of the signal conditioning path adds to the measured voltage as follows:

$$V_{AI0} = V_{S1} + V_{O}$$

$$V_{AI16} = -V_{S1} + V_{O}$$

Using these measurements it is possible to cancel out the offset error of the signal conditioning path.

$$0.5 \times (V_{AI0} - V_{AI16}) = 0.5 \times ([V_{S1} + V_{O}] - [-V_{S1} + V_{O}]) = 0.5 \times 2 V_{S1} = V_{S1}$$

The PXIe-4309 performs these mathematical operations in hardware when chopping is enabled. Refer to the *PXIe-4309 Specifications* for performance specifications of Chopping.

Table 2-7 shows the AI channel connections that are paired to create a chopping channel.

Table 2-7. Al Channel Connection Pairing for Chopping Channels

Chopping Channel	Al Channel
ai0	ai0, ai16
ai1	ail, ai17
ai2	ai2, ai18
ai3	ai3, ai19
ai4	ai4, ai20
ai5	ai5, ai21
ai6	ai6, ai22
ai7	ai7, ai23
ai8	ai8, ai24
ai9	ai9, ai25
ai10	ai10, ai26
ai11	ai11, ai27
ai12	ai12, ai28
ai13	ai13, ai29
ai14	ai14, ai30
ai15	ai15, ai31

Table 2-8 shows the maximum sample rate when chopping is enabled based on the maximum number of chopping channels per ADC.

Table 2-8. Maximum Sample Rate (Chopping Enabled)

Maximum Chopping Channels per ADC	Maximum Sample Rate	Aggregate Rate
1	10 kHz	
2	5 kHz	10 kHz

Use the DAQmx Channel Property Node (Analog Input»General Properties»Advanced» High Accuracy Settings» Chop» Enable) to enable or disable chopping. For a LabVIEW example to enable chopping, refer to ni.com/info and enter the Info Code chopexp.



Note Chopping, Auto Zero Once, and Auto Zero Every Sample are mutually exclusive configurations.

Timing and Triggering

This section contains information about timing and triggering for the PXIe-4309.

Sample Clock Timebase

The PXIe-4309 sample clock is derived from the onboard 100 MHz clock, which can be locked to the PXIe backplane 100 MHz clock to synchronize multiple modules.

External Clock

The PXIe-4309 sample clock can be driven from external sources by either PXI TRIG<0..7>, PFI<0..1>, PXI STAR, or PXIe DSTAR<A,B>.



Note When using the external sample clock, the PXIe-4309 does not support averaging, auto zero, or chopping.

Digital Triggering

The PXIe-4309 supports acquisition start in response to a digital trigger signal from one of the PXI Express backplane trigger lines or the PFI from the front connector. The trigger circuit can respond either to a rising or a falling edge.

PFI


Each PFI can be individually configured as the following:

- A timing input signal for AI.
- A timing output signal from AI.

Each PFI input has a programmable debouncing filter.

Figure 2-13 shows the PFI circuitry.

Figure 2-13. PFI Circuitry

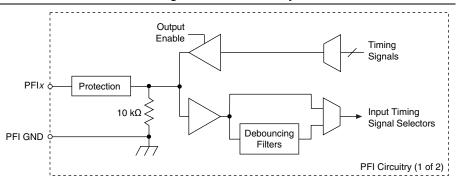

On the I/O connector, each terminal is labeled PFIx.

Refer to the PXIe-4309 Specifications for the PFI input and output specifications.

Using PFI Terminals as Timing Input Signals

Use PFI terminals to route external timing signals to different functions of the PXIe 4309. Each PFI input can be routed to any of the following signals:

- AI Sample Clock (ai/SampleClock)
- AI Start Trigger (ai/StartTrigger)
- AI Reference Trigger (ai/ReferenceTrigger)
- AI Pause Trigger (ai/PauseTrigger)
- AI Sample Clock Timebase (ai/SampleClockTimebase)



Note The polarity of PFI inputs is configurable for most functions.



Note Pause triggers are only sensitive to the level of the source, not the edge.

Exporting Timing Output Signals Using PFI Terminals

You can route any of the following timing signals to any PFI terminal configured as an output:

- AI Reference Trigger (ai/ReferenceTrigger)
- AI Sample Clock (ai/SampleClock)
- AI Start Trigger (ai/StartTrigger)
- AI Pause Trigger (ai/PauseTrigger)
- PXI STAR
- PXI Trig <0..7>

Connecting PFI Input Signals

All PFI I/O lines are referenced to PFI GND.

PFI Filters

The PXIe-4309 has programmable debouncing filters on each PFI, PXI Trig, PXI STAR, or PXI DSTAR. When filters are enabled, the PXIe-4309 uses the onboard oscillator as the filter clock to sample the filter input on each rising edge and to generate the filtered digital signal.



Note Enabling filters introduces jitter on the input signal. The maximum jitter is one period of the filter clock.

Figure 2-14 shows an example of a low-to-high transition for a custom filter with N set to 5.

Figure 2-14. Low-to-High Transition for a Custom Filter with N Set to 5 Filter Input _____ Filtered input goes high when terminal is sampled high on five consecutive filter clocks Filtered Digital Signal _

Table 2-9 shows the debouncing filters supported by the PXIe-4309.

N (Filter **Pulse Width** Clocks **Pulse Width** Guaranteed Needed to Guaranteed to Not Pass Filter Setting Filter Clock Pass Signal) to Pass Filter Filter None 100 MHz 9 90 ns 90 ns 80 ns (short) $5.12 \mu s$ 100 MHz 512 $5.12 \mu s$ 5.11 µs (medium) (N - 1)100 MHz N Custom 100 MHz 100 MHz

Table 2-9. Debouncing Filters

The filter setting for each input can be configured independently. Filters are disabled on power up.

I/O Protection

Each PFI signal is protected against overvoltage, undervoltage, and overcurrent conditions as well as ESD events; however, avoid these fault conditions by following these guidelines:

- When using PFI lines as an output, do *not* connect them to any external signal source, ground, or power supply; understand the current requirements of the load connected to these lines; do not exceed the specified current drive of the device. NI has several signal conditioning solutions for digital applications requiring high current drive.
- When using PFI lines as an input, do not drive the line with voltages outside of its normal operating range. The PFI lines have a smaller operating range than the analog input channels.
- Treat the device as you would treat any static sensitive device. Always properly ground yourself and the equipment when handling the device or connecting to it.

Synchronization

Some applications require tight synchronization between input and output operations on multiple modules. Synchronization is important to minimize skew between channels and to eliminate clock drift between modules in long-duration operations. You can synchronize the analog input operations on two or more PXIe-4309 modules to extend the channel count for your measurements. In addition, the PXIe-4309 can synchronize with 62XX (M Series) and 63XX (X Series) devices/modules, such as the PXI-6289 and PXIe-6363, respectively. Refer to the Reference Clock Synchronization section for more information.

Reference Clock Synchronization

With reference clock synchronization, master and slave modules generate their ADC sample clock from the shared 100 MHz reference clock from the PXIe backplane (PXIe CLK100). The backplane supplies an identical copy of this clock to each peripheral slot. In addition, multiple chassis can be synchronized by using a timing and synchronization board to lock the 100 MHz clock across chassis.

When acquiring data from multiple modules within the same NI-DAQmx task, NI-DAQmx will automatically handle all of the Reference Clock Synchronization details required to synchronize the modules within the task. This is known as a multi-device task.

To perform reference clock synchronization when using multiple NI-DAQmx tasks that are acquiring at the same rate, complete the following steps to synchronize the hardware.

- 1 Specify PXIe CLK100 as the reference clock source for all modules to force all the modules to lock to the reference clock on the PXIe chassis.
- 2. Choose one module to be the start trigger master.
- Configure the rest of the modules in your system to receive their start trigger from the start trigger master module. This ensures that all modules will begin returning data on the same sample.

- Set the synchronization type of the Start Trigger slaves at **DAQmx Trigger»Advanced» Synchronization**»**Synchronization Type** to Slave and that of the Master to Master.
- Start all of the start trigger slave module tasks. This sets them up to expect the start trigger from the master.
- Start the start trigger master module task. You can now acquire data.



Tip Consider using a multi-device task when synchronizing multiple devices at the same rate



Tip You can find example VIs in the NI Example Finder. Select Help»Find **Examples** to launch the NI Example Finder.

Accessory Auto-Detection

SC Express modules automatically detect compatible accessories or terminal blocks. The RSVD pins on the I/O connector provide power to the accessories as well as digital communication lines. This allows software to detect when accessories are inserted or removed. In addition, software can automatically identify the specific terminal block as well as access any calibration or scaling information associated with the terminal block.

MAX allows you to see which accessories are currently connected to your module. In MAX, expand Devices and Interfaces and locate your module. If a terminal block is connected to your module, it will be displayed beneath the module. Unsupported terminal blocks appear in MAX with an X next to them.

NI-DAQmx property nodes can be used to programmatically access information about connected accessories in your application. Refer to the NI-DAOmx Help for documentation about programmatically accessing accessory status.

SC Express Considerations

This chapter details the clock and trigger functionality available through the PXI Express chassis

SC Express Clock and Trigger Signals

PXIe_CLK100

PXIe CLK100 is a common, low-skew 100 MHz reference clock used for synchronization of multiple modules in a PXI Express measurement or control system. The PXIe backplane is responsible for generating PXIe CLK100 independently to each peripheral slot in a PXI Express chassis. For more information, refer to the PXI Express Specification at www.pxisa.org.

PXI CLK10

PXI CLK10 is a common, low-skew 10 MHz reference clock for synchronization of multiple modules in a PXI measurement or control system. The PXI backplane is responsible for generating PXI CLK10 independently to each peripheral slot in a PXI chassis. In PXIe chassis, the PXI CLK10 signal is in phase with PXIe CLK100.



Note PXI CLK10 cannot be used as a reference clock for SC Express modules.

PXI Triggers

A PXI/PXIe chassis provides eight bused trigger lines to each module in a system. Triggers may be passed from one module to another, allowing precisely timed responses to asynchronous external events that are being monitored or controlled. Triggers can be used to synchronize the operation of several different PXI peripheral modules.

In a PXI chassis with more than eight slots, the PXI trigger lines may be divided into multiple independent buses. Refer to the documentation for your chassis for details.

PXI STAR Trigger

In a PXI Express system, the Star Trigger bus implements a dedicated trigger line between the system timing slot and the other peripheral slots. The Star Trigger can be used to synchronize multiple modules or to share a common trigger signal among modules.

A system timing controller can be installed in the system timing slot to provide trigger signals to other peripheral modules. Systems that do not require this functionality can install any standard peripheral module in this system timing slot.

An SC Express module receives the Star Trigger signal (PXI STAR) from a System timing controller. PXI STAR can be used as a trigger signal for input operations.

An SC Express module is not a System timing controller. An SC Express module can be used in the system timing slot of a PXI system, but the system will not be able to use the Star Trigger feature.

PXIe DSTAR<A..C>

PXI Express devices can provide high-quality and high-frequency point-to-point connections between each slot and a system timing slot. These connections come in the form of three low-voltage differential star triggers that route between a PXI Express system timing controller and a peripheral device. Using multiple connections simplifies the creation of applications because of the increased routing capabilities.

Table 3-1 describes the three differential star (DSTAR) lines and how they are used.

Trigger Line	Purpose
PXIe_DSTARA	Distributes high-speed, high-quality clock signals from the system timing slot to the peripherals (input).
PXIe_DSTARB	Distributes high-speed, high-quality trigger signals from the system timing slot to the peripherals (input).
PXIe_DSTARC	Sends high-speed, high-quality trigger or clock signals from the peripherals to the system timing slot (output).

Table 3-1. PXIe DSTAR Line Descriptions

The DSTAR lines are only available for PXI Express devices when used with a PXI Express system timing module. For more information, refer to the PXI Express Specification at www.pxisa.org.



NI Services

NI provides global services and support as part of our commitment to your success. Take advantage of product services in addition to training and certification programs that meet your needs during each phase of the application life cycle; from planning and development through deployment and ongoing maintenance.

To get started, register your product at ni.com/myproducts.

As a registered NI product user, you are entitled to the following benefits:

- Access to applicable product services.
- Easier product management with an online account.
- Receive critical part notifications, software updates, and service expirations.

Log in to your MyNI user profile to get personalized access to your services.

Services and Resources

- Maintenance and Hardware Services—NI helps you identify your systems' accuracy and reliability requirements and provides warranty, sparing, and calibration services to help you maintain accuracy and minimize downtime over the life of your system. Visit ni.com/ services for more information.
 - Warranty and Repair—All NI hardware features a one-year standard warranty that is extendable up to five years. NI offers repair services performed in a timely manner by highly trained factory technicians using only original parts at a NI service center.
 - Calibration—Through regular calibration, you can quantify and improve the measurement performance of an instrument. NI provides state-of-the-art calibration services. If your product supports calibration, you can obtain the calibration certificate for your product at ni.com/calibration.
- **System Integration**—If you have time constraints, limited in-house technical resources, or other project challenges, National Instruments Alliance Partner members can help. To learn more, call your local NI office or visit ni.com/alliance.

- Training and Certification—The NI training and certification program is the most effective way to increase application development proficiency and productivity. Visit ni.com/training for more information.
 - The Skills Guide assists you in identifying the proficiency requirements of your current application and gives you options for obtaining those skills consistent with your time and budget constraints and personal learning preferences. Visit ni.com/skills-quide to see these custom paths.
 - NI offers courses in several languages and formats including instructor-led classes at facilities worldwide, courses on-site at your facility, and online courses to serve your individual needs.
- Technical Support—Support at ni.com/support includes the following resources:
 - Self-Help Technical Resources—Visit ni.com/support for software drivers and updates, a searchable KnowledgeBase, product manuals, step-by-step troubleshooting wizards, thousands of example programs, tutorials, application notes, instrument drivers, and so on. Registered users also receive access to the NI Discussion Forums at ni.com/forums. NI Applications Engineers make sure every question submitted online receives an answer.
 - Software Support Service Membership—The Standard Service Program (SSP) is a renewable one-year subscription included with almost every NI software product, including NI Developer Suite. This program entitles members to direct access to NI Applications Engineers through phone and email for one-to-one technical support, as well as exclusive access to online training modules at ni.com/self-paced-training. NI also offers flexible extended contract options that guarantee your SSP benefits are available without interruption for as long as you need them. Visit ni.com/ssp for more information.
- **Declaration of Conformity (DoC)**—A DoC is our claim of compliance with the Council of the European Communities using the manufacturer's declaration of conformity. This system affords the user protection for electromagnetic compatibility (EMC) and product safety. You can obtain the DoC for your product by visiting ni.com/certification.

For information about other technical support options in your area, visit ni.com/services, or contact your local office at ni.com/contact.

You also can visit the Worldwide Offices section of ni.com/niglobal to access the branch office websites, which provide up-to-date contact information, support phone numbers, email addresses, and current events.